

**WHEELED COACH INDUSTRIES**

2737 N. FORSYTH ROAD

WINTER PARK, FL 32792

(800) 932-7077



**ELECTRONIC MANUAL**



Introduction



Safety Standards



Operating Instructions



Maintenance Information



Warranty and Coverage



Warranty Claim Form



Parts and EMS Supplies



Service Centers



Work Order



Electrical Drawings



---

Corporate Office: 2737 North Forsyth Road, Winter Park, FL 32792 / Post Office Box 677339, Orlando, FL 32867-7339  
800-342-0720 ◊ (407) 677-7777 ◊ Fax: (407) 679-1337 ◊ <http://www.wheeledcoach.com>

#### To Our Customers:

Congratulations for choosing *Wheeled Coach*; you have chosen the best! Your vehicle represents the most innovative and sophisticated piece of equipment of its kind. We are proud to say that every effort to ensure a quality product was incorporated into the manufacturing process.

As with any well-designed system, proper operation and maintenance is vital to ensure continued trouble free service. We designed this manual to be used as the primary source of reference for general operation and maintenance of your ambulance.

Refer to this manual if you have any questions or problems with your vehicle. If you need further assistance, contact either a service center in your area or call our *convenient 24-hour service line at: 800-628-8178*. All *Wheeled Coach* distributors and district sales managers are ready and eager to help you in any way they can.



## **Chassis, Engine and Warranty Information**

---

### **Important Do Not Violate Your Manufacturer's Warranties**

Your vehicle comes with a manufacturer's information packet provided by the chassis and engine manufacturers. If you haven't already done so, take time to go through the packet to see what is included. Please refer to your manufacturer's packet for any questions concerning chassis, engine or components which came from the manufacturer as part of the chassis and engine assembly.

If you need assistance from Chevrolet or Ford, and the dealership cannot assist you, please call the following numbers:

**Chevrolet:** 1-800-243-8872

**Ford:** 1-800-392-3673      USA  
1-800-565-3673      Canada  
1-313-317-4282      Outside the USA or Canada

**Freightliner:** 1-800-385-4357

**Navistar International:** 1-800-385-4357



## SAFETY STANDARDS



### National Highway Traffic Safety Administration

One mission of the National Highway Traffic Safety Administration (NHTSA) is to monitor the safety of all types of automotive vehicles. Our goal in producing this manual is to ensure safe and dependable operation of our vehicles. Also, in support of the NHTSA's efforts to make our roads and highways safer, we would like to point out the following NHTSA procedures:

*If you believe that your vehicle has a defect which could cause a crash, injury, or death, you should immediately inform the NHTSA, in addition to notifying Wheeled Coach Industries.*

*If the NHTSA receives similar complaints, it may open an investigation, and if it finds that a safety defect exists in a group of vehicles, it may order a recall and remedy campaign. However, the NHTSA cannot become involved in individual problems between you, your dealer, or Wheeled Coach Industries.*

*To contact the NHTSA, you may call the Auto Safety Hotline, toll-free, at 1-888-327-4236 (202-366-0123 in the Washington, D.C. area), or write to NHTSA, Office of the Administrator (NAO-10) 400 Seventh Street S.W. Washington, D.C. 20590. You can also obtain other information about motor vehicle safety from the hotline.*



## **Safety Standards**

---

The following section reflects our interest in the safe operation of our vehicles. We urge all operators, drivers, and attendants to review these tips to make your job safer for yourself, your patients and your passengers.

### **Inspection and Maintenance Safety**

- ◆ Read chassis manufacturer's literature.
- ◆ Check brake shoes and pads for necessary replacement.
- ◆ Check cooling system under pressure to spot any high pressure leaks. Inspect radiator fan and shroud for proper operation.
- ◆ Inspect fuel system regularly.
- ◆ Check lug nuts for torque as specified in Owner's Guide.
- ◆ Check high idle operator for correct settings and operation.
- ◆ Use the chassis maintenance schedule for "Unique Operating Conditions."
- ◆ Avoid after-market fuel system components.
- ◆ Keep heat shields intact and in place.
- ◆ Keep undercoating at least 10 inches clear of exhaust system components.
- ◆ Keep front of engine clear; do not mount anything in front of radiator. Do not modify radiator fan or shroud.
- ◆ Replace all-season coolants at least once a year.
- ◆ Inspect Engine air intake and exhaust systems regularly.
- ◆ Use the instructions of your chassis manufacturer's for cleaning all glass (FMVSS 205)
- ◆ Seats not intended for passenger transport are marked accordingly in your ambulance (FMVSS 207)
- ◆ Child restraints have operating instructions attached, use them. (FMUSS 210)



## **Safety Standards**

---

- ◆ During the winter, run the air conditioners for 5 minutes each week and inspect for function, leaks, etc.
- ◆ Clean any combustible materials from the engine compartment each month.

## **Fuel Safety**

- ◆ Fill the tanks only to the automatic nozzle's first shut-off point. This will assure adequate vapor expansion space in the tanks.
- ◆ Fuel your vehicles on level ground to avoid overfilling tanks.
- ◆ Use approximately two gallons of fuel from each tank immediately after filling.
- ◆ **Do not remove** the fuel cap if there is a hissing sound coming from the cap, or if fuel can be seen seeping out. These are signs that the cap is venting pressure from the tank, and removing it could cause fuel to spray out.
- ◆ **Do not overfill transmission fluid.** Transmission fluid expands when heated and could cause a fire if it is spilled onto a hot exhaust system.

## **Operations Safety**

- ◆ Drivers should be trained in driving techniques for emergency vehicles.
- ◆ Use vehicle only within the Gross Vehicle Weight Rating (GVWR) and Gross Axle Weight Rating (GAWR).
- ◆ Avoid parking the vehicle over combustible materials.



## **Safety Standards**

---

- ♦ During movement of the vehicle, be sure:

*All passengers have their seat belts securely fastened.*

*All cots are secured with their cot mounts.*

*Patients are secured to their cots in accordance with cot manufacturer's instructions.*

*A fully charged fire extinguisher is on board.*

*Never drive faster than road conditions allow. It is far better to arrive at the scene dispatched to than create a second scene!*

*Warranty registration serves as the link between you and any potential recalls. Be sure to sign and return or you stand the chance of not being notified of potential safety defects in your ambulance (FMVSS Part 574)*



## OPERATING INSTRUCTIONS



## Operating Instructions

---

### Starting the Vehicle

1. The Master Switch on the driver's console must be in the **OFF** position before starting the engine. Move it to **OFF** if it is not.
2. The switches controlling all electrical equipment on both the driver's and attendant's consoles should be in the **OFF** position before starting the engine. Move them to **OFF** if they are not.

*Starting the engine with the electrical equipment drawing current places unnecessary stress on the charging system.*

3. The green indicator will come on with the ignition.

If you believe that one of the batteries is defective or discharged, we recommend that you keep the switch in the **OFF** position and not start the vehicle until the condition is corrected. Only if it is absolutely necessary should you start and use the vehicle with a defective battery and then only for the shortest time possible.

*A defective or discharged battery will drain charge from the good battery and could prevent the vehicle from starting.*

4. Start the engine. Refer to the manufacturer's manual for starting instructions.



## Operating Instructions

---

### **Turning on the Electrical Systems** (Patient Support, Emergency Signals, etc.)

1. With the engine running, move the Master Switch on the driver's control console to **ON**.
2. You may turn on the equipment switches in any sequence. However, switches must be turned on one at a time, with a brief pause between turning them on. (An optional sequencing device can be installed which automatically activates each pre-selected switch, one at a time, avoiding the undesirable surge which can result from multiple or simultaneous switching.)

*Never load the system by turning on more than one switch at a time. A sudden demand for power can cause the alternator belt to slip on its pulley and generate a damaging electrical spike in the voltage regulator. This will reduce the life of the belt and the charging system.*

### **Stationary Vehicle Operation**

If the vehicle is to remain stationary, engage the Electric Fast Idle Throttle.

*This throttle increases the alternator output which helps the system to support the electrical load. Some manufacturers require the emergency brake to be on to operate the throttle system. Please consult your chassis manufacturer's operating manual.*

**Watch the ammeter!** If it shows discharge, turn off as many electrical accessories as possible to prevent system failure. Do not allow the ammeter to remain in the discharge (-) region.

*The ammeter will read discharge (-) if there is too large a load for the alternator to handle. Reduce the load to prevent electrical system failure due to excessive battery drain.*



## **Operating Instructions**

---

### **Systems Operations**

NOTE: To aid in identifying switch location, switches will be labeled as (D), located on the Driver's console, or (A), located on the Attendant's console.

#### **Action Area Lighting**

IGN	ON
Master Power Switch (D)	ON
Action Area Light Switch (A)	ON

#### **Air Conditioning/Heater, Patient Compartment**

IGN	ON
Master Power Switch (D)	ON
Patient Heat/AC Switch (D) & (A)	ON

To adjust compartment temperature:

Adjust temperature control to desired level.

Select desired fan speed.

#### **Air Conditioner/Heater, Patient Compartment, Split System**

IGN	ON
Master Power Switch (D)	ON
Patient Heat/AC Switch (D) & (A)	ON

To adjust compartment temperature:

Adjust temperature control to desired level.

Select desired fan speed.



## **Operating Instructions**

---

### **Back-Up Alarm**

Placing vehicle in reverse gear automatically activates the Back-Up Alarm.

### **Beacon (see Light Bar)**

### **Control Console, Attendant**

All functions within the patient compartment are controlled at this console. The battery switch (if equipped) and the Master Power Switch must be ON for the console to be operable.

### **Control Console, Driver**

The instrument panel light intensity, switch pilot lights, and the voltmeter and ammeter lights are controlled by the vehicle's headlight switch.

### **Dome Lights, Left**

Depending on the compartment configuration, between three and five dual intensity dome lights on the left side (as one faces the front of the vehicle) of the compartment illuminate the primary patient area. The lights are automatically activated in the low intensity mode when either the side access door or the rear door is open. They may also be controlled manually from the attendant's control console.

IGN	ON
Master Power Switch (D)	ON
Dome Light, Left (HIGH) (A)	ON
Dome Light, Left (LOW) (A)	ON



## **Operating Instructions**

---

### **Dome Lights, Right**

Depending on compartment configuration, either two or three dual intensity dome lights illuminate the secondary patient area.

Ignition	ON
Master Power Switch (D)	ON
Dome Light, Right (High) (A)	ON
Dome Light, Right (Low) (A)	ON

### **Exhaust Blower**

The exhaust blower is used in conjunction with a static fresh air vent to provide ambient air changes within the patient compartment.

Ignition	ON
Master Power Switch (D)	ON
Blower Control Switch (A)	ON

### **Flood Lights, Side**

The left and right side floodlights provide external illumination for nighttime operations. They are switched separately and can be operated independently or simultaneously.

Ignition	ON
Master Power Switch (D)	ON
Right Floodlight Switch (D)	ON
Left Floodlight Switch (D)	ON

### **Floodlights, Rear (see Loading Lights, Rear)**

### **High-Beam Flasher System (optional equipment)**

When in use, this system pulsates the headlight high beam.



## **Operating Instructions**

---

Ignition	ON
Master Power Switch (D)	ON
High-Beam Flasher Switch	ON

### **Indicator Lights**

A set of indicator lights on the driver's console, gives information about the security of the doors.

Ignition	ON
----------	----

A flashing red light will appear on the driver's console if the patient compartment side access door and/or rear access is/are opened.

### **Light Bar**

The light bar is used as an additional visual alert signal to indicate the presence of an emergency vehicle. A light in the light bar switch indicates when the light bar has been activated. The rear light bar operation is identical to the front.

Ignition	ON
Master Power Switch (D)	ON
Light Bar Front (D)	ON
Light Bar Rear (D)	ON

See the light bar manufacturer's maintenance instructions for more information.

### **Loading Light(s), Rear**

The light(s) illuminate the area to the rear of the vehicle during nighttime operations.



## **Operating Instructions**

---

Ignition	ON
Master Power Switch	ON
Rear Load Light Switch (D)	ON

These lights are also activated when the vehicle is in reverse and when the rear entry doors are open.

### **Master Power Switch**

This is the first switch of the driver's switch panel and is labeled "Master". It controls all electrical equipment and devices related to patient services, including the attendant's console. Battery switch (if equipped) must be ON for the Master Switch to be active.

### **Primary/Secondary Warning Lights**

The vehicle is equipped with eight corner-mounted and two grill-mounted red flashing warning lights, a center forward-mounted clear flashing light, center rear amber flashing light, and two red intersection lights. In the primary mode, all warning lights are activated to provide a visual "clear the right-of-way" signal. In the secondary mode, the eight corner lights and rear amber light are activated to give warning that the vehicle has stopped in the right-of-way.

Ignition	ON
Master Power Switch (D)	ON
Warning Lights (Primary)	ON
Warning Lights (Secondary)	ON

### **Sequencer (optional equipment)**

The sequencer automatically "sequences" the onset of pre-set warning switches, to prevent surging caused by pre-loading the switches. All preset switches will come on upon activation of the master and sequencer switches.



## **Operating Instructions**

---

### **Siren**

For detailed operating instructions of the siren, see the operator's manual provided by the siren manufacturer. The vehicle ignition switch must be on RUN for the siren to operate.

### **Speaker, Patient Compartment (optional equipment)**

The rear speaker for stereos, if installed, is controlled from the Attendant's console.

### **Spotlights**

Hand-held spotlights are for directional nighttime illumination. Each spotlight has a momentary switch on its handle.

Ignition	ON
Master Power Switch (D)	ON
Spotlight Switch	ON

### **Strobe Low Power System**

This system allows warning lights to be used at low intensity in the case of weather conditions that create reflections back to the driver and interfere with night vision. Refer to the lighting manufacturer's operator's manual for more details.

### **Suction**

The patient aspiration system utilizes an electric suction system.

Ignition	ON
Master Power Switch (D)	ON
Electric Suction Switch (A)	ON



## Operating Instructions

---

### **Printed Circuit Board**

On vehicles equipped with the printed circuit board power distribution system, diagnostics are built in to indicate power in, power out, and the switching leg of all relays. This is accomplished with the use of three LED's at each relay, labeled for each leg of the circuit. If the "IN" LED is lit, there is power into the relay. If the "SW" LED is lit, the switch circuit is functioning properly. If the "OUT" LED is lit, the relay is functioning properly.

There are two terminal strips on the circuit board for interconnecting to the external circuitry. The terminal strips are used to help troubleshoot. Each position is labeled with the color, number, and function of the wire connected to that point.

Other than circuit breakers and relays, there are no other serviceable parts on the circuit board.

**Caution: You should check for loose connections, but be aware; you may damage the circuit board if you over-tighten the connections.**

When replacing circuit breakers, *never* exceed the maximum circuit breaker rating printed on the circuit board.



## MAINTENANCE INFORMATION



## Maintenance Information

---

### EASY CARE TIPS FOR YOUR COMMERCIAL VEHICLE PAINT FINISH

CONGRATULATIONS! Your vehicle has been painted with "Sikkens" brand paint, manufactured by AKZO NOBEL – the world's largest coating manufacturer. To maintain the factory appearance, be sure to follow these simple tips when caring for your finish.

For the first 90 days:

NO WAXING OR POLISHING  
NO ABRASIVE CLEANING METHODS OR MATERIALS

#### Proper Care

Use cool water with a mild detergent when washing.

Wash in a shaded area.

Wipe only with a cloth that is damp or wet.

Once clean water no longer beads up on the surface, use wax or polish specific for urethane paint and clear coat finishes.

#### Avoid:

Commercial washes that use high pressure, steam or caustic cleaners. Poorly maintained brush type vehicle washes. Exposure to common fluids (diesel fuel, gasoline, antifreeze or brake fluid) should be washed off immediately with mild detergent and cool water followed by a clean water rinse.

Note: Failure to follow maintenance steps can result in forfeiture of warranty coverage.



## **Maintenance Information**

---

### **Fluorescent Light Maintenance**

For preventive maintenance you should replace the fluorescent tubes when they turn black at the ends, when the light does not come to full brightness when first turned on, or if the light flickers when first turned on. Replacing the fluorescent tubes at regular intervals will prevent inverter ballast failure.

When a fluorescent light is turned on, it should be left on for a minimum of five minutes to allow the gasses in the fluorescent tubes to be activated fully. Failure to do this will shorten the life of the fluorescent tubes.

Please note that there is no such thing as a 12 volt D.C. fluorescent tube. All fluorescent tubes are for A.C. operation. The inverter ballast takes the 12 volt D.C. input and changes it to the A.C. requirement of the fluorescent tubes. Each different wattage fluorescent tube has a different A.C. requirement. Therefore inverter ballasts are lamp specific and not interchangeable for different wattages.



## Maintenance Information

---

This section should be carefully reviewed. It provides a maintenance schedule for those items that need the most frequent attention. The systems cited below are also described in the Systems Operating Section and the Maintenance and Service Information Section.

### Daily Maintenance Check

- ♦ Battery Water Level and/or Charge Condition Indicator
- ♦ Lighting Systems

*Check all bulbs daily, internal and external*

- ♦ Charging System Performance

*See Electrical Systems Operations. Continuously monitor the ammeter and voltmeter while vehicle is in operation. Refer any problems to a qualified technician.*

- ♦ Siren System Performance

*Check to see that both speakers are operating. The siren will work with one speaker but adequate sound levels cannot be obtained unless both speakers work properly. Please refer to system manual for further details.*

- ♦ Engine Fluids

*Be sure to follow the chassis manufacturer's recommendations for quality and quantity of fluids to be added and keep within the time schedule for fluid changes. Always maintain a 50/50 coolant to water ratio in the radiator*



## **Maintenance Information**

---

While the electronic unit itself does not require maintenance, the throttle actuator may require fine adjustments from time to time. The auto-throttle rpm setting must be checked every 2,000 miles; **the first 2,090 mile check is particularly important.**

Some auto-throttles are adjustable to compensate for possible changes in the alternator's output and vehicle loads. It may be necessary to increase or decrease the auto-throttle setting. If nothing is on besides the circuits which the vehicle needs to run in idle, the auto-throttle setting should never be greater than half the rpm required to reach the SAB Net Horsepower rating. Ford auto throttles should not exceed 1400 rpm. Please see manufacturer's operating manual for more information.

The **Wheeled Coach** ammeter should show charging (+) or (0) when items 1 through 9 listed on the yellow "Test Data" sticker are in operations. If it shows a negative (-) reading, have the system checked by a qualified technician.

### **Belts and Pulleys**

All belts are regarded as "normal wear" items and should be checked for signs of wear, cracking and abrasions whenever the vehicle is to be placed into service.

Whenever a belt is replaced, see that all residual rubber from the old belt is removed from the pulley grooves. Lacquer thinner or fine emery paper can be used for this purpose. If not removed, the glazing will cause the new belt to slip. Please refer to Original Equipment Manufacturer's (OBM) specifications for proper belt information.

### **Compartment and Access Door**

The basic latch system consists of two-stage rotary latch, paddle handle, and a steelstriker pin. Each of the larger doors has two rotary latches and incorporates a 5/32" plated rod(s) to couple the rotary latches.



## Maintenance Information

---

### Adjusting the Lead Door

Remove the inside panel and adjust the door as described above, using the second door edge as a guide. The second door should be aligned with the upper and lower edges of the frame.

*It is adjusted in the same manner as the other doors. When finished, replace all the panels. If the screw holes on the diamond plates do not line up with the door, replace the plates in accordance with the parts ordering procedures.*

### Vertical Adjustments

To allow for vertical adjustments, hinges are slotted on the door side. If a vertical adjustment is made, the striker pins may also need to be adjusted.



## Maintenance Information

---

### Ground Fault Circuit Interrupter (GFCI)

The shoreline AC power input line is protected by a ground fault circuit interrupter. Should an electrical short or ground fault occur on the ambulance while off-vehicle AC power is being used, the circuit breaker will trip and disconnect the AC power. To reconnect AC power, the circuit breaker must be set manually. Do not reset it until the cause of the short has been corrected.

We recommend testing the Ground Fault Circuit Interrupter (GFCI) each month. The test procedure and a chart on which to record data are provided on the unit itself. Any repair of the GFCI equipment should be undertaken only by a licensed electrician.

### Heat/Cool Air Conditioner

The patient compartment air conditioner has been especially designed for ambulance service. When properly maintained it has an indefinite life expectancy. You should therefore refer to the manufacturer's instructions in ***Equipment Manufacturer's*** section and also observe the following suggestions for keeping your unit in prime operating condition.

PAF operating info

Keep the condenser and evaporator fins clean. The dirt on the fins is dry so a soft brushing or wiping with a soft cloth should be adequate. Clean the fins on a regular basis but observe caution; the fins are made of very soft metal and can be easily damaged. See the ***Equipment Manufacturer's*** section.

Check all water hoses and fittings. Look for leaks and signs of deterioration such as weak points. Bubbles, cracks, swelling, or hoses which feel spongy or soft to the touch. Defective parts should be replaced immediately. Hoses should be replaced every 24 months as part of the coolant system tune-up which consists of:

- ◆ Draining and flushing the radiator
- ◆ Refilling the system with fresh coolant. Maintain a 50/50 mixture ratio of coolant and water
- ◆ Replacing all water hoses



## Maintenance Information

---

Check the compressor drive belt tension. Adjust the tension or replace the belt as necessary per OEM manuals.

Two vacuum shut-off valves in the engine compartment stop the flow of hot water into the rear heater core.

### Replacing a Hose

1. Drain the radiator completely, or at least below the level at which you are working. The vehicle should be on level ground.
2. Loosen the hose clamps and move them away from the connection. Rusted clamps can be cut with metal snips. Use worm-driven stainless steel replacement clamps which are sturdy and can be easily adjusted with a screwdriver.

***Caution: Silicone hoses require a special non-serrated clamp to prevent hose damage.***

3. Be gentle with hose connections. If a hose does not twist off easily, cut the hose lengthwise to remove it from the nipple.
4. Clean the nipple with a wire brush and file smooth any rough or sharp edges.
5. Place the new clamps over the new hose and slip the hose over the nipple. Check to see that the hose is free of obstructions.
6. Position and tighten the clamps around the hose, centered over the nipple.
7. Close the drain cock and refill the radiator with a 50/50 mixture of coolant and water. Start the engine and check the connection for leaks. Check again after several days of vehicle use and keep a close watch on the engine temperature during this time in case of further tightening in necessary.



## **WARRANTY AND COVERAGE**



## **Warranty Policy**

---

### **Standard Warranty Periods**

Conversion – 12 months or 12,000 miles

Modular Body – 15 Years

Limited Lifetime Electrical Warranty

Warranty cards are used to register vehicle ownership. This registration is a federal requirement. Warranty work **will not be** authorized without a signed warranty card on file at Wheeled Coach, Orlando. You can fill out the warranty registration card online by going to our website at:

<http://www.wheeledcoach.com/services/warranty/warranty-registration-form.htm>

### **Prior Authorization**

Wheeled Coach's warranty policy requires prior authorization for all warranty repairs. This allows Wheeled Coach to work with the service center to assure the customers receive the best service. Prior authorization also allows Wheeled Coach to keep costs to a minimum by avoiding unnecessary repair times and part orders. Any unauthorized repairs **may be rejected** or may **receive payment reductions** based on Wheeled Coach's review of the repair.

### **Parts**

Wheeled Coach reserves the right to replace any parts used for warranty repairs. Should this happen, Wheeled Coach will ship the replacement part via ground freight and subtract the cost from the repair bill.

Wheeled Coach requires the return of all parts valued over \$50.00. Occasionally, Wheeled Coach may require the return of items costing less than \$50.00. Wheeled Coach will issue call tags for all items returned, to avoid charges to the customer. Should Wheeled Coach require the return of an item, the customer's account won't be credited until the part is received.



## **Warranty Policy**

---

All warranty parts must be ordered through the Warranty Department. Any parts ordered as a part sale will not be credited as warranty. Vendor warranty items will be invoiced if not returned.

Wheeled Coach does not provide warranty coverage for vendor products such as sirens, flashers, speaker drivers, etc. Wheeled Coach will assist customers who are within the Wheeled Coach warranty period, in obtaining warranty coverage from the vendor. Wheeled Coach is not responsible for the administration or payment of vendor warranty. Additionally, Wheeled Coach is not responsible for repair charges not covered by the vendor.

### **Filing a Claim**

#### **Fax Submission**

The best way to submit a warranty is to fax the completed claim form to the warranty department. Faxed warranties are received quicker (and therefore process quicker), provide the sender with proof of submission, and decrease the chance of error by providing written documentation. Fax your completed claim to the warranty department at: 800-241-5177. Any fax sent elsewhere is subject to delay and possible loss.

#### **Phone Submission**

Phone submissions will be accepted, however, a completed warranty form ***must*** be submitted to process a claim.

#### **Service Centers**

Many service centers will process warranty claims with us. This will save time, so when looking for a service center, pick one that will work with Wheeled Coach.



## Warranty Policy

---

### **How to Get Paid**

#### **Original Bill Required**

Warranty payments won't be issued without an original invoice. Wheeled Coach won't accept photocopies or pay off of statements.

#### **Method of Payment**

Properly submitted warranty claims should be paid within thirty (30) days of receipt of original invoice. Warranty checks will be made out to the appropriate party as requested at time of submission.

All warranty payments will be withheld until past due parts accounts are satisfied. ***There will be no exception to this policy.*** All payment arrangements on parts accounts must be directed to Accounting or the Parts Department.



## LIMITED LIFETIME CABINET CONSTRUCTION WARRANTY

Wheeled Coach Industries, Inc., (hereinafter Wheeled Coach) ("Warrantor") warrants to the original retail purchaser, upon expiration of the attached twelve (12) months standard conversion vehicle warranty. The following parts or components of the patient compartment cabinets of the vehicle will remain free from defects in material and workmanship:

- That the wood or non-wood material used for the construction of the cabinets shall not delaminate.
- That the wooden dowels used for the construction of the cabinets shall not allow the cabinet sections to separate.

This Cabinet Construction Warranty commences upon the expiration of the original Wheeled Coach twelve (12) months standard conversion vehicle warranty and continues for the lifetime of the vehicle for the original owner on the original chassis. For the purpose of the Limited Lifetime Cabinet Construction Warranty, a lifetime is defined by Wheeled Coach as: 10 years from the expiration of the original Wheeled Coach twelve (12) months standard conversion vehicle warranty for only the only original retail purchaser/owner. This warranty is valid only in the United States and Canada, and any vehicle sold for use outside of the United States and Canada is sold "AS IS".

If the owner discovers within this period a defect in the workmanship or material of any of the conditions listed above, it must promptly notify Wheeled Coach in writing. In no event shall such notification be received by Wheeled Coach later than one month after the date this Cabinet Construction Warranty expires. Within a reasonable time after such notification, Wheeled Coach will correct any defect in workmanship or material by repairing or replacing same, at Wheeled Coach's sole option. Such repair, including both parts and labor is at Wheeled Coach's expense. All warranty service is subject to Wheeled Coach's prior examination and written approval and will be performed by Wheeled Coach or at service centers designated by Wheeled Coach. All transportation to and from the designated service center will be at the owner's expense and is not included as a cost of repair covered by this warranty. These remedies are the owner's exclusive remedies for breach of warranty.

### Warranty Exclusions and Limitations

This Limited Lifetime Cabinet Construction Warranty DOES NOT COVER:

1. Any product, components or parts not manufactured by Wheeled Coach or covered by the warranty of another manufacturer, whether installed by Wheeled Coach or parties other than Wheeled Coach, including but not limited to hinges, cabinet hardware, Lexan sliding or hinged door sections, frames around said door sections, cabinet or door hold open or closed devices, handles and similar equipment or items,
2. Damage caused by use of the vehicle for purposes other than those for which it was designed
3. Damage caused by accident or the negligence of the owner or any third party or by disasters such as fire, flood, wind and lightning,
4. Routine maintenance and maintenance items
5. Damage caused by the owner's failure to provide normal maintenance
6. Damage caused by unauthorized or improper installation of attachments, repairs or modifications
7. Damage during shipment
8. Any other abuse or misuse by the owner.
9. Exported and/or vehicles shipped to an international destination are excluded from this policy.

**Limitation on Damages:** Wheeled Coach shall not be liable for any special, incidental, or consequential damages based upon breach of warranty, breach of contract, negligence, strict tort, or any other legal theory. Such damages include, but are not limited to, lost wages, lost vehicle rental expenses, loss of savings, profits, or revenue, loss of use of the vehicle or any associated equipment, cost of capital, cost of any substitute vehicle, equipment, facilities or services, downtime, the claims of third parties including customers, and injury to property. This limitation does not apply to claims for personal injury. Some states do not allow limits on warranties, or on remedies for breach in certain transactions. In such states, the limits in this paragraph and in the disclaimer of warranty may not apply.

**Warranty Registration:** This warranty is conditioned upon receipt by Wheeled Coach of a completed and signed warranty registration card within two weeks of delivery. It is the obligation of the purchaser to sign the customer registration card and return it to Wheeled Coach within the two weeks following delivery. The warranty registration card must be on file for any warranty claim to be considered. This warranty is not transferrable.

**Warranty Claim Procedure:** If the owner discovers, within the limited lifetime warranty period, a defect in material or workmanship, the owner must contact Wheeled Coach in writing or call Wheeled Coach's customer service department toll free at 855-661-9232 Monday through Friday from 8:00 to 5:00 pm (Eastern Time) prior to any warranty work being performed. Wheeled Coach will provide the purchaser with a warranty claim form and further instructions on how to proceed with such warranty claim. All warranty correspondence must be sent to Wheeled Coach, PO Box 677339, Orlando, Florida 32867-7339. In the case of products, components or parts covered by the warranty of another manufacturer, Wheeled Coach will pass on to the owner the existing manufacturer's warranty to the extent reasonably possible. However, purchasers are responsible for checking the original component manufacturer's warranty regarding its coverage. In no event shall such notification be received by Wheeled Coach later than one month following the expiration of the warranty period. All warranty work is subject to Wheeled Coach's examination and approval. If authorized by Wheeled Coach, the obligation of Wheeled Coach is limited to correcting any defect in material or workmanship, with either new or used replacement parts, at Wheeled Coach's option, to the original purchaser. If approved, such repair, including both parts and labor, is at Wheeled Coach's expense.

**Third Party Representations:** Wheeled Coach does not authorize any person to create for Wheeled Coach any other obligations or liability in connection with its custom module, and Wheeled Coach is not responsible for any representation, promise or warranty made by a dealer, component or vehicle manufacturer, or other person beyond what is expressly stated in this limited warranty.

**Warrantor's Rights:** Warrantor reserves the unrestricted right to change the parts and designs of its products from time to time without notice and with no obligation to make corresponding changes in the product previously manufactured.

**Disclaimer of Warranty:** The foregoing warranties are in lieu of all other warranties, express or implied. The implied warranties of merchantability and fitness for a particular purpose are expressly disclaimed.

**Warranty Termination:** The warranty shall be void, and Wheeled Coach shall be relieved from any and all obligations hereunder, if:

1. The motor vehicle is sold or any of the components are removed.
2. The chassis is remounted, unless the remount work is performed by Wheeled Coach or by a service facility approved, in advance and in writing, by Wheeled Coach.
3. Owner fails to comply with the warranty registration requirements described above.
4. Policy expiration date

**Time Limit for requesting Arbitration:** Any action for breach of warranty must be commenced within three months following the expiration of the warranty period.

**Legal Remedies:** Any claim or controversy arising out of or relating to this limited warranty, or breach thereof, shall be settled by arbitration administered by the American Arbitration Association in Orange County, Florida in accordance with the Commercial Arbitration Rules of the American Arbitration Association. The determination of the arbitrator(s) shall be in writing and shall include an explanation of the basis for the determination. The determination of the arbitrator(s) shall be final and binding and judgment upon such determination may be entered in any court having jurisdiction.

**No Other Warranty:** Unless modified in a written statement signed by both parties, this agreement is understood to be the complete and exclusive agreement between the parties, superseding all prior agreements, oral or written, and all other communications between the parties (including without limitation any terms and conditions contained in any purchase order or sales invoice issued pursuant to the sale of the vehicle) relating to the subject matter of this agreement. No employee of Marque or any other party is authorized to make any warranty in addition to those made in this agreement.



## Limited Electrical Warranty

Wheeled Coach Industries, Inc., (hereinafter Wheeled Coach) ("Warrantor") warrants to the original retail purchaser for a period of twelve (12) months or twelve thousand (12,000) miles from the date of delivery of the completed new custom module to the end user, regardless of subsequent ownership. This product shall be free of substantial defects in materials and workmanship, which are attributable to Wheeled Coach and which arise during the course of normal use and service. Wheeled Coach will correct any defect in covered parts or workmanship, with either new or used replacement parts, at Wheeled Coach's option. Covered parts are limited to custom module electrical systems and components such as electrical harness, harness installation, wires (but only to the extent that wires are broken, chafed, or pinched), electrical connections, terminal blocks, junction posts, and related components. Printed circuit boards are covered for a limited lifetime. A lifetime is defined by Wheeled Coach Limited Electrical Warranty as; 10 years from the expiration of the original Wheeled Coach standard conversion warranty which is for an unlimited mileage for a period of twelve (12) months from the date of delivery for only original retail purchaser/owner. All warranty service is subject to Wheeled Coach prior examination and (written) approval and will be performed by Wheeled Coach or at service centers designated by Wheeled Coach. This warranty is valid only in the United States and Canada, and any vehicle sold for use outside of the United States and Canada is sold "AS IS."

### Component Manufacturer Warranty

Original component manufacturers may provide their own warranties. Purchasers should check the original component manufacturer's warranty regarding its coverage. Wheeled Coach may assist the original retailer in submitting claim to the original component manufacturer within original manufacturer term period mentioned below.

- ❖ **Weldon (Vehicle Multiplex system)** (hereinafter VMUX) correctly is warranted against mechanical, electrical and physical defects for the period defined in the table below per module. The period is defined as the date of manufacture from Weldon; each module carries its own date of manufacture. Without (VMUX Certified Supplier) VCS, Weldon offers a 30-day parts-only warranty on all V-MUX hardware. Written notice of a defect must be received by Weldon within four (4) years after initial product shipment; unless an extended warranty is purchased. Note: Extended warranties are offered by purchasing tamper proof stickers from Weldon with the years of the warranty on them. Extended warranties extend the base warranty, a node with a 7 year sticker has an additional 3 years added to the 4 year standard. Damaged or removal of the label will result in no extended warranty. Weldon will correct by repair or replacement, at its option, equipment or parts which fail because of mechanical, electrical or physical defects, provided that the goods have been properly handled and stored prior to installation and properly installed and properly operated after installation. End users requiring warranty repairs on the V-MUX® system should contact the original equipment manufacturer's customer service or other appropriate department for service/warranty repairs. Weldon Repair Parts Labor Period is defined as the labor provided by Weldon's service department to repair/replace the returned device. This does not reflect the labor to remove/replace any such devices from the vehicle, nor indicate Weldon accepts any responsibility for such removal, replacement or troubleshooting of said devices. Please read below warranty exclusions below for items not covered by Weldon. Customer is responsible for verifying manufacturer warranties.
- ❖ **Weldon Exclusions and Limitations** – Weldon excludes and/or limits:
  1. Cost to remove defective part, installation of repaired product, labor or consequential damages of any kind, and the exclusive remedy being to require such parts to be furnished.
  2. Weldon's liability under no circumstances will exceed the contract price of goods claimed to be defective
  3. Any returns under this guarantee are to be on a transportation charges prepaid basis.

### WARRANTY PERFORMANCE

Owner's exclusive remedy under this warranty shall be repair or replacement of defective parts, free of charge to the Owner. Owner shall bear all expenses arising out of or relating to transporting the product to the appropriate Warranty Service Location.

### Warranty Exclusions and Limitations

This Limited Ambulance Warranty **DOES NOT COVER**:

1. Damage to products that have been improperly installed, altered or modified by any party other than Warrantor.
2. Chassis electrical system and related components installed by the chassis supplier.
3. Electrical items which are not manufactured by Wheeled Coach, but which are installed by Wheeled Coach personnel. These items include, but are not limited to, the following: battery chargers; siren; inverters; generators; light bars and similar OEM equipment. These excluded items are typically covered by separate warranties supplied by the OEM manufacturer of the component. When necessary Wheeled Coach personnel may assist the purchaser in pursuing warranty assistance with the OEM manufacturer of these excluded-components.
4. The cab and chassis (owner is responsible for verify warranty provided by Cab & Chassis manufacturer.) Any collateral damage from cab and chassis related electrical failure.
5. Replacement of routine maintenance items, such as voltage regulator, flashers, bulbs, batteries, switches, etc.
6. Non Electrical Components and Items not specified hereinabove as covered items.
7. Damage caused by misuse, neglect, negligence or accident. Usage of this product in a manner, which is inconsistent with design intentions or inconsistent with owner's manual directions, will invalidate this warranty in regard to damage caused by or relating to such inconsistent usage.
8. Periodic tightening and cleaning of connection terminals as required by customary routine maintenance.
9. Damage caused by owner's failure to provide normal maintenance. Failure to follow operating instructions.
10. Damage caused by floods, fire, wind, hail, lightning, and/or any natural disaster.
11. Damage caused by replacement or modification of original parts and/or components this includes but is not limited to improper or unauthorized installation of repairs, and/or attachments. Installation of any "aftermarket" devices including ancillary equipment used in emergency service calls or the modification of any existing system or component originally installed by Wheeled Coach without Wheeled Coach's express written approval and any problems resulting from such installation or modification.
12. Wheeled Coach shall not be liable for incidental, consequential, direct, indirect, or other damages (such as, but not limited to, lost wages, inconvenience, lodging, travel costs (expense of transporting the product to an authorized service center for service, and all expense arising from or relating to such transportation), or lost vehicle rental expenses) that result from breach of written warranty or any implied warranty.
13. Damage to a Wheeled Coach vehicle that is leased or rented to a second party

**Limitation on Damages:** Wheeled Coach shall not be liable for any special, incidental, or consequential damages based upon breach of warranty, breach of contract, negligence, strict tort, or any other legal theory. Such damages include, but are not limited to, lost wages, lost vehicle rental expenses, loss of savings, profits, or revenue, loss of use of the vehicle or any associated equipment, cost of capital, cost of any substitute vehicle, equipment, facilities or services, downtime, the claims of third parties including customers, and injury to property. This limitation does not apply to claims for personal injury. Some states do not allow limits on warranties, or on remedies for breach in certain transactions. In such states, the limits in this paragraph and in the disclaimer of warranty may not apply.

**Warranty Registration:** This warranty is conditioned upon receipt by Wheeled Coach of a completed and signed warranty registration card within two weeks of delivery. It is the obligation of the purchaser to sign the customer registration card and return it to Wheeled Coach within the two weeks following delivery. The warranty registration card must be on file for any warranty claim to be considered. This warranty is not transferrable.

**Warranty Claim Procedure:** If the owner discovers within this period a defect in material or workmanship it must promptly notify Wheeled Coach in writing. In no event shall such notification be received by Wheeled Coach later than 30 days after said warranty expires. All warranty work is subject to Wheeled Coach's examination and approval. If authorized by Wheeled Coach, the obligation of Wheeled Coach is limited to repairing or replacing, per the warranty schedule, to the original purchaser. The Purchaser must contact Wheeled Coach in writing or call Wheeled Coach's customer service department toll free at 855-661-9232 Monday through Friday from 8:00 to 5:00 pm (Eastern Time) prior to any warranty work being performed. Wheeled Coach will provide the purchaser with a warranty claim form and further instructions on how to proceed with such warranty claim. All warranty correspondence must be sent to Wheeled Coach Industries, PO Box 677339, Orlando, Florida 32867-7339. In the case of products, components or parts covered by the warranty of another manufacturer, Wheeled Coach will pass on to the owner the existing manufacturer's warranty to the extent reasonably possible. All claimed defective parts are to be returned by the original purchaser, transportation expenses for returned parts pre-paid to Wheeled Coach within 30 days of receipt of the replacement part or parts sent by Wheeled Coach. If the claimed defective part(s) are not received by Wheeled Coach within the 30 day grace period, or if the claimed defect cannot be authenticated, the original purchaser will be charged the full value of the replacement parts plus shipping and handling.

**Third Party Representations:** Wheeled Coach does not authorize any person to create for Wheeled Coach any other obligations or liability in connection with its custom module, and Wheeled Coach is not responsible for any representation, promise or warranty made by a dealer, component or vehicle manufacturer, or other person beyond what is expressly stated in this limited warranty.

**Warrantor's Rights:** Warrantor reserves the unrestricted right to change the parts and designs of its products from time to time without notice and with no obligation to make corresponding changes in the product previously manufactured.

**Warranty Termination:** The warranty shall be void, and Wheeled Coach shall be relieved from any and all obligations hereunder, if:

1. The motor vehicle is sold or any of the components are removed.
2. The chassis is remounted, unless the remount work is performed by Wheeled Coach or by a service facility approved, in advance and in writing, by Wheeled Coach.
3. Owner fails to comply with the warranty registration requirements described above.
4. Policy expiration date

**Time Limit for requesting Arbitration:** Any action for breach of warranty must be commenced within three months following the expiration of the warranty period.

**Legal Remedies:** Any claim or controversy arising out of or relating to this limited warranty, or breach thereof, shall be settled by arbitration administered by the American Arbitration Association in Orange County, Florida in accordance with the Commercial Arbitration Rules of the American Arbitration Association. The determination of the arbitrator(s) shall be in writing and shall include an explanation of the basis for the determination. The determination of the arbitrator(s) shall be final and binding and judgment upon such determination may be entered in any court having jurisdiction.

**No Other Warranty:** Unless modified in a written statement signed by both parties, this agreement is understood to be the complete and exclusive agreement between the parties, superseding all prior agreements, oral or written, and all other communications between the parties (including without limitation any terms and conditions contained in any purchase order or sales invoice issued pursuant to the sale of the vehicle) relating to the subject matter of this agreement. No employee of Wheeled Coach or any other party is authorized to make any warranty in addition to those made in this agreement.





## Modular Structure Limited Lifetime Warranty

Wheeled Coach Industries, Inc., (hereinafter Wheeled Coach) ("Warrantor") warrants to the original retail purchaser only, that the module structure that is the subject of this sale is structurally sound and free from all structural defects in material and workmanship and further warrants the module structure will remain free of structural damage due to rusting caused by electrolysis. The custom module structure limited warranty is in effect for the lifetime of a new vehicle. For the purpose of the lifetime custom module limited warranty, a lifetime is defined by Wheeled Coach as; 15 years from the date of original retail owner's purchase/in-service date from Wheeled Coach or the period of time the ambulance is in continuous front line service with the original retail purchaser. In the event of a module remount this custom module structural warranty shall remain in effect provided the remount work is completed within the defined lifetime period, and remount work is completed by Wheeled Coach or a facility Wheeled Coach authorizes. This limited warranty covers repairs or replacement of any part of your new Wheeled Coach custom module structure (hereinafter Covered Parts) in which a defect in materials or workmanship appears during normal use, maintenance or service within the limited warranty period, subject to the limitations and exclusions. The repair or replacement does not extend the life of the limited warranty except where state or provincial law otherwise provides for an extension during the time that the Covered Part is being repaired or replaced under this limited warranty. Covered Parts are limited to custom module systems and components such as custom module body, exterior doors, and interior aluminum cabinets. Wheeled Coach will not replace the ambulance or repurchase the ambulance. All warranty service is subject to Wheeled Coach's prior examination and (written) approval and will be performed by Wheeled Coach or at service centers designated by Wheeled Coach. This warranty is valid only in the United States and Canada, and any vehicle sold for use outside of the United States and Canada is sold "AS IS."

### **WARRANTY PERFORMANCE**

Owner's exclusive remedy under this warranty shall be repair or replacement of defective parts, free of charge to the Owner. Owner shall bear all expenses arising out of or relating to transporting the product to the appropriate Warranty Service Location.

### **Warranty Exclusions and Limitations**

#### This Limited Ambulance Warranty DOES NOT COVER:

1. All paint, striping, hardware, hinges, latches, moldings, windows, appointments and accessories.
2. Any and all structural defects resulting from any form of accident or operator negligence.
3. Material or parts manufactured by other entities, including but not limited to sealants, hardware, moldings, windows, fenderettes and other appointments and accessories.
4. Damage caused by misuse, neglect, negligence or accident. Usage of this product in a manner, which is inconsistent with design intentions or inconsistent with owner's manual directions, will invalidate this warranty in regard to damage caused by or relating to such inconsistent usage.
5. Damage caused by owner's failure to provide normal maintenance. Failure to follow operating instructions.
6. Damage caused by floods, fire, wind, hail, lightning, any natural disaster, overloading, chemical fallout, industrial contamination, and/or vandalism.
7. Damage caused by replacement or modification of original parts and/or components this includes but is not limited to improper or unauthorized installation of attachments, repairs, modifications or alterations.
8. Damage caused by replacement of original parts or components with unauthorized substitutes.
9. Damage during shipment.
10. Wheeled Coach shall not be liable for incidental, consequential, direct, indirect, or other damages (such as, but not limited to, lost wages, inconvenience, lodging, travel costs, or lost vehicle rental expenses) that result from breach of written warranty or any implied warranty.
11. Damage to a Wheeled Coach vehicle that is leased or rented to a second party.
12. Damage to products that have been improperly installed, altered or modified by any party other than Warrantor.
13. General tightening, lubrication of latches, catches, and electrical components.
14. Interior cabinets.
15. Fiberglass components.
16. Installation of any "aftermarket" devices including ancillary equipment used in emergency service calls or the modification of any existing system or component originally installed by Wheeled Coach without Wheeled Coach's prior express written approval and any problems resulting from such installation or modification.
17. Custom module if altered by welding, cutting or splicing, or improper drilling without Wheeled Coach's prior written approval.
18. Custom module components damaged as a result of corrosion, including, but not limited to, acid rain, blisters, deicing chemicals, road salt additives and/or acidic exposure.

**Limitation on Damages:** Wheeled Coach shall not be liable for any special, incidental, or consequential damages based upon breach of warranty, breach of contract, negligence, strict tort, or any other legal theory. Such damages include, but are not limited to, lost wages, lost vehicle rental expenses, loss of savings, profits, or revenue, loss of use of the vehicle or any associated equipment, cost of capital, cost of any substitute vehicle, equipment, facilities or services, downtime, the claims of third parties including customers, and injury to property. This limitation does not apply to claims for personal injury. Some States do not allow limits on warranties, or on remedies for breach in certain transactions. In such States, the limits in this paragraph and in the disclaimer of warranty may not apply.

**Warranty Registration:** This warranty is conditioned upon receipt by Wheeled Coach of a completed and signed warranty registration card within two weeks of delivery. It is the obligation of the purchaser to sign the customer registration card and return it to Wheeled Coach within the two weeks following delivery. The warranty registration card must be on file for any warranty claim to be considered. This warranty is not transferable.

**Warranty Claim Procedure:** If the owner discovers within this period a defect in material or workmanship, the owner must promptly notify Wheeled Coach in writing. In no event shall such notification be received by Wheeled Coach later than 3 months from the expiration of the Modular Structure Limited Lifetime Warranty. All warranty work is subject to Wheeled Coach's examination and approval. If authorized by Wheeled Coach, the obligation of Wheeled Coach is limited to repairing or replacing, per the warranty term period, to the original purchaser. The Purchaser must contact Wheeled Coach in writing or call Wheeled Coach's customer service department toll free at 855-661-9232 Monday through Friday from 8:00 to 5:00 pm (Eastern Time) prior to any warranty work being performed. Wheeled Coach will provide the purchaser with a warranty claim form and further instructions on how to proceed with such warranty claim. All warranty correspondence must be sent to Wheeled Coach, PO Box 677339, Orlando, Florida 32867-7339. In the case of products, components or parts covered by the warranty of another manufacturer, Wheeled Coach will pass on to the owner the existing manufacturer's warranty to the extent reasonably possible. All claimed defective parts are to be returned by the original purchaser, transportation expenses for returned parts pre-paid to Wheeled Coach within 30 days of receipt of the replacement part or parts sent by Wheeled Coach. If the claimed defective part(s) are not received by Wheeled Coach within the 30 day grace period, or if the claimed defect cannot be authenticated, the original purchaser will be charged the full value of the replacement parts plus shipping and handling.

**Third Party Representations:** Wheeled Coach does not authorize any person to create for Wheeled Coach any other obligations or liability in connection with its custom module, and Wheeled Coach is not responsible for any representation, promise or warranty made by a dealer, component or vehicle manufacturer, or other person beyond what is expressly stated in this limited warranty.

**Warrantor's Rights:** Warrantor reserves the unrestricted right to change the parts and designs of its products from time to time without notice and with no obligation to make corresponding changes in the product previously manufactured.

**Disclaimer of Warranty:** The foregoing warranties are in lieu of all other warranties, express or implied. The implied warranties of merchantability and fitness for a particular purpose are expressly disclaimed.

**Warranty Termination:** The warranty shall be void, and Wheeled Coach shall be relieved from any and all obligations hereunder, if:

1. The motor vehicle is sold or any of the components are removed.
2. The chassis is remounted, unless the remount work is performed by Wheeled Coach or by a service facility approved, in advance and in writing, by Wheeled Coach.
3. Owner fails to comply with the warranty registration requirements described above.
4. Policy expiration date

**Time Limit for requesting Arbitration:** Any action for breach of warranty must be commenced within three months following the expiration of the warranty period.

**Legal Remedies:** Any claim or controversy arising out of or relating to this limited warranty, or breach thereof, shall be settled by arbitration administered by the American Arbitration Association in Orange County, Florida in accordance with the Commercial Arbitration Rules of the American Arbitration Association. The determination of the arbitrator(s) shall be in writing and shall include an explanation of the basis for the determination. The determination of the arbitrator(s) shall be final and binding and judgment upon such determination may be entered in any court having jurisdiction.

**No Other Warranty:** Unless modified in a written statement signed by both parties, this agreement is understood to be the complete and exclusive agreement between the parties, superseding all prior agreements, oral or written, and all other communications between the parties (including without limitation any terms and conditions contained in any purchase order or sales invoice issued pursuant to the sale of the vehicle) relating to the subject matter of this agreement. No employee of Wheeled Coach or any other party is authorized to make any warranty in addition to those made in this agreement.





**Wheeled  
Coach®** OUT IN FRONT

Modular Warranty Certificate  
Body Paint Warranty

Wheeled Coach Industries, Inc., (hereinafter Wheeled Coach) ("Warrantor") warrants to the original retail purchaser, under normal use and service, each new Wheeled Coach modular body paint job is free of all material and workmanship defects for a prorated period of five (5) years from the date of delivery. All warranty service is subject to Wheeled Coach's prior examination and (written) approval. All warranty service will be performed by Wheeled Coach or at service centers designated by Wheeled Coach. This warranty is valid only in the United States and Canada, and any vehicle sold for use outside of the United States and Canada is sold "AS IS." The paint applied by Wheeled Coach is limited to the original user and limited to the exterior painted surface of the module according to this warranty schedule:

0-36 months	100%
37-48 months	50%
49-60 months	25%

The warranty provided herein shall cover and extend to the following properties of the paint system according to the warranty schedule:

- Loss of adhesion of the paint system resulting in rust
- Cracking of paint system
- Fading or loss of gloss

Repair or replacement of covered parts by a Wheeled Coach Authorized Service Center is the exclusive remedy under this limited warranty. Wheeled Coach will not replace the ambulance or repurchase the ambulance from you. The repair of a vehicle does not extend the life of the warranty except where state or provincial law otherwise provides for an extension during the time covered repair is under warranty.

**Modular Body Paint Warranty Exclusions and Limitations**

The following conditions are specifically excluded from the Wheeled Coach Modular Body Paint Warranty:

- Exterior finish on any fiberglass components, including but not limited to hardware moldings, windows, mirrors, lights and other appointments and accessories
- Hazing, chalking, or loss of gloss caused by improper care, abrasive polishes, cleaning agents, heavy duty pressure washing, or aggressive mechanical wash systems;
- Paint deterioration caused by abuse, accidents, acid rain, hail, blisters, deicing chemicals, road salt additives, film degradation due to rust or corrosion originating from substrate, chemical fallout or acts of nature;
- Damage resulting from crevice corrosion
- Accidents, scratches, chips, bruises, and gloss reduction due to normal vehicle use and maintenance
- Re-application of custom finishes, custom decals, hand lettering, exotic finishes or any other finish other than standard refinish procedure;
- Scratches, chips, UV Paint fade, or gloss reduction due to normal wear and tear
- All paint hardware, hinges, latches, windows, appointments and accessories
- Gold leaf, decals, graphics, stickers, name plates, or striping except that which is affected by repair
- Paint failures resulting from the improper application of striping, tapes, decals or custom painted designs I numerals
- Repairs done to previously refinished areas unless stripped to bare metal and appropriate substrate
- Any product finishes, component finishes or finishes of parts not manufactured by Wheeled Coach, including the chassis
- Defects resulting from normal and customary wear as a result of operating the vehicle
- Defects caused by installation of any devices and/or modification of any existing system or component originally installed by Wheeled Coach without Wheeled Coach's prior express written approval and any problems resulting from such installation or modification
- Repairs or alterations done over previously refinished areas, unless stripped to bare metal or appropriate substrate
- Paint if altered by welding, cutting or splicing, or improper drilling without Wheeled Coach's prior written approval.
- Damage that occurs during shipment
- Damage caused by abrasion of external foreign objects
- Claims presented without proper Warranty documentation or pictures
- Wheeled Coach shall not be liable for incidental, consequential, direct, indirect, or other damages (such as, but not limited to, lost wages, lodging, travel costs, or lost vehicle rental expenses) that result from breach of written warranty or any implied warranty.

**Limitation on Damages:** Wheeled Coach shall not be liable for any special, incidental, or consequential damages based upon breach of warranty, breach of contract, negligence, strict tort, or any other legal theory. Such damages include, but are not limited to, lost wages, lost vehicle rental expenses, loss of savings, profits, or revenue, loss of use of the vehicle or any associated equipment, cost of capital, cost of any substitute vehicle, equipment, facilities or services, downtime, the claims of third parties including customers, and injury to property. This limitation does not apply to claims for personal injury. Some states do not allow limits on warranties, or on remedies for breach in certain transactions. In such states, the limits in this paragraph and in the disclaimer of warranty may not apply.

**Warranty Registration:** This warranty is conditioned upon receipt by Wheeled Coach of a completed and signed warranty registration card within two weeks of delivery. It is the obligation of the purchaser to sign the customer registration card and return it to Wheeled Coach within the two weeks following delivery. The warranty registration card must be on file for any warranty claim to be considered. This warranty is not transferrable.

**Warranty Claim Procedure:** The original purchaser shall notify Wheeled Coach within thirty (30) days of the first signs of the claimed defect or perforation. All warranty work is subject to Wheeled Coach's examination and approval. If authorized by Wheeled Coach, the obligation of Wheeled Coach is limited to repairing or replacing, per the warranty schedule, to the original purchaser. The Purchaser must contact Wheeled Coach in writing or call Wheeled Coach's customer service department toll free at 855-661-9232 Monday through Friday from 8:00 to 5:00 pm (Eastern Time) prior to any warranty work being performed. Wheeled Coach will provide the purchaser with a warranty claim form and further instructions on how to proceed with such warranty claim. All warranty correspondence must be sent to Wheeled Coach, PO Box 677339, Orlando, Florida 32867-7339. In the case of products, components or parts covered by the warranty of another manufacturer, Wheeled Coach will pass on to the owner the existing manufacturer's warranty to the extent reasonably possible.

**Third Party Representations:** Wheeled Coach does not authorize any person to create for Wheeled Coach any other obligations or liability in connection with its custom module, and Wheeled Coach is not responsible for any representation, promise or warranty made by a dealer, component or vehicle manufacturer, or other person beyond what is expressly stated in this limited warranty.

**Warrantor's Rights:** Warrantor reserves the unrestricted right to change the parts and designs of its products from time to time without notice and with no obligation to make corresponding changes in the product previously manufactured.

**Disclaimer of Warranty:** The foregoing warranties are in lieu of all other warranties, express or implied. The implied warranties of merchantability and fitness for a particular purpose are expressly disclaimed.

**Warranty Termination:** The warranty shall be void, and Wheeled Coach shall be relieved from any and all obligations hereunder, if:

1. The motor vehicle is sold or any of the components are removed.
2. The chassis is remounted, unless the remount work is performed by Wheeled Coach or by a service facility approved, in advance and in writing, by Wheeled Coach.
3. Owner fails to comply with the warranty registration requirements described above.
4. Policy expiration date

**Time Limit for requesting Arbitration:** Any action for breach of warranty must be commenced within three months following the expiration of the warranty period.

**Legal Remedies:** Any claim or controversy arising out of or relating to this limited warranty, or breach thereof, shall be settled by arbitration administered by the American Arbitration Association in Orange County, Florida in accordance with the Commercial Arbitration Rules of the American Arbitration Association. The determination of the arbitrator(s) shall be in writing and shall include an explanation of the basis for the determination. The determination of the arbitrator(s) shall be final and binding and judgment upon such determination may be entered in any court having jurisdiction.

**No Other Warranty:** Unless modified in a written statement signed by both parties, this agreement is understood to be the complete and exclusive agreement between the parties, superseding all prior agreements, oral or written, and all other communications between the parties (including without limitation any terms and conditions contained in any purchase order or sales invoice issued pursuant to the sale of the vehicle) relating to the subject matter of this agreement. No employee of Marque or any other party is authorized to make any warranty in addition to those made in this agreement.



# U.S. AMBULANCE CORPORATION



## Standard Conversion Warranty

Wheeled Coach Industries, Inc., (hereinafter Wheeled Coach) ("Warrantor") warrants to the original retail purchaser for an unlimited mileage for a period of twelve (12) months from the date of delivery. This product shall be free of substantial defects in materials and workmanship, which are attributable to Wheeled Coach and which arise during the course of normal use and service. Wheeled Coach will correct any defect in material or workmanship, with either new or used replacement parts, at Wheeled Coach's option. Such repair including both parts and labor, is at Wheeled Coach's expense. All warranty service is subject to Wheeled Coach's prior examination and (written) approval and will be performed by Wheeled Coach or at service centers designated by Wheeled Coach. This warranty is valid only in the United States and Canada, and any vehicle sold for use outside of the United States and Canada is sold "AS IS."

### **WARRANTY PERFORMANCE**

Owner's exclusive remedy under this warranty shall be repair or replacement of defective parts, free of charge to the Owner. Owner shall bear all expenses arising out of or relating to transporting the product to the appropriate Warranty Service Location.

### **Warranty Exclusions and Limitations**

This Limited Ambulance Warranty **DOES NOT COVER:**

1. Damage to the soft trim and appearance items if such damage is due to normal use, wear and tear, or exposure to elements.
2. Accessories, components, and/or parts not manufactured by Warrantor, which items include (but are not limited to): the chassis and its component parts, tires, tire balancing, aftermarket suspension components, wheel alignment, inverters, sirens, light bulbs, light bars, battery chargers, generator, air conditioners, radios, power converters, batteries, other electronic, sealants, hardware, moldings, windows, fender rings, and/or attachments and accessories. The manufacturers of these products may provide warranties covering the performance of their particular products.
3. Normal maintenance such as lubrication, batteries, tires, filter and oil replacement, belts and hoses, brake lining and adjustment, and vehicle alignments. Normal wear parts such as electrical accessories, voltage regulator, flashers, switches, etc.
4. Scratches or dents on the body, windows, and window shields not reported within 7 days of delivery.
5. Vinyl painted graphics in any form not originally installed by Wheeled Coach.
6. Damage to products that have been improperly installed, altered or modified by any party other than Warrantor.
7. Damage caused by misuse, neglect, negligence or accident. Usage of this product in a manner, which is inconsistent with design intentions or inconsistent with owner's manual directions, will invalidate this warranty in regard to damage caused by or relating to such inconsistent usage.
8. Damage caused by owner's failure to provide normal maintenance. Failure to follow operating instructions.
9. Damage caused by floods, fire, wind, hail, lightning, and/or any natural disaster.
10. Damage caused by replacement or modification of original parts and/or components. This includes but is not limited to improper or unauthorized installation of repairs, and/or attachments.
11. Wheeled Coach shall not be liable for incidental, consequential, direct, indirect, or other damages (such as, but not limited to, lost wages, inconvenience, lodging, travel costs, or lost vehicle rental expenses) that result from breach of written warranty or any implied warranty.
12. Damage to a Wheeled Coach vehicle that is leased or rented to a second party
13. Custom module if altered by welding, cutting or splicing, or improper drilling without Wheeled Coach's prior written approval
14. Custom module systems and components damaged as a result of corrosion, including, but not limited to, salt and/or acidic exposure

**Limitation on Damages:** Wheeled Coach shall not be liable for any special, incidental, or consequential damages based upon breach of warranty, breach of contract, negligence, strict tort, or any other legal theory. Such damages include, but are not limited to, lost wages, lost vehicle rental expenses, loss of savings, profits, or revenue, loss of use of the vehicle or any associated equipment, cost of capital, cost of any substitute vehicle, equipment, facilities or services, downtime, the claims of third parties including customers, and injury to property. This limitation does not apply to claims for personal injury. Some states do not allow limits on warranties, or on remedies for breach in certain transactions. In such states, the limits in this paragraph and in the disclaimer of warranty may not apply.

**Warranty Registration:** This warranty is conditioned upon receipt by Wheeled Coach of a completed and signed warranty registration card within two weeks of delivery. It is the obligation of the purchaser to sign the customer registration card and return it to Wheeled Coach within the two weeks following delivery. The warranty registration card must be on file for any warranty claim to be considered. This warranty is not transferrable.

**Warranty Claim Procedure:** If the owner discovers within this period a defect in material or workmanship, the owner must promptly notify Wheeled Coach in writing. In no event shall such notification be received by Wheeled Coach later than 13 months from the date of delivery to the original purchaser. All warranty work is subject to Wheeled Coach's examination and approval. If authorized by Wheeled Coach, the obligation of Wheeled Coach is limited to repairing or replacing, per the warranty schedule, to the original purchaser. The Purchaser must contact Wheeled Coach in writing or call Wheeled Coach's customer service department toll free at 855-661-9232 Monday through Friday from 8:00 to 5:00 pm (Eastern Time) prior to any warranty work being performed. Wheeled Coach will provide the purchaser with a warranty claim form and further instructions on how to proceed with such warranty claim. All warranty correspondence must be sent to Wheeled Coach, PO Box 677339, Orlando, Florida 32867-7339. In the case of products, components or parts covered by the warranty of another manufacturer, Wheeled Coach will pass on to the owner the existing manufacturer's warranty to the extent reasonably possible. All claimed defective parts are to be returned by the original purchaser, transportation expenses for returned parts pre-paid to Wheeled Coach within 30 days of receipt of the replacement part or parts sent by Wheeled Coach. If the claimed defective part(s) are not received by Wheeled Coach within the 30 day grace period, or if the claimed defect cannot be authenticated, the original purchaser will be charged the full value of the replacement parts plus shipping and handling.

**Third Party Representations:** Wheeled Coach does not authorize any person to create for Wheeled Coach any other obligations or liability in connection with its custom module, and Wheeled Coach is not responsible for any representation, promise or warranty made by a dealer, component or vehicle manufacturer, or other person beyond what is expressly stated in this limited warranty.

**Warrantor's Rights:** Warrantor reserves the unrestricted right to change the parts and designs of its products from time to time without notice and with no obligation to make corresponding changes in the product previously manufactured.

**Disclaimer of Warranty:** The foregoing warranties are in lieu of all other warranties, express or implied. The implied warranties of merchantability and fitness for a particular purpose are expressly disclaimed.

**Warranty Termination:** The warranty shall be void, and Wheeled Coach shall be relieved from any and all obligations hereunder, if:

1. The motor vehicle is sold or any of the components are removed.
2. The chassis is remounted, unless the remount work is performed by Wheeled Coach or by a service facility approved, in advance and in writing, by Wheeled Coach.
3. Owner fails to comply with the warranty registration requirements described above.
4. Policy expiration date

**Time Limit for requesting Arbitration:** Any action for breach of warranty must be commenced within three months following the expiration of the warranty period.

**Legal Remedies:** Any claim or controversy arising out of or relating to this limited warranty, or breach thereof, shall be settled by arbitration administered by the American Arbitration Association in Orange County, Florida in accordance with the Commercial Arbitration Rules of the American Arbitration Association. The determination of the arbitrator(s) shall be in writing and shall include an explanation of the basis for the determination. The determination of the arbitrator(s) shall be final and binding and judgment upon such determination may be entered in any court having jurisdiction.

**No Other Warranty:** Unless modified in a written statement signed by both parties, this agreement is understood to be the complete and exclusive agreement between the parties, superseding all prior agreements, oral or written, and all other communications between the parties (including without limitation any terms and conditions contained in any purchase order or sales invoice issued pursuant to the sale of the vehicle) relating to the subject matter of this agreement. No employee of Wheeled Coach or any other party is authorized to make any warranty in addition to those made in this agreement.



## **WARRANTY CLAIM FORM**



# Warranty Claim Form

Please fax back to: **800-241-5177**

Form# F - 433  
Revision Date: 07.18.2013  
Page 1 of 1

**\*\*\*This is NOT an authorization for repair. All Claims are subject to management review before approval.\*\*\***

*Form must be fully completed and returned to Wheeled Coach for review and processing. Please print clearly.*

Date: \_\_\_\_\_

## Customer Section

Name: \_\_\_\_\_

Phone: \_\_\_\_\_

Address: \_\_\_\_\_

City: \_\_\_\_\_ State/  
Province: \_\_\_\_\_

Country: \_\_\_\_\_ Zip/Postal Code: \_\_\_\_\_

Will payment be made payable to the customer?  Yes or  No

Will parts be shipped to customer?  Yes or  No

## Dealer or Repair Facility Section

Company Name: \_\_\_\_\_

Contact Name: \_\_\_\_\_

Address: \_\_\_\_\_

City: \_\_\_\_\_ State/  
Province: \_\_\_\_\_

Country: \_\_\_\_\_ Zip/Postal Code: \_\_\_\_\_

Email Address: \_\_\_\_\_

Phone: \_\_\_\_\_

Fax: \_\_\_\_\_

Will payment be made payable to the dealer/repair facility?  Yes or  No Will parts be shipped to dealer/repair facility?  Yes or  No

## Vehicle Information

**\*\*\*Only 1 vehicle per claim form\*\*\***

Work Order #: \_\_\_\_\_

Last 8 VIN #: \_\_\_\_\_

Current Mileage: \_\_\_\_\_

## Warranty Instructions

**\*\*\*Please read and check on each item.\*\*\***

Digital Pictures must be provided to support claim.

Defective parts must be returned before new parts and/or payment is released.

Estimates from repair shops and sub-contractors must be faxed with this form.

Invoices must be submitted within 30 days of approval date.

## Warranty Claim Information

Please check if this is a claim for:  Labor Only  Warranty Parts Only  Both

QTY	PART NAME OR #	DESCRIPTION OF PROBLEM
<b>LABOR</b>		
HRS:	LABOR RATE:	
QTY	PART NAME OR #	DESCRIPTION OF PROBLEM
<b>LABOR</b>		
HRS:	LABOR RATE:	

**I ACKNOWLEDGE RECEIPT OF THE PARTS AND SERVICES DESCRIBED HEREIN.**

**X**

*Customer Signature*

We are dedicated to the highest level of customer service. For any questions regarding your claim, please call our offices at 800-628-8178 between the hours of 8 am - 5 pm EST.



## PARTS AND EMS SUPPLIES



## **Parts and EMSupply Depot**

---

### **Parts**

Wheeled Coach's Customer Service (Parts Sales) is organized to meet the needs of our customers with maximum efficiency. We offer replacement parts of the same design and quality as the original equipment on your Wheeled Coach vehicle. A parts list for your vehicle is located in the back of the manual under the parts section. [www.wheeledcoach.com/parts](http://www.wheeledcoach.com/parts)

### **Terms**

Open account terms will be extended to all tax supported organizations, state, local and deferral government agencies and any other customer who completes the credit approval process.

***Open account terms must be strictly adhered to. Payment is expected within 30 days of merchandise delivery.***

All payments must be made in U.S. dollars and there will be a \$20.00 charge for all returned checks. Customers without an open account status must enclose a certified check or money order for the entire amount plus shipping charges.

### **Prices**

Prices are subject to change without notice. **Current prices are available upon request.**



## **Parts**

### **Shipping/Freight**

Wheeled Coach reserves the right to determine the most efficient means of shipping. All orders shipped F.O.B. Origin (shipping point). Shipping charges will be pre-paid and added to the invoice on all open account purchases. Whenever possible, merchandise will be shipped by UPS or Parcel Post.

### **Delivery**

All shipments leave the point of origin within 48 hours (manufactured items within 7 working days), ARO subject to stock on hand. In out of stock situations, items are automatically back-ordered and shipped as soon as possible.

### **Minimum Order**

The minimum order for open account purchases is \$25.00. Orders which fall below this minimum may be subject to a \$5.00 handling charge.

### **Business Hours**

Our normal business hours are Monday through Friday, 8:00a.m. to 5:00p.m. EST. We are closed on major holidays.

### **Returns**

Should you choose to return an item for any reason within 30 days of purchase, simply obtain a return authorization from our Customer Service Department at 800-628-8178. Merchandise must be in original packaging and in merchantable condition. The following items are non-refundable.

- ◆ Discontinued items
- ◆ Items purchased for competitive bid
- ◆ Custom/Special Order Items
- ◆ Sterile Products



## **Parts**

For proper credit, please follow the procedures below:

1. Contact Customer Service to obtain a Return Goods Authorization (RGA) number. This number must be marked clearly on the outside of the package.
2. Write your reason for returning the merchandise on the packing list. If the item is being returned because of defect, please briefly describe the nature of the defect.
3. Items being returned due to customer preference or error may be charged a 15% restocking fee.
4. Ship the merchandise back in the original packaging and include a copy of the invoice or packing list to:

**Wheeled Coach  
Attn: Shipping Department  
2737 N. Forsyth Rd.  
Winter Park, FL. 32792**

## **Claims**

Wheeled Coach carefully inspects, documents and weighs all orders prior to shipment. If there is a problem with damaged shipments, the customer should first inform the carrier. We will assist in any way in claims made against the carrier. All claims should be made within five (5) days of receipt of merchandise. In situations where the customer received a short shipment, the Customer Service Department should be notified at 800-628-8178 (for vehicle parts). In the event we made an error, we will promptly make the adjustment at no charge to the customer.



## **Parts**

### **Ordering Information**

Wheeled Coach Customer Service has several methods for placing orders. You are encouraged to use the method most convenient for you.

#### **By Mail**

Fill in the "bill" to section on the order form. This is the name of the person, business, or organization responsible for paying the charges for the order. To avoid billing errors, please use the same name that you have used on previous orders. Reference your customer number when possible for faster processing. The ship to address should be filled in with full name, address and telephone number where the order is to be delivered. We are unable to ship to Post Office Boxes. Fill in the description, item number (from the catalog), quantity, size/color, unit price and total. Please provide as much information as possible and be sure to make note of a daytime phone number where you can be reached should we have questions. If you phone in an order and later confirm it in writing, please indicate that you are only confirming an order. Otherwise you could be billed for both orders.

*Since a special order item (manufactured part) made for a vehicle, not currently in production, requires the use of exact measurements, these requests must be submitted in writing and are non-refundable.*

#### **By Phone**

Have information on the order form along with your charge card and/or customer number ready. This will help us take your order more efficiently.

#### **By Fax**

Wheeled Coach has 24-hour customer service fax line, dedicated to receiving your orders. The number is 800-241-5177.



## SERVICE CENTERS

**WHEELED COACH**  
SERVICE CENTER LISTING

<b>ALABAMA</b>						
<b>NAME OF BUSINESS</b>	<b>ADDRESS</b>	<b>CITY</b>	<b>STATE</b>	<b>ZIP CODE</b>	<b>PHONE NUMBER</b>	<b>FAX NUMBER</b>
Peach State Ambulance, Inc	130 Peach State Ct	Tyrone	GA	30290	770-486-7410	770-631-1930
<b>ALASKA</b>						
<b>NAME OF BUSINESS</b>	<b>ADDRESS</b>	<b>CITY</b>	<b>STATE</b>	<b>ZIP CODE</b>	<b>PHONE NUMBER</b>	<b>FAX NUMBER</b>
Rocky Mountain Ambulance	312 South Main	Paragonah	UT	84760	435-477-3430	435-477-3796
<b>ARKANSAS</b>						
<b>NAME OF BUSINESS</b>	<b>ADDRESS</b>	<b>CITY</b>	<b>STATE</b>	<b>ZIP CODE</b>	<b>PHONE NUMBER</b>	<b>FAX NUMBER</b>
Vital Link, Inc	1033 EMS Dr	Batesville	AR	72501	870-612-8319	870-793-3159
<b>ARIZONA</b>						
<b>NAME OF BUSINESS</b>	<b>ADDRESS</b>	<b>CITY</b>	<b>STATE</b>	<b>ZIP CODE</b>	<b>PHONE NUMBER</b>	<b>FAX NUMBER</b>
Emergency Vehicle Group	2883 East Coronado Street	Anaheim	CA	92806	714-238-0110	714-238-0120
Canyon County Fire Service	2079 Fresh Aire St	Flagstaff	AZ	86001	928-310-8584	
<b>CALIFORNIA</b>						
<b>NAME OF BUSINESS</b>	<b>ADDRESS</b>	<b>CITY</b>	<b>STATE</b>	<b>ZIP CODE</b>	<b>PHONE NUMBER</b>	<b>FAX NUMBER</b>
Emergency Vehicle Group	2883 East Coronado Street	Anaheim	CA	92806	714-238-0110	714-238-0120
Folsom Lake Ford	1 Plaza Dr	Folsom	CA	96530	916-353-2047	916-353-2128
S & C Motors Inc	35 Dolores St	San Francisco	CA	94103	415-553-4400	415-715-6988
Peck Road Truck	2450 Kella Ave	Whittier	CA	90601	562-692-7267	562-692-8390

**WHEELED COACH**  
SERVICE CENTER LISTING

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
<b>COLORADO</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Peak Motor Coach	880 Brickyard Circle Unit #140	Golden	CO	80403	303-295-0111	303-295-0362
Alpine Automotive	508 Pine River Lane	Estes Park	CO	80517	970-586-2316	
Color on Wheels	7192 N Washington St	Denver	CO	80229	303-287-7322	
Coyne's Maintenance	411 W Emerson	Holyoke	CO	80734	970-854-2310	
John Baldwin Chevrolet Buick	112 N Sunflower Dr	Holyoke	CO	80734	970-854-2211	
L.A.W.S.	7124 S Dillon Ct	Englewood	CO	80112	720-870-4965	
Larkspur Fire Protection District	9414 S Spruce Mtn Rd	Larkspur	CO	80118	303-681-3284	303-681-3201
Summit County Fleet	0218 County Road 1003	Frisco	CO	80443	970-668-3590	
Western Colorado Automotive	3190 Mesa Ave.	Grand Junction	CO	81504	970-523-1172	
Wireless Advance Communications	3901 W Service Rd	Evans	CO	80620	970-330-9415	970-330-3632
<b>CONNECTICUT</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
North Eastern Rescue Vehicles	10 Dwight Park Drive	Syracuse	NY	13209	315-414-9999	315-437-3007
<b>DELAWARE</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
DPC Emergency Equipment	335 Strauss Ave	Marydel	DE	19964	302-492-1245	302-492-1340

## **WHEELED COACH SERVICE CENTER LISTING**

**WHEELED COACH**  
SERVICE CENTER LISTING

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Feld Equipment Company, Inc.	113 North Griffith Road	Carroll	IA	51401	800-568-2403	712-792-6658
<b>KANSAS</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Nebkan, LLC	1203 Fredrich Drive	Salina	KS	67401	785-577-4455	
<b>KENTUCKY</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Muster Emergency Vehicles Inc.	135 E Third St	Calhoun	KY	42327	270-273-9220	270-273-9230
<b>LOUISIANA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Peach State Ambulance, Inc	130 Peach State Ct	Tyrone	GA	30290	770-486-7410	770-631-1930
Phillip Bonura Auto Service	1618 Veterans Blvd.	Metairie	LA	70005	504-833-8584	504-835-6048
<b>MARYLAND</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
DPC Emergency Equipment	335 Strauss Ave	Marydel	DE	19964	302-492-1245	302-492-1340
<b>MASSACHUSETTS</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
North Eastern Rescue Vehicles	10 Dwight Park Drive	Syracuse	NY	13209	315-414-9999	315-437-3007
<b>MICHIGAN</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Emergency Vehicles Plus	670 East 16th Street	Holland	MI	49423	800-320-9749	734-439-1280
<b>MINNESOTA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Everest Emergency Vehicles	1133 Rankin Street Suite D-1	St. Paul	MN	55116	800-889-6143	
<b>MISSISSIPPI</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER

**WHEELED COACH**  
SERVICE CENTER LISTING

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Peach State Ambulance, Inc	130 Peach State Ct	Tyrone	GA	30290	770-486-7410	770-631-1930
North Mississippi Medical Center	830 South Gloster St	Tupelo	MS	38801	662-841-3022	662-377-7430
<b>MISSOURI</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Feld Equipment Company, Inc	113 North Griffith Road	Caroll	IA	51401	800-568-2403	712-792-6658
<b>MONTANA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Eagles Nest Garage	201 Lewis St	Belt	MT	59412	406-277-3500	
Mobile Service Dustin Ballanntyne	2735 Aspen Dr	Big Sky	MT	59716	406-995-2100	
Fire Equipment Repair Company	3317 McGirl Rd	Billings	MT	59105	888-320-3473	406-256-3943
G & J Diesel	1739 Main St	Billings	MT	59105	406-248-6366	406-259-3556
Fox Repair LLC	768 Old Corvallis Rd	Corvallis	MT	59828	406-961-4202 406-961-8491	406-961-8491
Ron's Automotive	1900 Gold Ave	Helena	MT	59601	406-443-0550	
Ken's Equipment Repair	1224 US Hwy 10 West	Livingston	MT	59047	406-222-8005	
Riteway Service	324 N Main St	Livingston	MT	59047	406-222-3213	
Mac's Frontierland	3016 Valley Dr East	Mile City	MT	59301	406-234-2456	
AT&T Repair	56605 Hwy 93 South	Polson	MT	59860	406-883-2970	
Harmon's	3545 Sodbuster Rd	Turner	MT	59542	406-379-2584	
<b>NEBRASKA</b>						
Feld Equipment Company, Inc	113 North Griffith Road	Caroll	IA	51401	800-568-2403	712-792-6658
<b>NEVADA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Rocky Mountain Ambulance	312 South Main	Paragonah	UT	84760	435-477-3430	435-477-3796
<b>NEW JERSEY</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Campbell Supply Company	110 Sylvania Place	South Plainfield	NJ	7080	732-287-8884	908-791-3101

**WHEELED COACH**  
SERVICE CENTER LISTING

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
<b>NEW MEXICO</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Pete's Equipment Repair	1400 Broadway NE	Albuquerque	NM	87102	505-242-6969	505-242-4156
Emergency Service	2601 Pajarito Rd SW	Albuquerque	NM	87105	505-934-8728	
<b>NEW YORK</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
North Eastern Rescue Vehicles	10 Dwight Park Drive	Syracuse	NY	13209	315-414-9999	315-437-3007
<b>NORTH CAROLINA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Select Custom Apparatus	6100 North main Street	Falkland	NC	27827	252-757-3787	252-757-1639
<b>NORTH DAKOTA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Ressler Chevy, Inc	1500 Second Street NE	Mandan	ND	58554	701-663-9851	701-667-8837
Westlie Ford Linc/Merc	500 S Broadway	Minot	ND	58702	701-857-1631	701-857-1610
<b>OHIO</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
911 Fleet and Fire Equipment	11 Lendale Drive	Florence	KY	41042	877-605-2378	859-371-0131
Ambulance Maintenance	2512 Norwood Ave	Cincinnati	OH	45212	513-221-7969	513-221-4135
Bus Service Inc	3153 Lamb Ave	Columbus	OH	42319	614-471-2877	614-471-8801
<b>OKLAHOMA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Nebkan, LLC	1203 Fredrich Drive	Salina	KS	67401	785-577-4455	
<b>PENNSYLVANNIA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Pfund Superior Sales	221 Chester Dr	Lower Burrell	PA	15068	724-339-1600	724-339-7509
Campbell Supply Company of PA	85 Griffith Road	Muncy	PA	17756	732-287-8884	570-546-4984

**WHEELED COACH**  
SERVICE CENTER LISTING

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
<b>RHODE ISLAND</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
North Eastern Rescue Vehicles	10 Dwight Park Drive	Syracuse	NY	13209	315-414-9999	315-437-3007
<b>SOUTH CAROLINA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Peach State Ambulance, Inc	130 Peach State Ct	Tyrone	GA	30290	770-486-7410	770-631-1930
Code 3 Outfitters	1640 Priceville Rd	Gilbert	SC	29054	803-622-6331	803-892-5925
Spartanburg County Fleet Services	290 Broadcast Dr	Spartanburg	SC	29323	864-596-3586	864-596-3588
<b>SOUTH DAKOTA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Feld Equipment Company, Inc	113 North Griffith Road	Carroll	IA	51401	800-568-2403	712-792-6658
<b>TENNESSEE</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Peach State Ambulance, Inc	130 Peach State Ct	Tyrone	GA	30290	770-486-7410	770-631-1930
<b>TEXAS</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Southwest Ambulance Sales	1715 Wall Street Suite 101	Garland	TX	75041	903-440-2358	903-440-2358
Absolute Vehicle Service	213 Bennington St	Houston	TX	77022	713-490-5416	713-490-5417
Boggus Motors	2521 S Expressway 77/83	Harlingen	TX	78550	954-423-2580	956-423-3117
Municipal Trucks and Equipment	1120 Produce Rd	Hidalgo	TX	78557	956-843-9661	956-843-7861
Pollard Ford Friendly Ford	3301 South Loop 289	Lubbock	TX	79423	806-797-3441	806-795-3414
<b>UTAH</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Rocky Mountain Ambulance	312 South Main	Paragonah	UT	84760	435-477-3430	435-477-3796
<b>VIRGINIA</b>						
NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER

**WHEELED COACH**  
SERVICE CENTER LISTING

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Select Custom Apparatus	2742 Mary Linda Ave, NE	Roanoke	VA	24012	540-344-2229	540-344-2203

**WEST VIRGINIA**

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Pfund Superior Sales	221 Chester Dr	Lower Burrell	PA	15068	724-339-1600	724-339-7509

**WISCONSIN**

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Everest Emergency Vehicles	1133 Rankin Street Suite D-1	St. Paul	MN	49423	800-889-6143	

**WYOMING**

NAME OF BUSINESS	ADDRESS	CITY	STATE	ZIP CODE	PHONE NUMBER	FAX NUMBER
Peak Motor Coach	880 Brickyard Circle Unit #140	Golden	CO	80403	303-295-0111	303-295-0362
High County V.E.C.T.R.C.	6030 Big Horn Ave	Sheridan	WY	82801	307-674-4006	
Owl Creek Manufacturing Co.	1080 Raintree Dr	Riverton	WY	82501	307-856-6346	



## WORK ORDER



## ELECTRICAL DRAWINGS

Kempsville Volunteer Rescue  
 Kevin Lipscomb  
 5145 Ruritan Court  
 Virginia Beach, VA 23462

Select Custom Apparatus, Inc.  
 Todd Pittman  
 1 Unit  
 Vin# 1GB6G5CL7F1165544,  
 SHA

Quote No: CF031915-  
 03/19/2015 14:47:18

WO#: 353031

Page 1

PART NO	DESCRIPTION	QTY	ID	REF. NO
	<b>== Chev, 159"WB, 170" Module, #3170C - 620.004 04/09/13</b> == Eng: Ref only previous unit #333075.	1	WC	
00-00-0247	***** OUTSIDE GRAPHICS- 24/7 ***** <b>CHASSIS:</b>	1	WC	00000247
00-30-5659	<b>Chev, 2015, 6.6L Diesel Cutaway,159" WB,170"</b> <b>Module,#3170C 2015-4</b> (with DEF) G4500	1	WC	00-30-5659
	170X95 Auto throttle OEM <b>Y537</b> Alternators OEM Rear Swaybar OEM			
1A-10-0100	<b>CHASSIS ADD-ONS:</b> Lettering, By Outside Vendor (SO-LETTR-353031) SR #215193F RFP 03/10/15	1	WC	9980011
	Location: 24 / 7 Graphics # of Days: 2 Unit # 924	1	WC	
1A-20-6370	<b>MIRRORS:</b> <b>Mirrors, Velvac #715917,2020XG,C3,Remote/Heated,</b> <b>Black,ILOS</b> <b>With switch on driver's side bracket.</b> <b>Y556 w21 w23</b>	1	WC	1A-20-6370
1F-30-1600	<b>SUSPENSION:</b> Suspension Upgrade, Front, Timbren, GMFK15WC, G4500	1	WC	1F-30-1600
1F-30-1700	Upgrade Handling Kit , G4500, Bilstein shocks	1	WC	1F-30-1700
1K-40-3520	<b>ALTERNATOR(S):</b> <b>Alternator, Add'l, Chevrolet G3500,4500,Pentax PX422V10,</b> <b>IATS</b> <b>Y150</b>	1	WC	1K-40-3520
1P-50-5055	<b>THROTTLE CONTROLS:</b> <b>High Idle Controls,OEM Switch in Steering Wheel,Chev T3/T7</b>	1	WC	9980003
5U-10-0001	<b>HEATING/AIR CONDITIONING:</b> <b>Special Circuit for AC system IATS</b> (SO-AC-353032) Ref (SO-AC-323049) SR #215106F	1	WC	
	Special Circuit for rear A/C. Rear HVAC fan is to be wired to ignition and not battery timer. Rear fan to shut off when ignition is turned off. Special Circuit for rear A/C. Rear HVAC fan is to be	0501 0503 0505 B513 B515		

PART NO	DESCRIPTION	QTY	ID	REF. NO
	wired to ignition and not battery timer. Rear fan to shut off when ignition is turned off. ILOS		Y119@0520	
5U-12-5000	Hoses, Heater, No Max, to Rear, STD	1	WC	9980003
5U-70-0407	Heat/AC, C3, Combo, w/Coolbar with thermostat in action area.  12vdc External Condenser, in Coolbar mounted on front of module.	1	WC	5U-70-0407
	Painted: Chevy white: FLNA4738 (2015 Chevy- FLNA41383)			
	Note: Rear Heat/AC Fan shall cut off when ignition switch is cut off. (Do Not Run Thru 5 Minute Timer)			
	<b>WHEELCOVERS:</b> == Type 3 Chev - 170" Module - #3170C - 620.004 04/09/13	1	WC	
	==	1	WC	
35-05-0159	Type 3 - 170" Module - Standard Model Series	1	WC	35-05-0159
35-05-0170	Type 3 - 170" Module General Body Construction Spec	1	WC	35-05-0170
35-05-0180	Type 3 - 170" Module Vehicle Body Structure Spec	1	WC	35-05-0180
35-05-2000	INTERIOR CABINETS -Duralite Construction, Wood	1	WC	35-05-2000
	<b>MODULE BODIES:</b>	1	WC	
	BS NOTE install front stone guards & corner caps & rear kickplate with heavy rivets ILOS fasteners. Customer Request			
	Note std compt lights			
35-10-0066	Aisle Space, 50 inches, STD, Type 3, Std Model Only (49-50" actual)	1	WC	35-10-0066
	Install Black(non-logo !!!) anti-slip on curbside and rear stainless steel thresholds, ILOS !!!!!!!			
35-10-0101	10 Inch Extended Cab, Chev 165"-170", Std	1	WC	35-10-0101
35-10-0320	Module Skins, 0.125 inch	1	WC	35-10-0320
35-10-0921	Headroom, 72 in. Type 3 - 170" module, #3170C/F	1	WC	35-10-0921
35-11-0100	Drop Skirt, 3 inches forward of rear wheels, Both sides	1	WC	35-11-0100
37-00-0350	Coating, Scorpion X02, Per Compt, (1), Std light gray	6	WC	3F-70-0520
37-00-0360	Coating, Scorpion, Rear Bumper Supports	1	WC	3F-70-0525
37-00-0500	Compartments, Std, Floor 3" Drop Down from Door Opening  compts B1,B2,D,& E are STD drop down floor	4	WC	37-00-0500
37-00-0600	Compartments, Sweepout - Compartments "A" & "F" are standard sweep out	2	WC	37-00-0600
37-00-0710	Duraseam Doors, with Hidden Jambs - Magnetic door switches, Standard	1	WC	37-00-0710
39-A0-0472	A, Compt, Split, Hoseline, 170" Module, 72" H.R, T3	1	WC	39-A0-0472
39-B1-0153	B1, Compt, Std, 170" Module, 72" HR, #3170C/F	1	WC	39-B1-0153
39-B2-0172	B2, Compt, Std, 170" Module, 72" H.R, T3	1	WC	39-B2-0172
39-D0-0353	D, Compt, 3/4, 170" Module, 72" H.R, #3170C/F	1	WC	39-D0-0353

PART NO	DESCRIPTION	QTY	ID	REF. NO
	<b>with single door !!!</b>			
39-E0-0103	E, Compt, Std, 170" Module, #3170C/F With Double Doors  Dog Leg for tilt out sharps container	1	WC	39-E0-0103
39-F0-0153	F, Compt, 3/4, 170" Module, 72" H.R, #3170C/F  Eng reference 333075 SL inlet is forward of "F"	1	WC	39-F0-0153
3A-10-4000	<b>FUEL FILL HOUSING:</b> Housing, Fuel Fill, Cast, F3/C3/C4500/CM+(7)	1	WC	1000774
3F-10-1103	<b>EXTERIOR DIVIDERS:</b> Divider, Vertical, B2, 3/16" Thick, Std Model only compt "B2", approximately 8.25" from wall 3  NOTE: divider to be 14" deep	1	WC	3F-10-1103
3F-10-5500	Divider, Vertical, Custom, (1), Mounted:  Compartment "D" - Install two (2) adjustable vertical dividers on wall #2. The dividers are to be 3/16 inch thick, 40 inches tall and 7.5 inches deep. Install the first divider 11 inches from wall #1 (just behind the spring shackle dogleg) and the second divider 22 inches from wall #1. Shelf track to run the width of compartment from the shackle box back.  *Note: Rubber finish and dri-deck on the floor per the work order.	2	WC	1100333
	Note: Mount divider's on adjustable track only, do not mount to floor or shelf. Edge Trim on Dividers			
	*See photos submitted with the order file.			
3F-20-1000	<b>SHELF- COMPARTMENT "A"/RF CABINET:</b> Shelf, Adjustable,(1), RF Cabt Custom (1) RF Compt "A"/cabt "E2" area, centered.	1	WC	1101600
3F-20-1500	Shelf, Additional, Adjustable, (1), RF Cab Custom  Locate in RF compt "A"/cabinet "E1", installed mid-height. w/ aluminum shelf track, std.	1	WC	1101601
3F-30-1000	<b>SHELF- COMPARTMENT "B"/"B1":</b> Shelf, Custom, Adj, Cmpt "B", (1)  Locate in compt "B-2", Mid height between vertical divider and wall #3.  Install 20" of adjustable aluminum shelf track to allow for 10" of adjustment up or down	1	WC	1101602
	<b>SHELF- COMPARTMENT "D":</b>	1	WC	

PART NO	DESCRIPTION	QTY	ID	REF. NO
3F-40-3001	Shelf, Custom, Adj, Cmpt "D", (1), Std Model only at bottom of cabt "H" inside/outside access opening.  w/ aluminum shelf track the entire height of the compt.  Note: Second (IATS) shelf to be installed above this shelf.	1	WC	3F-40-3001
3F-40-3500	Shelf, Custom, Adj, Cmpt "D", (1), Addtl'  Install above existing shelf in compt D, install dri deck on shelf per wo.	1	WC	3F-40-3500
3F-50-0200	<b>SHELF- COMPARTMENT "E":</b> Shelf, Fixed, Cmpt E, for Electrical storage, No Divider  consolidate electrical equipment on the ceiling and on shelf as close to wall #1 & wall #2 as possible. Install an expanded alum cover to be as minimal as possible and still allow for protection and ventilation.  Inverter to be on ceiling. Suction pump and airhorn compressor to be on shelf below inverter.	1	WC WC	1101621  E1, E2 LOOK!!
3F-50-2000	Shelf, Custom, Adj, Cmpt "E", (1)  Install mid height w/ aluminum shelf track.	1	WC	1101606
3F-70-0400	<b>COMPARTMENT / SHELF, MISC:</b> Bracket, Stairchair Mount, State Location :  Per Engineering design Compt: "D", to hold a Stryker Model 6252 Stair Chair  *Stair chair dimensions (folded) Height - 37.5 inches Width - 20.5 inches Depth - 8 inches  *Add additional mounting screws across the top. See drawing on server for all mounting screw locations .  * Includes a Seat Belt across the stair chair installed approx 30" above the bottom of stair chair box.  **As installed on previous unit - 323049  **See Reference photos submitted with the order file.	1	WC WC	2350076
3F-70-1200	Upgrade, Shelf/Div to 3/16", W/1" Lips, ILOS  for both compt D dividers	2	WC	1101618
3F-70-1700	Reinforcement, Door, SCBA Bracket  Reinforce compt door "D" for Stryker stair chair bracket.	1	WC	1109966

PART NO	DESCRIPTION	QTY	ID	REF. NO
	<b>WINDOWS:</b>			
3P-10-5201	Window,Upper,CS/Slider,Rr/Fixed,Priv Tint,PAN,Std Model only	1	WC	3P-10-5201
3P-11-2050	Window, Mod, C/S Over S/B, Privacy Tint, 12" x 36"	1	WC	3P-11-2050
	<b>MODULE ENTRY DOOR HANDLES:</b>			
3U-30-8010	Handles, Patient Entry, Trimark, Chrome "SafePass" ILOS	1	WC	3U-30-8010
	<b>EXTERIOR COMPT DOOR HANDLES:</b>			
3U-40-0512	Handle, Module Compt, Trimark, Chrome, ILOS	7	WC	3U-40-0512
3U-40-0514	Latch, Paddle, Trimark, For Trailing Door compt "E"	1	WC	3U-40-0514
	<b>ELECTRICAL LOCKS SWITCHING:</b>			
3U-60-1900	Switch, Stealth, Cab Doors (not Dodge) In grille on passenger's side.	1	WC	1800024
	For cab doors only !			<i>y529</i>
	<b>DOOR HOLD OPENS:</b>			
40-10-1611	Hold Opens,Rr Doors,Cast,Grabber (2), 5.5", Std Model only - Installed so doors will open as wide as possible - Mount at top of doors.	1	WC	40-10-1611
40-10-7206	Hold Open, Gas Strut,Comp Door Only,110 Degree,30lb,ILOS  SR# 215106F Note: For all compartment doors except battery compartment.  Note: Install so that doors open wide enough to give full access to inside width of the compartment	7	WC	40-10-7206
40-10-74SR	Hold Open, Gas Strut, Dampening, C/S Door, ILOS (SO-STRUT-353032)Ref (SO-STRUT-313070)SR #215106F	1		
	<b>DOOR PANELS:</b>			
45-10-45SR	Panels, Entry, Stainless Steel, Brushed (SO-PANEL-353032) SR #213779F  Panels, Entry, Stainless Steel, Brushed. NO reflective chevrons. Center section to be stainless steel.	1	WC	
	<b>FENDER FLARES:</b>			
4A-10-1001	Flare, Fender, Rear, Black Rubber, Std Model only	1	WC	4A-10-1001
	<b>BUMPERS:</b>			
4F-10-9502	Bumper,Rear,w/ Skids & Flip Up,C3,w/LED DOT lights,Std Model  <i>y3(x2)</i>	1	WC	4F-10-9502
	w/ Gator Grip on flip up step, STD			
	<b>RUNNING BOARDS:</b>			
4K-10-5611	Running Boards, Chev Type 3,w/Gator Grip,#3153C/65C/70C only	1	WC	4K-10-5611
	<b>SKIRTRAILS</b>			
4U-10-9210	Skirtrails,Extruded Rubber w/Refl Tape,170-175" Module	1	WC	4U-10-9210
4U-10-9500	Tape, Reflective, White, For Extruded Rubber Skirtrails	1	WC	4U-10-9500
	<b>STONE GUARDS (FRONT):</b>			
4U-11-0001	Rivots ILOS Hardware (SO-RIVOT-353031)	1	WC	

PART NO	DESCRIPTION	QTY	ID	REF. NO
	Use RIVOTS ILOS hardware to install front stone guards, front corner caps, rear corner caps and rear kickplate.			
4U-11-5000	Stone Guard, Front, Dia Plate, 17"	1	WC	1100141
	Install with rivets ILOS fasteners			
	<b>KICKPLATES (w/STD. REAR STONEGUARDS):</b>	1	WC	
50-10-2105	Kickplate, W/Recess, No Light Holes,Duraseam	1	WC	1100137D
	Install with rivets ILOS fasteners			
	<b>MODULE ADD-ONS:</b>	1	WC	
55-10-3000	Holder, License Plate, Cast, C30002, STD - in rear kick plate, centered	1	WC	1600055
	<i>y3 (x2)</i>			
	<b>MODULE STEPWELL:</b>	1	WC	
55-90-6020	Stepwell, Double Step, ILOS, Chev Type 3	1	WC	55-90-6020
	<b>MUD FLAPS:</b>	1	WC	
5A-10-1001	Mud Flaps, Rear, W/Logo, Std Model only	1	WC	5A-10-1001
	<b>PAINT/BELTS:</b>	1	WC	
	Paint shop heads up see painted opti com bracket.			
	OEM White, Dodge/Chevy white: FLNA4738 (2015 Chevy- FLNA41383)			
5F-10-0800	Paint, OEM White Sikkens Crossover Codes Dodge/Chevy white: FLNA4738 (2015 Chevy- FLNA41383)	1	WC	5F-10-0800
5F-10-1200	Paint Belt, Single, 20% Upch Green Metallic - FLNA 90374.	1		
	Delete the paint belt from the rear of the ambulance, including the rear doors.			
	*The paint belt should wrap around the corner of the rear upright extrusions but stop there.			
	**See Reference photos submitted with the order file. Match Previous WO# 313070			
	<b>DECAL PACKAGES INSTALLED:</b>	1	WC	
	No roof star			
	<b>ELECTRICAL SYSTEM: Circuitboard</b>	1	WC	
	This unit has std compt lights, NO elec locks on module and NO clock.			
	<b>BATTERY SWITCHES:</b>	1	WC	
60-22-1020	Battery Sw, 5min Timer, C3/CM+, Fig 5B Activated thru OEM ignition switch timer function upon ignition "off".	1	WC	60-22-1020
	<i>y19</i> <i>y571</i> <i>y573</i> <i>y575</i> <i>w350</i>			
	Momentary rocker switch, on driver's side of cab console, to function as timer shutoff and also reactivate timer.			

PART NO	DESCRIPTION	QTY	ID	REF. NO
	(1) OEM battery under hood (2) Delco 1150 batteries in RF Lower Compartment			
	ONLY Module power is turned on/off. Chassis power is NOT turned off.			
60-30-5000	Battery, Delco, 1150, (1) Each	2	WC	1200230
	COAX CABLES/ANTENNAS:	1	WC	
65-10-3001	Coax Cable, RG58/U, (1), Additional, Std Model only	1	WC	65-10-3001
	Module roof port #1 to ECC - with 6' pigtailed		COAX #1	
65-30-1000	Coax Cable, RG58/U, (1) Each, Standard, Terminate:	1	WC	1850001
	Module roof port #2 to ECC - with 6' pigtailed		COAX #2	
65-30-1500	Coax Cable, RG58/U, (1) Each, Additional, Terminate:	3	WC	1850002
	(1) From module roof port #3 in bulkhead cabinet at CB.		COAX #3	
	(1) From module roof port #4 in bulkhead cabinet at CB.		COAX #4	
	(1) From module roof port #5 in bulkhead cabinet at CB.		COAX #5	
	All with 6' pigtailed			
	12VDC POWER SOURCES & OUTLETS:			
6A-22-200E	Outlet, Cigar Lighter, STD, (2), A/A	1	WC	
	*Ignition/Shoreline Hot* ILOS same as 333075	1	WC	6A-22-200E
	(1) in A/A		0107 (x2)	
	(1) in the RF cabt "E-3", wall #2, upper right.			
6A-22-2500	Outlet, Cigar Lighter, IATS, (1)	1	WC	1801258
	*Ignition/Shoreline Hot*			
	Install on extended A/A wall, forward of 110VAC Outlet, same level as SSCOR canister.			
6A-22-5100	Outlet, Dual USB Port, 5VDC, 2.1Amp output, Standard	1	WC	6A-22-5100
	Passenger side of the ECC, per Elec Eng Ign sw hot		Y110	
6A-22-5150	Outlet, Dual USB Port, 5VDC, 2.1Amp output, IATS	1	WC	6A-22-5150
	Location : A/A Wall, Ign sw hot		Y110	
	POWER SOURCE, 12 VDC			
6A-23-0910	Power Source, 12VDC, 20A, Ignition/Shoreline Hot	1	WC	6A-23-0910
	(1) pre-wire coil and tagged in action area		0554	
	(1) pre-wire coil and tagged behind driver's seat.		0554	
6A-23-1200	Power Source, 12 VDC, 15 Amp (+)	1	WC	1800005
	Ign hot, terminating behind A/A with 6 foot tails (hot and ground).		W152	
6A-23-1201	Power Source, 12 VDC, 15 Amp (+), Std Model only	1	WC	6A-23-1201
	Terminating in the (ECC) w/6 foot tails (hot and ground).		W150	
			14gblk	

PART NO	DESCRIPTION	QTY	ID	REF. NO
	Switch: Constant hot.			
6A-23-3001	Power Source, 12 VDC, 30 Amp (+-), Std Model only Constant Hot, to terminate behind the pass seat with 6 foot tails (hot and ground)	1	WC	6A-23-3001 <i>B 150</i>
6F-80-5400	SHORELINE INLETS: Shoreline Indicator Light, 110VAC, LED, Green near shoreline inlet	1	WC	1200601 <i>12/3 SO Yellow</i>
6F-90-4101	Super Auto Eject, 20 amp, ILOS, White, Std Model Only Install Super Auto Eject on the drivers side just forward of comp F in 10 inch extension.	1	WC	6F-90-4101 <i>12/3 SG Y121</i>
	**See Reference photos submitted with the order file.			
	110 VAC OUTLETS:	1	WC	
6K-40-1000	Outlet, 110 VAC, Duplex, (2), STD (1) outlet in action area (1) outlet in the RF Cabinet "E-3", wall #2, upper right.	1	WC	1800001 <i>12/3 SO #1</i> <i>12/3 SO #4</i>
6K-40-2000	Outlet, 110 VAC, Duplex, IATS (1) installed on extended A/A wall, centered above wheel well, same height as SSCOR suction canister. (1) installed in bulkhead wall behind passenger's seat. (1) installed in "H" Cabinet	3	WC	1800219 <i>12/3 SO #2</i> <i>12/3 SO #5</i> <i>12/3 SO #3</i>
	ENGINE BLOCKHEATER:	1	WC	
	OEM Per ECN 1804			
	INVERTERS:	1	WC	
7F-10-9000	Inverter/Charger, Vanner 20-1050CUL-DC Compt E Ceiling 1050 watt inverter to power all 110 VAC receptacles w/ switch in the A/A panel.	1	WC	1296422 <i>Y110, Y539</i> <i>0554</i>
	CONSOLES:	1	WC	
7U-10-7000	Console, A/A, Wood, Angled	1	WC	1301700
7U-11-0050	Console, ECC, Chevy Cutaway, STD, W/Map Storage, 2007+ Painted black.	1	WC	7U-11-0050
7U-11-72SR	Console, Custom, Dove Gray Mica, Chevy, IATS (SO-CC-353032) Ref (SO-CC-323049) SR #215106F CC Dove Gray Mica, Chevy, IATS Note: Slanted area at front of console for customer installed two way radios with holes cut for (1) Motorola XTL2500 & (1)Motorola CDM1250, two (2) cup holders in middle section, and map book storage at rear with	1		

PART NO	DESCRIPTION	QTY	ID	REF. NO
	binning strips and four (4) lexan dividers.			
D0-10-6500	<b>GAUGES/METERS:</b> Alarm, Low Voltage, Audio/Visual Light in cab console and Buzzer in cab.	1	WC Y552 Y553	D0-10-6500
D0-10-7000	Voltmeter, OEM, In Dash	1	WC	
D5-10-1000	<b>ENGINE HOUR METER:</b> Meter, Engine Hour, Datcon 100-223 12-24 mounted on driver's side of cab console.	1	WC W11a BH14Cn	1530040
DA-11-3000	<b>SEQUENCING/LOAD MANAGEMENT:</b> Sequencer,Kussmaul,Load Manager, W/ Switch W/Sequencer switch and override on console	1	WC Y104 Y105 Y107 Y109	1530049
DF-20-5600	<b>SIREN SPEAKERS:</b> Speakers,(2),Whelen, SA-315, In Grille area, C3/C7	1	WC Y70 Y71	DF-20-5600
DK-10-1400	<b>SIRENS:</b> Siren, Whelen, WS-295-SLSA1, Standard	1	WC Y70 Y73	1700297
	Switch Siren - Sequencer Hot (only)		Y105 R B2	
DP-10-9100	<b>AIR HORNS:</b> Airhorn, Buell, 5440 comp, 10 & 12 inch horns	1	WC B10 BIK10	1700533
	*Install compressor on fixed shelf in compartment "E", below inverter, as close to wall #1 & wall #2 as possible. Install an expanded alum cover to be as minimal as possible and still allow for protection and ventilation.*		Y5 Y110	
DP-10-9400	Brackets, Airhorn, Thru Bumper, C3/C7	1	WC	DP-10-9400
DP-11-1000	<b>SWITCHING OPTIONS FOR AIR HORNS &amp; SIRENS:</b> Switch, Siren/Horn Thru Horn Ring	1	WC 02 03	1700622
DP-11-5000	Valve, w/ Lanyard, For Airhorns:	1	WC	1530033
	Locate in the center of the headliner so it is a equal reach from the driver or passenger.		BIK10 (sw) Y5 wires drop at LF OF LINER	
	**See Reference photos on file.			
DU-10-3000	<b>SWITCHES ADDITIONAL:</b> Circuit, Auto Reset, For: See opticom momentary switch	1	WC Y59, W11	1200024
DU-10-6000	Switch, Momentary, Rocker, (1) Momentary rocker, IATS opti com switch, Additional momentary switch to activate the opticom, ign hot and in gear,once they shift it back into park it will automatically shut the opti com off.	1	WC	1539630
F0-11-1500	<b>AUDIBLE ALARMS:</b> Alarm, Back-up, No Cutoff Switch	1	WC	1700708
F3-10-01SR	<b>VEHICLE EXTERIOR LIGHTING</b> Opticom, GTT LED on Cool Bar (SO-OPTI-353032)SR #215106F OPTICOM,GTT #794 LED	1	WC	

PART NO	DESCRIPTION	QTY	ID	REF. NO
	GTT- Opticom-Com Emitter with OEM opticom mount bracket, Eng note include std opti com switch. <i>w10 w11</i>			
	Mount bracket to the top of the cool bar, PAINT bracket white to match the cool bar.			
FF-10-6001	Cab Headliner Removal Fee, Type 1/3  for LED cab dome light and lanyard installation.	1	WC	FF-10-6001
	<b>WHELEN FLANGES FOR 9E SERIES LIGHTS:</b>	1	WC	
FU-38-0100	Flange, Super 9L, #90FLANGC, (1) Each	12	WC	FU-38-0100
FU-38-0500	Flange, 97 Series, (1) Each	8	WC	1550023
	<b>WHELEN SUPER LED 9L FLASHING LIGHTS T1,3,9:</b> Front warning lights Ref 313070	1	WC	
FU-40-0107	Light, 9L, Super LED, Whelen, Red,Clear Lens,IATS, 90RR5FCR Eng Ref 313070	2	WC	FU-40-0107
	One (1) on each side of the center white light, cool bar.  Flash: Action Scan SW: Pri/Sec Flasher: Internal			<i>Y57(x2)</i>
	<b>WHELEN SUPER LED 9L STEADY BURN LIGHTS T1,3,9:</b>	1	WC	
FU-40-0505	Light, 9L, Super LED, Whelen, Red,Clear Lens,ILOS, 90RR5SCR (2) each side of module, standard locations Flash: Comet Flash	8	WC	FU-40-0505
	(2) on front in upper outer corners Flash: Comet Flash, alternating with center white light			<i>05 (x4) 05A (x4)</i>
	(2) on rear in upper outer corners Flash: Comet Flash, alternating with center amber light			
	SW: Primary/Secondary Flasher: Whelen AFM-560			
FU-40-0507	Light, 9L, Super LED, Whelen, Red,Clear Lens,IATS, 90RR5SCR (2) on rear of module in line with UPPER portion of rear windows to be visible when rear doors are opened.	2	WC	FU-40-0507
	Flash: Through AFM 560 flasher SW: Primary/Secondary			<i>013 (x2)</i>
	These will also be IATS brake lights to steady burn anytime brakes are applied			
FU-40-0510	Light, 9L, Super LED, Whelen,Amber,Clear Lens,ILOS, 90AA5SCR (1) on center, rear.	1	WC	FU-40-0510
	Flash: Comet Flash, alternating with upper, outer, red lights. SW: Primary Only Flasher: Whelen AFM-560			<i>09</i>

PART NO	DESCRIPTION	QTY	ID	REF. NO
FU-40-0514	Light, 9L, Super LED, Whelen, Clear, ILOS, 90CC5SCR  Install centered on cool bar.	1	WC	FU-40-0514
	09			
	Flash: Comet Flash, alternating with upper, outer, red lights SW: Primary Only Flasher: Whelen AFM-560			
	<b>WHELEN 6E FLANGES:</b>	1	WC	
H2-10-0050	Flange, 6E Series, (1) each	2	WC	H2-10-0050
	<b>WHELEN SUPER LED 7L LIGHTS W/FLASHERS:</b>	1	WC	
H5-58-0320	Light, 7L, Super LED, Whelen, Red, Clear Lens, ILOS,70R02FCR  (2) as front intersection lights.	2	WC	H5-58-0320
	Y55 (x2)			
	Flash: Action Scan SW: Primary Only Flasher: Internal			
H5-58-0325	Light, 7L, Super LED, Whelen, Red, Clear Lens, IATS,70R02FCR  (2) On module over rear wheel wells as intersectors. Flash: Action Scan SW: Primary Only Flasher: Internal	4	WC	H5-58-0325
	Y55 (x2)			
	(2) Angled ends of cool bar top, Eng see IATS scene lights also. Flash: Action Scan SW: Primary Only Flasher: Internal			
	Y55 (x2)			
	<b>7 SERIES SUPER LED STEADY BURN LIGHTS:</b>	1	WC	
H5-58-0520	Light, 7L, Super LED, Whelen, Red, Clear Lens, ILOS,70R02SCR  (2) as grille lights.	2	WC	H5-58-0520
	09 (x2)			
	Flash: Comet Flash, synchronized w/ upper center white light SW: Primary Only Flasher: Whelen AFM-560			
	<b>WHELEN FLANGES 7 SERIES LIGHTS</b>	1	WC	
H5-59-1000	Flange, Whelen, For 7E lights (1), Mods - For 7 Series Lights	6	WC	H5-59-1000
H5-59-1010	Housings, Whelen, 7E Grille Lights, Pair, Chev T3/CM - #7GCV6LRC	1	WC	H5-59-1010
H5-59-1012	Housings, Whelen, 7E Intersection Lights, Pair, Chev T3/CM - #7ICV6LRC	1	WC	H5-59-1012
	<b>CORNER CAP LIGHTS</b>	1	WC	
HF-11-6049	Corner Cap Lts,Multi-LED,(2)Amber(2)Red w/Flashers,Std Model - High intensity flashers; thru separate switch	1	WC	HF-11-6049
	Y35,Y37,Y39			
	<b>CLEARANCE LIGHTS:</b>	1	WC	
HF-12-9100	Light, Clearance, Amber LED, ILOS	3	WC	1570073-1
	Y3			
	Install on the front, center of module.			
HF-12-9601	Light, Clearance, Red LED, Std Model only - On rear as clearance lights.	3	WC	HF-12-9601
	Y3			

PART NO	DESCRIPTION	QTY	ID	REF. NO
HK-20-110M	<b>SCENE AND LOAD LIGHTS</b> Light, LED, Whelen, Clear, #9SC0ENZR, 24 LED, ILOS, (1) W Chrome Flange (2) each side as scene lights      W201    W251 W203    W253 (2) on rear as load lights      W205	1 6	WC WC	HK-20-110M
HK-20-12SR	Light, LED, Whelen, Clear, 70C0ELZR, (1) IATS (SO-LOAD-353032) Ref (SO-LOAD-313070) SR #215106F (2) on angled ends of cool bar, on bottom under warning lights. Switched with Left & Right Flood Light Switches.      W201 W203	2		
HK-20-7100	Scene Lights, C/S, Activated W/C/S Door Details: Curb side scene light will activate when C/S entry door is open.      W202 @ W203	1	WC	1600727
HK-20-7400	Rear Scene Lts Activate W/Rev Gear. Details: The rear scene lights on both sides of the module will activate when the vehicle is in reverse gear.      026 @ W251 W253 * The rear load lights wired to reverse is standard.	1	WC	1600721
<b>FLASH PATTERN WARNING LIGHTS:</b>				
HU-11-2500	Brake Override Circuit      013 For lights at rear windows, wired per Elec Eng.	1 1	WC WC	1200900
<b>WARNING LIGHT FLASHERS</b>				
HU-24-7500	Flasher, Whelen AFM-560, ILOS      Y13 Y15 Set flash pattern to Comet Flash      05 05A 07 09	1 1	WC WC	1521496
<b>HEADLIGHT FLASHER:</b>				
J0-42-2200	Flasher, Headlights, Whelen SSFPOS	1	WC	1550956
J0-42-7100	Circuit, Neutral Power, C1/C3/F1/F3/D1, 2001+	1	WC	1530072
J0-42-8200	Circuit, VA State Law Headlight Flashers      Y3,016,01B,Y21	1	WC	1550001
<b>HARNESS LAYOUT:</b>				
J5-12-1000	Halogen, " E " Spec, Whelen      020 021 023	1	WC	H960125
<b>STOP/TURN/TAIL LIGHTS:</b>				
JA-10-1200	Light, License Tag, (2), Trucklite #15205 LED, ILOS      Y3	1	WC	JA-10-1200
JA-13-8010	Tailights, Whelen, LED Stop/Tail, Halogen/Backup, w/ WC logo - Brake- Alert Flash before the steady Burn - Turn Arrow – sequential arrow	1	WC	JA-13-8010
JA-14-360E	Turn Signal, Whelen, Amber Arrow, 6LED, IATS, Pair Mounted on front of module, one each side, just below 9L red LED lights in upper outer corners.      W21 W23 Flash Pattern: Sequential Flash	1	WC	JA-14-360E
JA-14-70SR	Turn/Brake/Turn Indicators, LED, Trucklite (SO-IND-353032) Ref (SO-IND-333075) SR #215106F      W17 W21 W23 Locate Trucklite LED Turn/Brake/Turn Indicator Lights, between drop at rear of liner	1		

PART NO	DESCRIPTION	QTY	ID	REF. NO
	two (2) glove box holders in glove box cabinet, located above rear patient entry doors.			
JK-10-5500	<b>MODULE OVERHEAD HEADLINER:</b> Headliner, Flat, Expanded PVC, Std	1	WC	JK-10-5500
JP-10-0001	<b>12VDC INTERIOR LIGHTING:</b> Dome Light, Cab, Whelen Red/White LED IATS (SO-LT-353032)Ref (SO-LT-313070) SR #215106F  Ignition sw hot, switches on light head used to turn light on & off.  Whelen # 60CREGCS, 6 inch round surface mount light super LED Red/White.  Installed in the headliner forward of the OEM dome light, centered.  Look>> do not to interfere with the air horn lanyard.	1	WC	Y110 wires drop at front of liner
JP-10-0612	Lights, Dome, LED, Whelen, Chrome Flng, (4) S/S, (3) C/S, ILOS	1	WC	JP-10-0612
JP-10-7015	Light, Stepwell, LED, Whelen #TOCACCCR, 2in, ILOS	1	WC	JP-10-7015
JP-10-8030	Light, Intertek ZY-PIR38 12v, LED, In A/A ILOS	1	WC	JP-10-8030
JP-16-1206	<b>12VDC INTERIOR LIGHTING - FLUORESCENT:</b> Lights, Fluor, Thnlite #766,18IN,(3),W/Sw, Std Model Only (must change to surface mount if "ducted" A/C system)  Recessed, 18"  Mount in liner positions F, G, and H	1	WC	JP-16-1206
	W301 W351 With switch in A/A. Y108 → Y353			
JR-50-1401	<b>LIGHTS &amp; MISCELLANEOUS:</b> Spotlight, Blue Eye, 400,000 CP, Std Model only	1	WC	JR-50-1401
JR-60-2003	Timer, Momentary Sw,15 Minute, Constant Hot, Std Model onl	1	WC	JR-60-2003
	Y200 Y202 Mount switch on C/S wall at the head of the squad bench in the standard location. Switch to allow activation and also deactivation. Wire to fluorescent lights in liner.			Y351 to ECL
JR-70-0005	Light, in Circuit Board Area, for Electrical Troubleshooting - with integral switch; wired constant hot.	1	WC	JR-70-0005
LP-10-0500	<b>SUCTION PUMPS:</b> Suction Pump, Standard	1	WC	2300363
	W501 *Install suction pump on fixed shelf in compartment "E" below inverter, as close to wall #1 & wall #2 as possible. Install an Expanded metal cover to be as minimal as possible and still allow for protection and ventilation.*			
LU-10-3000	<b>INSULATION:</b> Insulation, Severe Weather Package one inch block foam in modular floor, walls and roof covered with reflective insulation.	1	WC	1000416

PART NO	DESCRIPTION	QTY	ID	REF. NO
	All edges of reflective insulation on walls and ceiling must be taped, std.			
N0-10-9101	<b>GRABRAIL(S)/GRABHANDLES:</b> Grabrail,Overhead,117",Handicap Style,Yellow,Std Model Only  Yellow Powder Coat Finish	1	WC	N0-10-9101
N0-10-9703	Handrails,(3),Custom "L" Shape,1",Yellow,Entry Drs,Std Mode  - Yellow Powder Coat Finish Stainless Steel Grab Handles	1	WC	N0-10-9703
N0-11-15SR	Grabhandle, 15 Inch, Yellow, 90 Degree (SO-GRAB-353032)Ref(SO-GRAB-321002)SR #215106F  Yellow Powder Coat Finish  Locate on curbside wall at head of squad bench.	1		
	<b>O2 CYLINDER RACKS:</b>	1	WC	
NA-10-1500	Rack, Versatile, "M" Or "H", ILOS	1	WC	2309971
	<b>STRAPS, O2:</b>	1	WC	
NA-30-2800	Strap, Cylinder, Rachet Style, (3) Each, ILOS	1	WC	2350804
	<b>O2 SYSTEMS COMPLETE:</b>	1	WC	
NF-50-1001	O2 Sys,LF,(2)O2/(1)Vac Port A/A,(1)O2 C/S, Std Model only Ohio style (2) oxygen outlets in action area (1) vacuum port in action area (1) oxygen outlet curbside wall at head of squad bench	1	WC	NF-50-1001
	<b>O2 SYSTEM PARTS:</b>	1	WC	
P0-59-7500	Regulator, Oxygen Cylinder, Preset 50PSI, Installed	1	WC	P0-59-7500
	<b>SUCTION CONTAINERS:</b>	1	WC	
PA-10-7001	Regulator/Holder, SSCOR 22000, Std Model only	1	WC	PA-10-7001
	<b>COT MOUNTS:</b>	1	WC	
PF-26-3100	Cot Mount, Stryker #6377, Single Position !!!!! Stryker Cot: MX - Pro Eng Same as 333075  Install (1) set of Cot Mount plate in the center position only !!!!! **Measurements are provided.**  Measurement to be approx 85" from the closed rear doors to the center of the most forward antler plate.  *Offset from base wall to be standard for center position.	1	WC	2411700
PF-26-31SR	Hook, Stryker Cot, 6060-36-18 (SO-HOOK-353032)Ref(SO-HOOK-323049)SR215106F WC part #AS81339  Note: Please through bolt rear safety hook. Note: See drawing/PDF in wo folder on server for placement of rear safety hook.	1		

PART NO	DESCRIPTION	QTY	ID	REF. NO
PP-20-1000	IV Hook, Perko, (1), W/Straps, STD  curbside and streetside standard locations	1 2	WC WC	2350508
PP-20-1500	IV Hook, Perko, (1), W/Straps, IATS  Locate Curbside and Streetside (IATS) locations.	2	WC	2351508
PU-10-1000	<b>KICK PLATES STAINLESS STEEL, ON CABINETS:</b> Kickplate, Brushed S/S, Face Of Squad Bench	1	WC	
PU-10-2000	Kickplate, Brushed S/S, On Face of Basewall	1	WC	2410410
R5-10-8000	<b>MISC. INSTALLED HARDWARE/EQUIPMENT:</b> Wrench, Oxygen, Tethered, Mods  in oxygen compt.	1	WC	2410412
R5-11-2010	Net, Squad Saver, at Head of Squad Bench, 2" Straps, Black  attached to floor at head of squad bench, CS wall and overhead liner with removable latches.	1	WC	2350810 R5-10-2010
RA-10-2000	<b>DRI-DECK/TURTLETILE:</b> Dri Deck, Black, All compt. Floors and Shelves	1 1	WC WC	
RK-13-3500	<b>STRAPS, INSTALLED:</b> Seat Belt, Across Compt Opening (with #HA24658 zinc footman loops)  (1) Install in compt "B2" across compt opening, install 30 inches down from ceiling and 17" from wall #2  (1) Install in compt "D" across vertical dividers, mid-height between the floor and the first shelf. Note. >> Install seat belt on wall #1 and wall #3, (8 inches) from wall #2.  (1) Install on compt. "D" door to hold in Stair Chair. Locate strap across the stair chair location approx 30 inches from the bottom of the door. Note: dimensions FOLDED of the Stryker MD# 6252 stair chair are: Height - 37.5 inches, Width - 20.5 inches, Depth - 8 inches  **See Reference photos submitted with the order file.	1 3	WC WC	1101018 2350027
RU-05-0000	<b>MODULE INTERIOR</b> Interior Cabinets - Duralite Construction, Wood  Attention: use exterior grade flexible caulk for interior cabts to avoid cracking !!	1 1	WC WC	RU-05-0000
RU-05-01SR	Loctite on all interior threaded fasteners (SO-LT-353031) SR# 213354F	1		


PART NO	DESCRIPTION	QTY	ID	REF. NO
	Apply Loctite or equal thread lock, to all interior machine thread screws to include all grab handles overhead and on doors and all hardware on cabinets (hinges and all restocking hardware)			
	<b>MICA COLOR:</b>	1	WC	
RU-10-3010	Mica, Gray Pampas #4168-01, High Gloss	1	WC	RU-10-3010
	<b>SUB FLOORING:</b>	1	WC	
RZ-99-0025	Flooring, Wood, Type 3/7/9, Std	1	WC	RZ-99-0025
	<b>FLOORING:</b>	1	WC	
T0-11-8203	Flooring, Lonplate II, Gunmetal Grey #424, Std Model only	1	WC	T0-11-8203
	<b>SEAMLESS ATTENDANT SEATS:</b>	1	WC	
T5-11-2D00	Seat, EVS, W/ Childseat, Dusty Jade Green, Swivel Base, ILC	1	WC	1926024
	Vacuum formed EVS1850S w/ 2 point seatbelt. With EVS swivel seatbase.			
	*With the seat moved all the way toward the cot position (forward if you are sitting in the seat facing the rear doors) Install the EVS seat 99" from the closed rear doors as measured to the center of the slide adjustment bar on the front of the seat. Offset from the Basewall...to allow for the seat to swivel from facing the back doors to facing the front of the vehicle.*			
	*See photo submitted with the order file*			
	<b>SEAMLESS UPHOLSTERY:</b>	1	WC	
T5-25-0400	Upholstery, Seamless, Dusty Jade Green, ILOS	1	WC	1926004
	<b>BULKHEAD CABINET:</b>	1	WC	
TF-51-9500	Cabinet, Bulkhead, 10" Ext , 72Hr, Type3	1	WC	C971817
	NO broom closet behind attendant seat.			
TF-51-95SR	Bulkhead Door, Upper, Double Wood, Side Hinged, ILOS (SO-BHDR-353032)Ref(SO-BHDR-303153)SR #215106F	1		
	With piano hinges.			
	Note: Electrical access doors to have "M1" latches, same as previous unit.			
	<b>O2 DOOR:</b>	1	WC	
TK-10-5010	Lexan Flaps, Hinged Left Side, Fixed Right Side, Std	1	WC	TK-10-5010
	<b>GLOVE HOLDERS:</b>	1	WC	
TN-20-1001	Cabinet, Glove Box Holder,(3), Above C/S door, Std Model only (narrows cabt "E1" on Pureair 1 w/low headroom)	1	WC	TN-20-1001
	<b>BULKHEAD SLIDING DOOR:</b>	1	WC	
TP-10-6000	Door, Sliding, Bulkhead, STD with sliding Lexan window	1	WC	2100178
	<b>RIGHT FRONT CABINET:</b>	1	WC	
TU-19-5100	Cabt, RF, Combo AC, Virginia Style (SO-RF-353032) Ref (SO-RF-313070) same as 333075  Upper Cabinet "RF/E1" - Inside Access Only. Upper Cabinet "RF/E2" - Inside/Outside Access.	1	WC	C990563


PART NO	DESCRIPTION	QTY	ID	REF. NO
	<p>Lower Cabinet "RF/E3" - Inside/Outside Access.</p> <p>Upper "RF" cabinet to be divided into two (2) equally spaced cabinets - "RF/E1" and "RF/E2"</p> <p>"RF/E1" &amp; "RF/E2" shall have separate access doors. Minimal partition and style between left and right side.</p> <p>"RF/E1 cabinet door (upper left cabinet) Two (2) drug boxes shall be stored length wise (front to rear) in "RF/E1" cabinet. The boxes shall be stored on top of each other separated by a shelf. The boxes are Flambeau PM 2072. The dimensions of drug boxes are (10.75" high x 19" long x 11.5" wide). Hinge the "E1" door on the right side, with hold open cable as long as possible.</p> <p>"RF/E2 cabinet door (upper right cabinet). Hinged on right side with chrome southco latch.</p> <p>"RF/E3" cabinet doors (lower cabinet). Hinged on right and left sides with leading door on right side.</p> <p>Note: Shelves w/shelf track shall be in cabinets - "RF/E1" and "RF/E2".</p>			
TU-20-8200	<p><b>RF CABINET DOORS W/COMBO A/C</b></p> <p>Door, Virginia Style, Split, "E1" &amp; "E2"</p> <p>"RF/E1" and "RF/E2" Cabinets - Both hinged on right side.</p>	1	WC	RFDR5
TU-21-SR01	<p>Doors, 1/2" Clear Lexan, Dual, RF Cabinet, ILOS (SO-RFDR-353032)Ref (SO-RFDR-313070)SR #215106F</p> <p>One (1) set for cabinet "E3" (both independent)</p>	1		
V0-10-2000	<p><b>HARDWARE,LATCHES AND LOCKS:</b></p> <p>Latch, Plunger Roller, Installed</p> <p>For "RF" Cabinet Doors.</p>	1	WC	2109714
V0-12-6200	<p>Hasp, Padlock, Marine, Installed, IATS (SO-HASP-353032) Ref (SO-HASP-301186)</p> <p>For "RF/E1" Cabinet Door.</p>	1		
V0-14-2075	<p>Latch, Southco M1-64-8, 2IN S/S, Non-Locking, (1), ILOS (self latching) Flush Mount. Non-Locking</p> <p>One (1) For "RF/E1" Cabinet Door. One (1) For "RF/E2" Cabinet Door.</p>	2	WC	V0-14-2075
V0-14-2076	<p>Latch, Southco M1-63-8, 2IN S/S, Non-Locking, (1), IATS</p> <p>Two (2) for RF cabinet "E3" doors (1/2" lexan)</p>	2	WC	SO-LOCK2-351015

PART NO	DESCRIPTION	QTY	ID	REF. NO
V0-14-2077	Latch, Southco M1-64-8, 2IN S/S, Non-Locking, (1), IATS (self latching) Flush Mount. Non-Locking  One (1) For Upper Bulkhead Lead Door. One (1) For Electrical Cabinet Door.	2	WC	V0-14-2077
V5-20-3400	<b>SQUAD BENCH:</b> Squad Bench, Standard, No Bio-Waste See tilt-out waste same as 333075  Aluminum trim at back of bench and CS wall ILOS caulk	1	WC	C940472
V5-20-34SR	Tilt-Out Waste Container At Head of S/B, Facing Aisle (SO-SB-353032) Ref (SO-SB-323049) SR #215106F  Per Eng design WC std container	1		
V5-64-1000	<b>LID 170 INCH BOX:</b> Lid, Squad Bench, Split 50/50	1	WC	170SBLD2
V5-80-1005	<b>HARDWARE:</b> Handle, Trimark, Squad Bench, STD	1	WC	V5-80-1005
V5-82-0500	Strut, Gas, 30 LB, Installed, ILOS	2	WC	2109715
VA-40-3400	<b>CURBSIDE SPLINT CABINET 160-170 inch Module</b> CS Splint,7in Deep,2 Sect, Clear Lexan Sliders Note: Will not meet KKK requirement for space between cushion on squad bench and bottom of splint cabt.	1	WC	C951217A
VF-20-4700	<b>CURBSIDE REAR CABINET:</b> Cabinet,CurbSide Rear,Standard,170",72"HR <b>CABINET H, BASEWALL 170", T1/T3</b>	1	WC	2107001
VH-H2-0300	H, Cabinet, Inside/ Outside Access <b>CABINET I, BASEWALL 170", T1/T3</b>	1	WC	VH-H2-0300
VH-I2-0400	I, Cabinet, Deleted <b>CABINET M, BASEWALL 160-170", T3</b>	1	WC	VH-I2-0400
VH-M1-0300	M, Cabinet, deleted <b>CABINET O, BASEWALL 160-170", T1/T3</b>	1	WC	VH-M1-0300
	DELETED			
VH-P2-0200	P, Cabinet, full height  CPR side seat moved toward rear of unit in place of cabt "O"	1	WC	VH-P2-0200
VH-Q0-SR00	<b>CABINET CPR, BASEWALL, 149"-170", T1 &amp; T3</b> CPR, side seat, Custom (SO-CPR-353032) Ref (SO-CPR-303098) SR #215106F  Move CPR side seat rearward in place of cabinet "O", telemetry work area & cabinet "K".  Eng: Dupl previous unit #333075	1	WC	VH-Q0-SR00
VH-Q0-SR01	CPR Side Seat, Custom Storage for O2 Cylinders, "D" (SO-CYL-353032) Ref (SO-CPR-303211)SR #215106F  Custom CPR Side Seat, Under Seat Storage - Increase depth/height and width (Wall #1 to Wall #3) to hold (3) "D" size	1		


PART NO	DESCRIPTION	QTY	ID	REF. NO
	<p>oxygen cylinders. Cylinders to rest in a "W" trough. The cylinders are to be secured with (2) quick release straps across them, equally spaced at 1/4 and 3/4 the length (height) of the tank. Cylinder dimensions are 21" height x 4.5" diameter. This storage area shall be manufactured so that you may take the cylinder in and out without turning it sideways.</p> <p>If CPR seat needs to be widened slightly, take space from cabt "P".</p> <p>See Drawing Mark Up Submitted With The Order File</p> <p>Eng: Dupl previous unit #333075</p>			
VI-J2-0300	<p><b>CABINET J, STREETSIDE UPPER 170", T1/T3</b> J, Cabinet, 72" H.R.</p> <p><b>CABINET K, STREETSIDE UPPER 160-170", T1/T3</b></p> <p>DELETED</p>	1 1 1	WC WC WC	VI-J2-0300
VI-K1-0200	K, Cabinet, Deleted	1		VI-K1-0200
VI-L1-0250	<p><b>CABINET L, STREETSIDE UPPER 160-170", T1/T3</b> L, Cabinet, L1, L2, L3 Combine 2 Equal Cabinets, 3/4F, 72" H</p> <p>Extended AA L1 &amp; L2 to occupy from the CPR seat forward.</p>	1	WC	VI-L1-0250
VK-15-3000 VK-16-6000	<p><b>BASEWALL DOORS:</b> Cab H, Clear Lexan Slider Cab P, Clear Lexan Slider</p> <p><b>TELEMETRY AREA WORK SURFACE</b></p> <p>DELETED, CPR seat moved toward rear of vehicle</p>	1 1 1 1	WC WC WC WC	2100001L 2100027L
VP-10-5918 VP-10-59SR	<p>A/A Tray, 160"-170", W/O Bio,Poured, Gray, Extended Tilt-Out Waste Container, Beneath A/A Work Surface, Facing Aisle (SO-TILT-353031)Ref (SO-AA-323049) SR #215106F</p> <p>Per Eng design WC std container</p> <p>Note: This is to be a waste container not a sharps. This is different from Princess Anne Unit. They want both tilt outs to be waste containers.</p>	1 1	WC	VP-10-5918
VP-15-7000	<p><b>ACTION AREA 165/170 INCH BOX:</b> Cab,A/A, Extended, Without Biowaste, 160-170" In ILOS (SO-AA-353031) Ref (SO-AA-333075)</p>	1	WC	
VU-14-4000	<p>Shelf, Interior, Adj, in Cabinet: (1) each in cabts "J", "L1", "L2", &amp; "P"</p>	4	WC	2101503
VU-14-4001	<p>Shelf Track, Cabt, #HA24663, Upgrade, Std Model only ("C" channel/exterior shelf style track)</p> <p>cabts "J", "L1", "L2", &amp; "P"</p>	4	WC	VU-14-4001

PART NO	DESCRIPTION	QTY	ID	REF. NO
VU-15-1000	Cab "J" Clear Lexan Slider	1	WC	2100060
VU-15-7000	Cab "L1" Clear Lexan Slider	1	WC	2100066L
VU-16-2000	Cab "L2" Clear Lexan Slider	1	WC	2100071L
	<b>SHIP LOOSE ITEMS:</b>	<b>1</b>	<b>WC</b>	
YY-00-0001	Bracket, Spare Tire Mounting, Ship Loose	1	WC	YY-00-0001
YY-10-6400	Ship Loose, W.C. Standard Items	1	WC	2609998
YY-13-2000	Fire Extinguisher. (1) 5 LB. ABC, IATS, Ship Loose	1	WC	2600601
YY-13-3000	Spare Tire, OEM only, Ship Loose	1	WC	1000800
YY-13-3001	Door Panel Screws (100) (SO-HDW-353031)	1		
	<b>WARRANTIES</b>	<b>1</b>	<b>WC</b>	
ZZ-10-0101	Warranty, Conversion, 12 Month	1	WC	0589998-14
ZZ-10-0601	Warranty, Limited Electrical	1	WC	0589993-14
ZZ-10-0701	Warranty, Cabinet Construction	1	WC	ZZ100701-14
ZZ-10-1900	Warranty, Lifetime, Mod Structural	1	WC	0589976
ZZ-10-2000	Warranty, Limited Lifetime, Mod Paint	1	WC	0589980
	<b>EMS SHIPLOOSE</b>	<b>1</b>	<b>WC</b>	
ZZ-ZZ-0200	End Of Order	1	WC	9980002

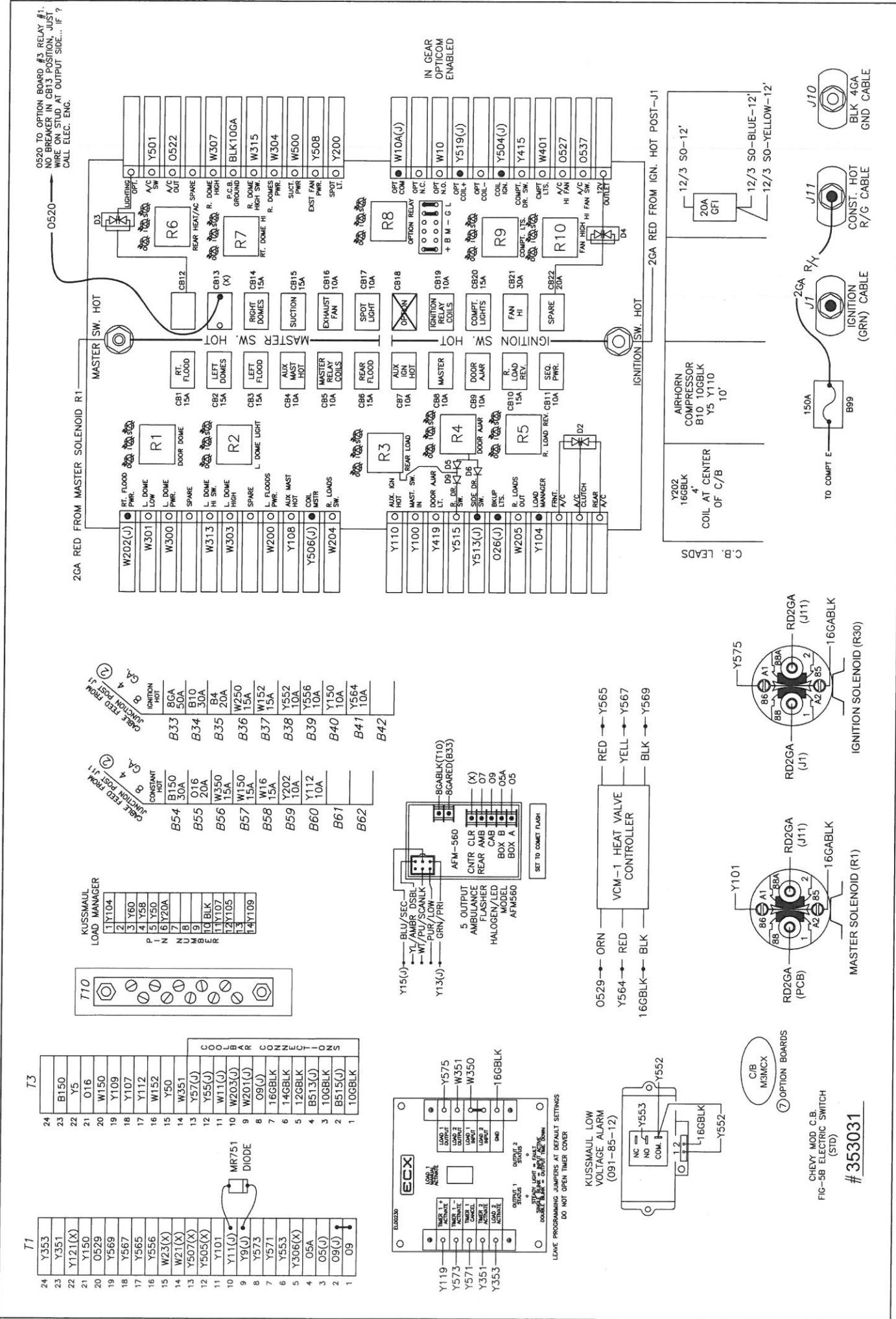


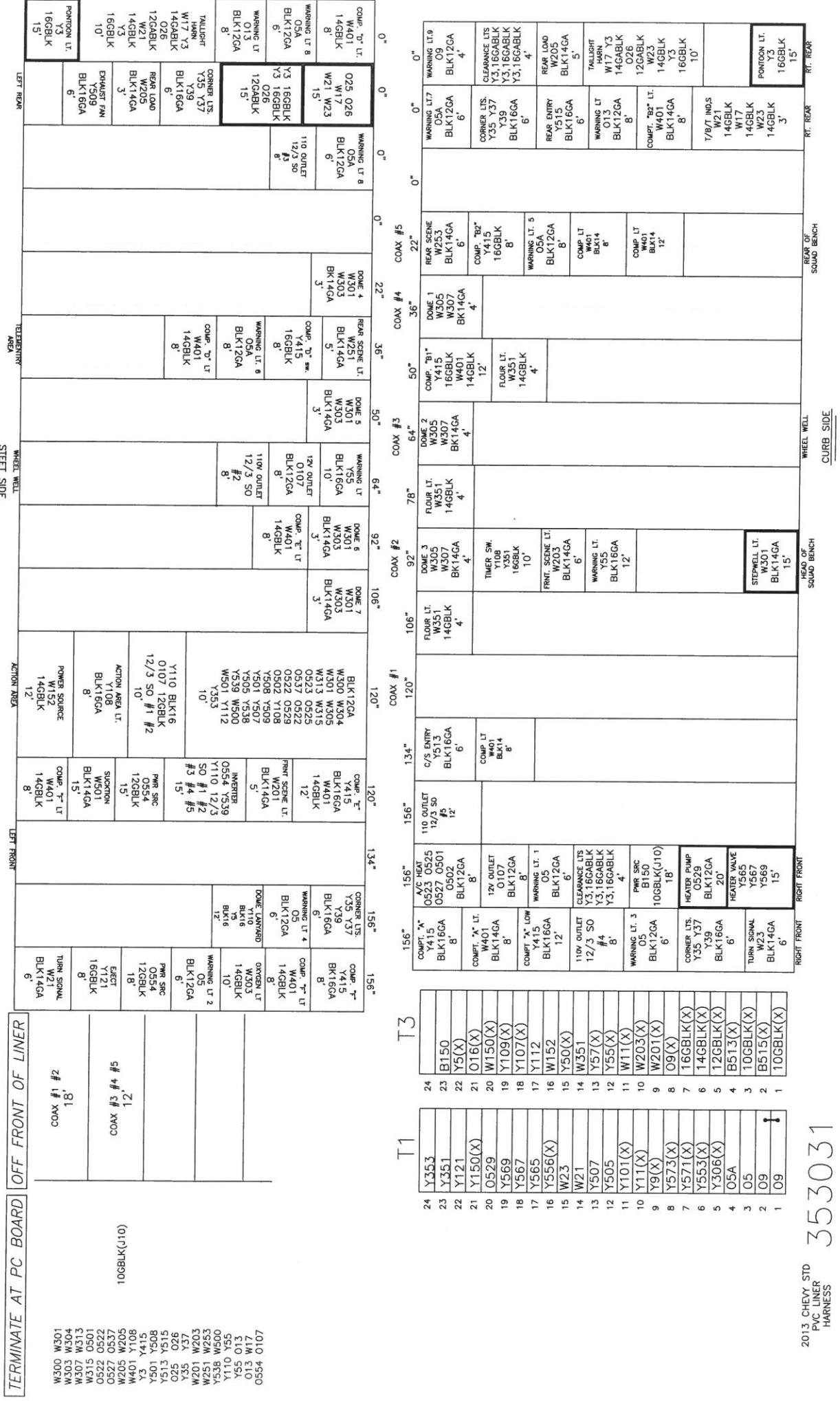

### OPTION BOARD # 1

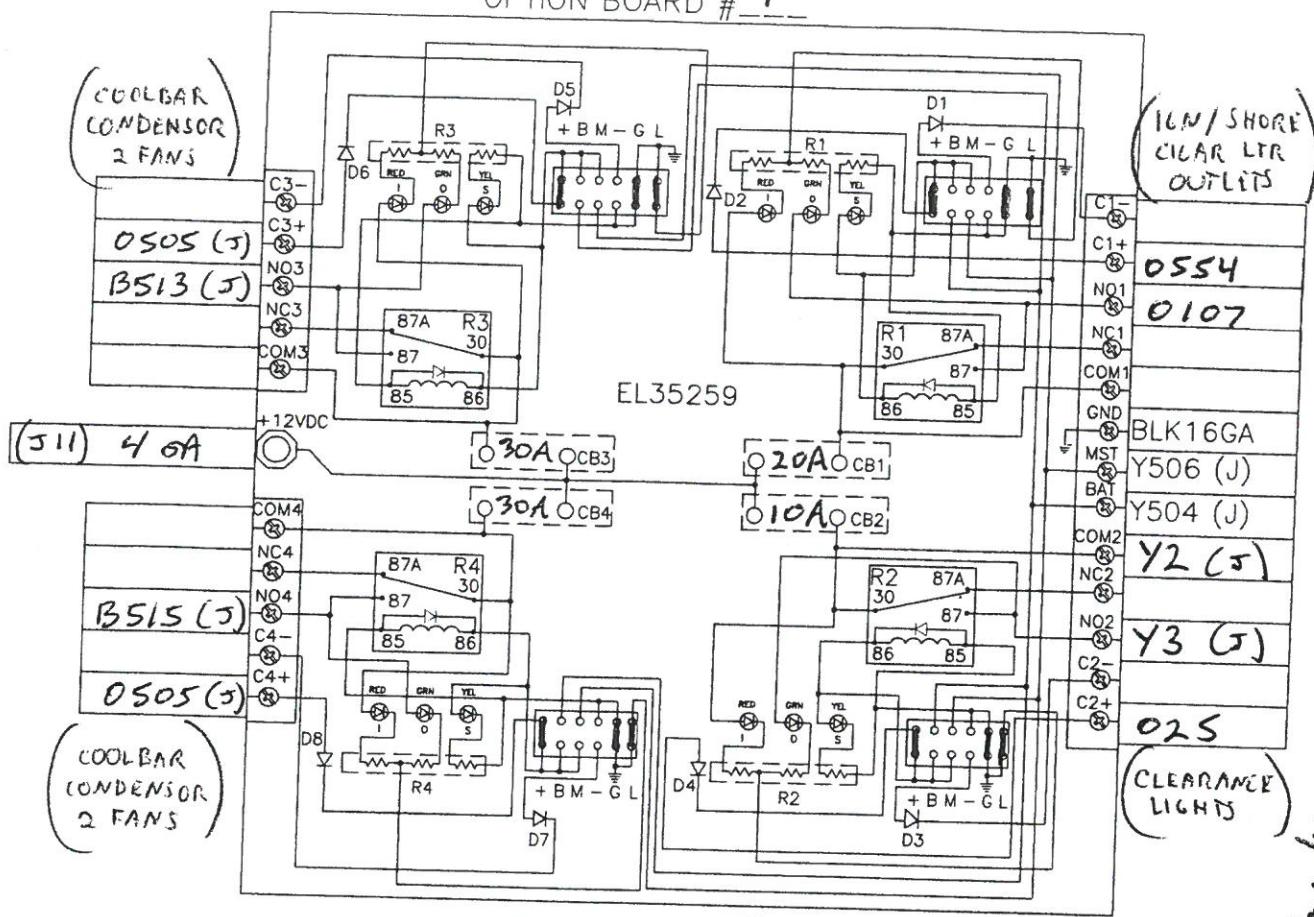



### OPTION BOARD # 2

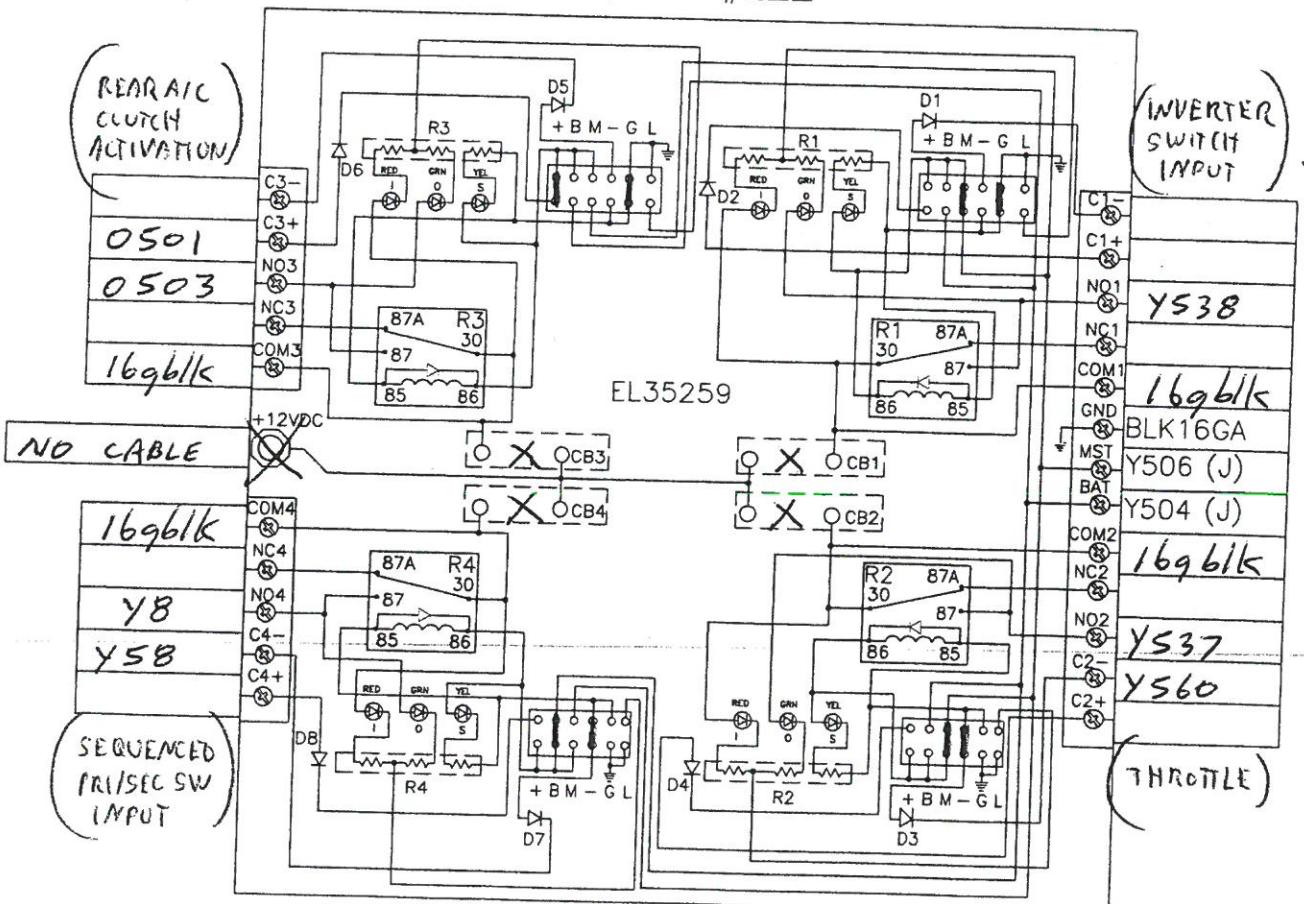

William S. Oberender  
4 / 20 / 15  
Date

William S. Oberender  
Electrical Designer

353031

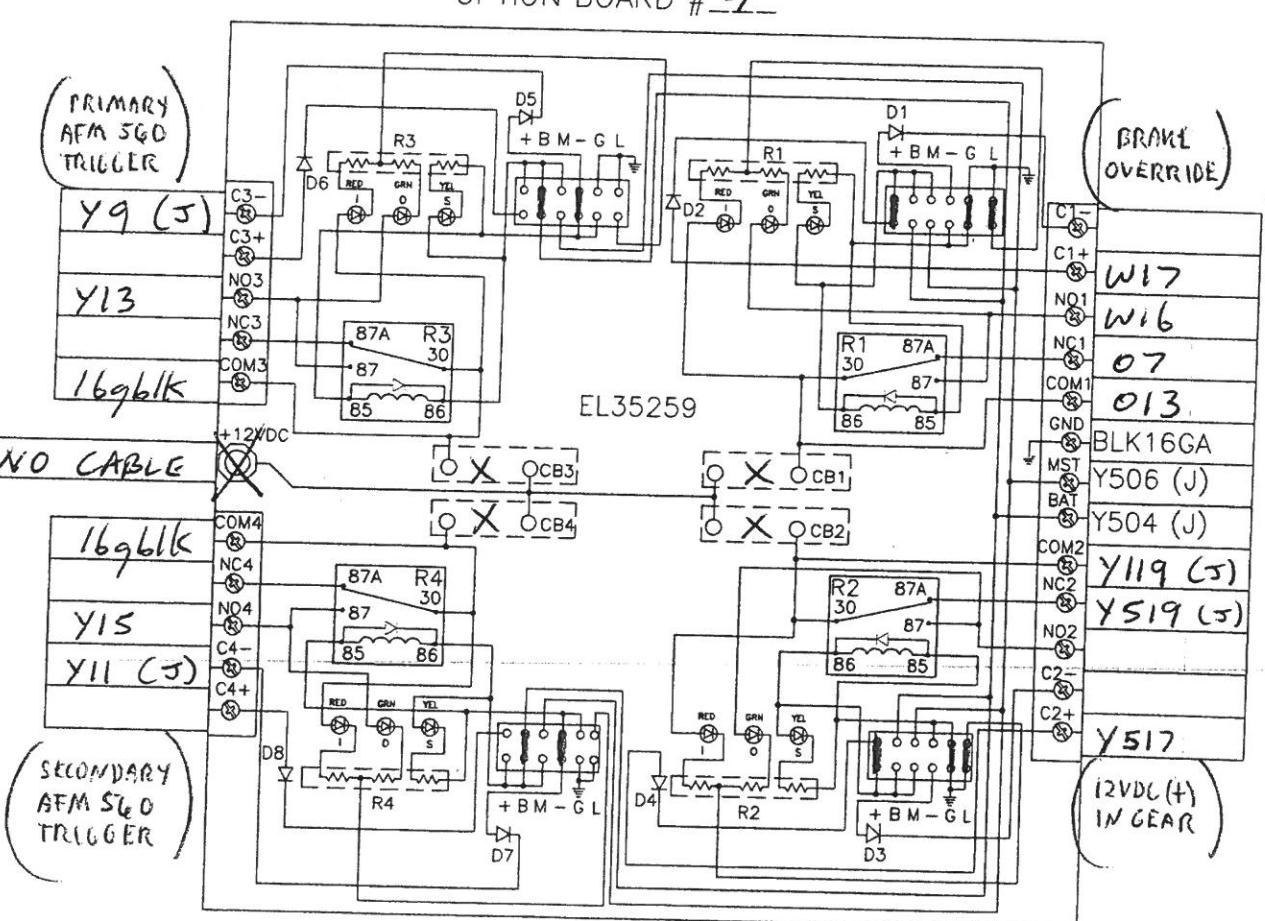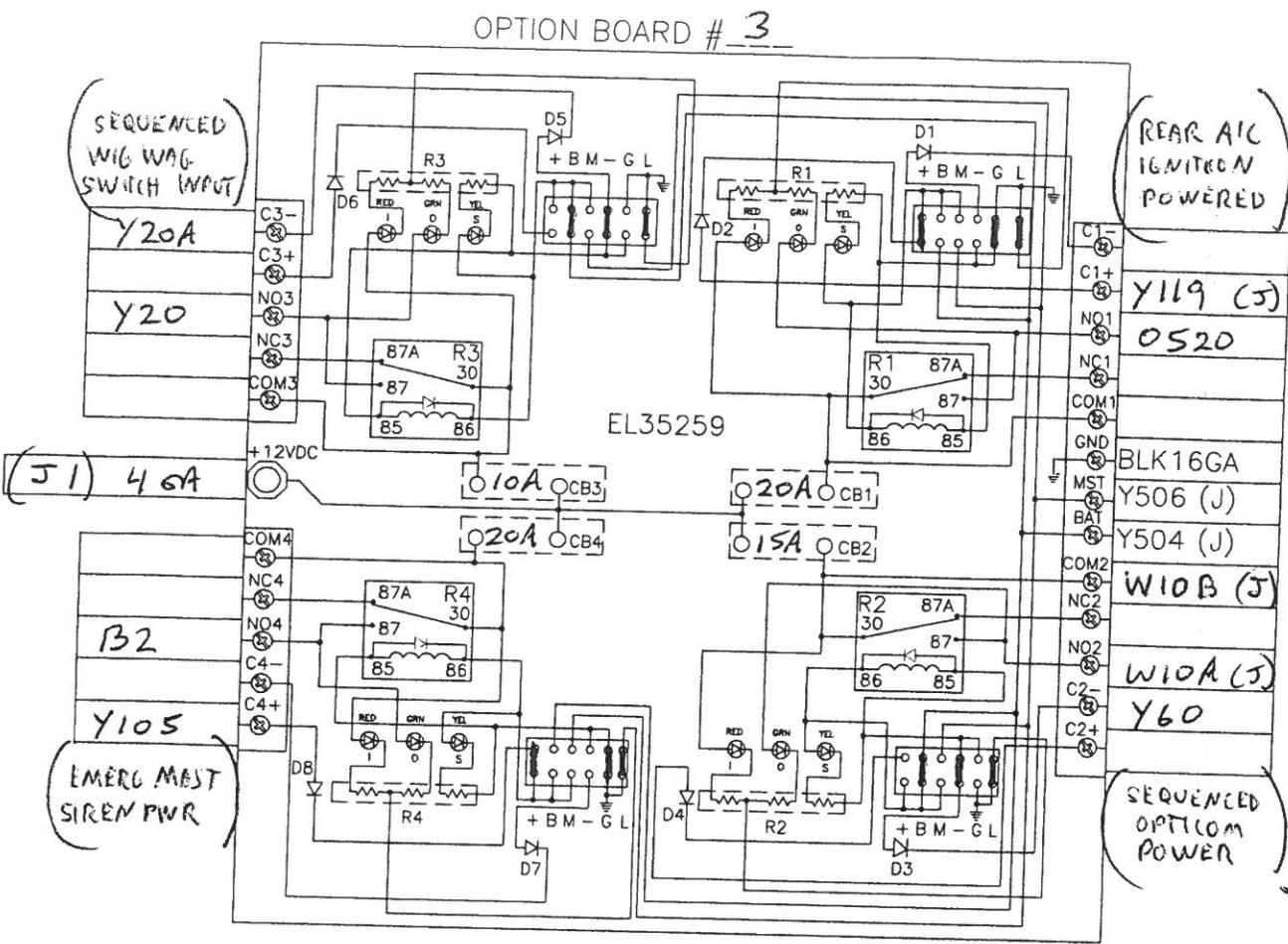


William Oberender  
Electrical Designer

4 / 20 / 15


William L. Oberender

353031

### OPTION BOARD #5

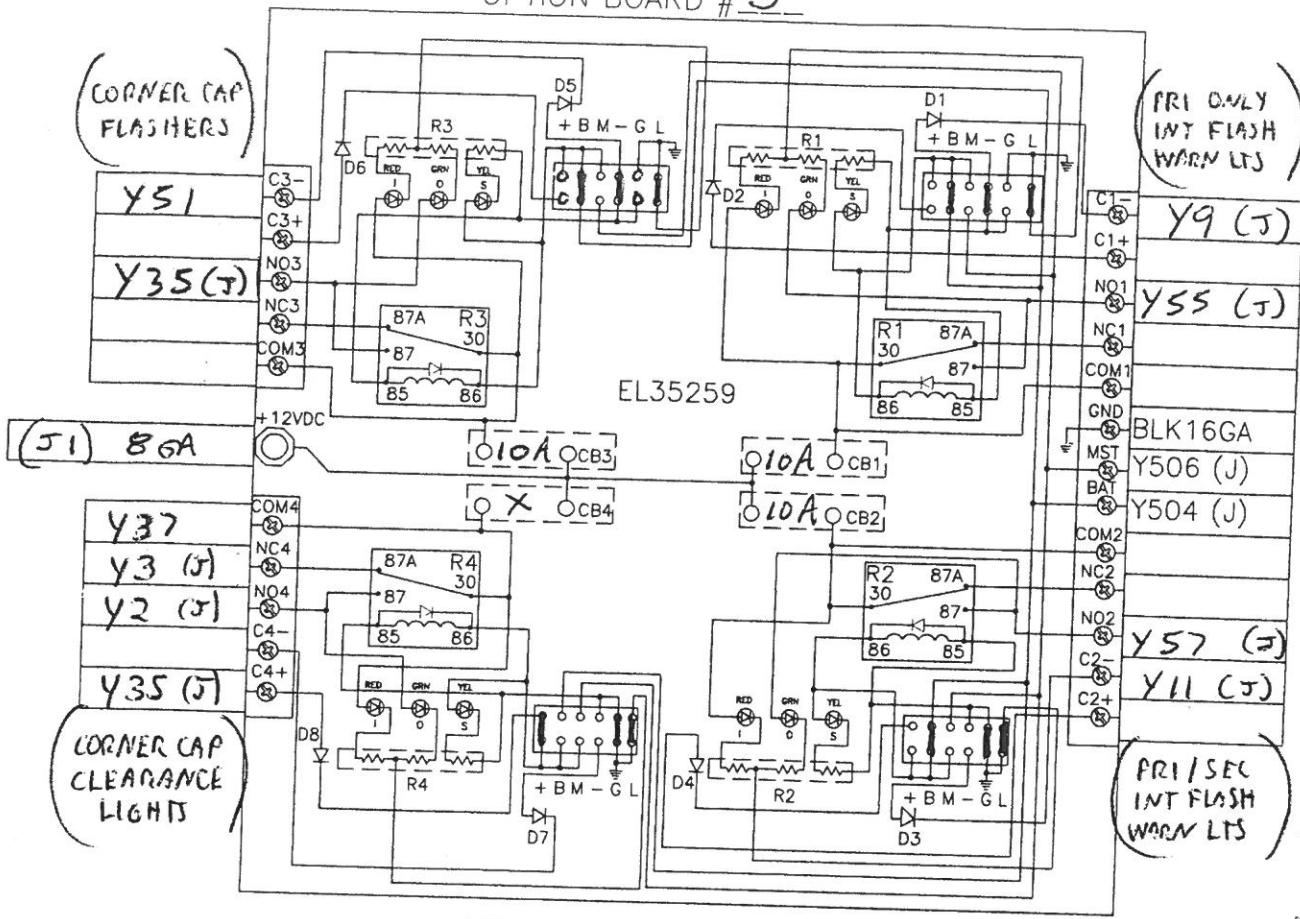

### OPTION BOARD #6

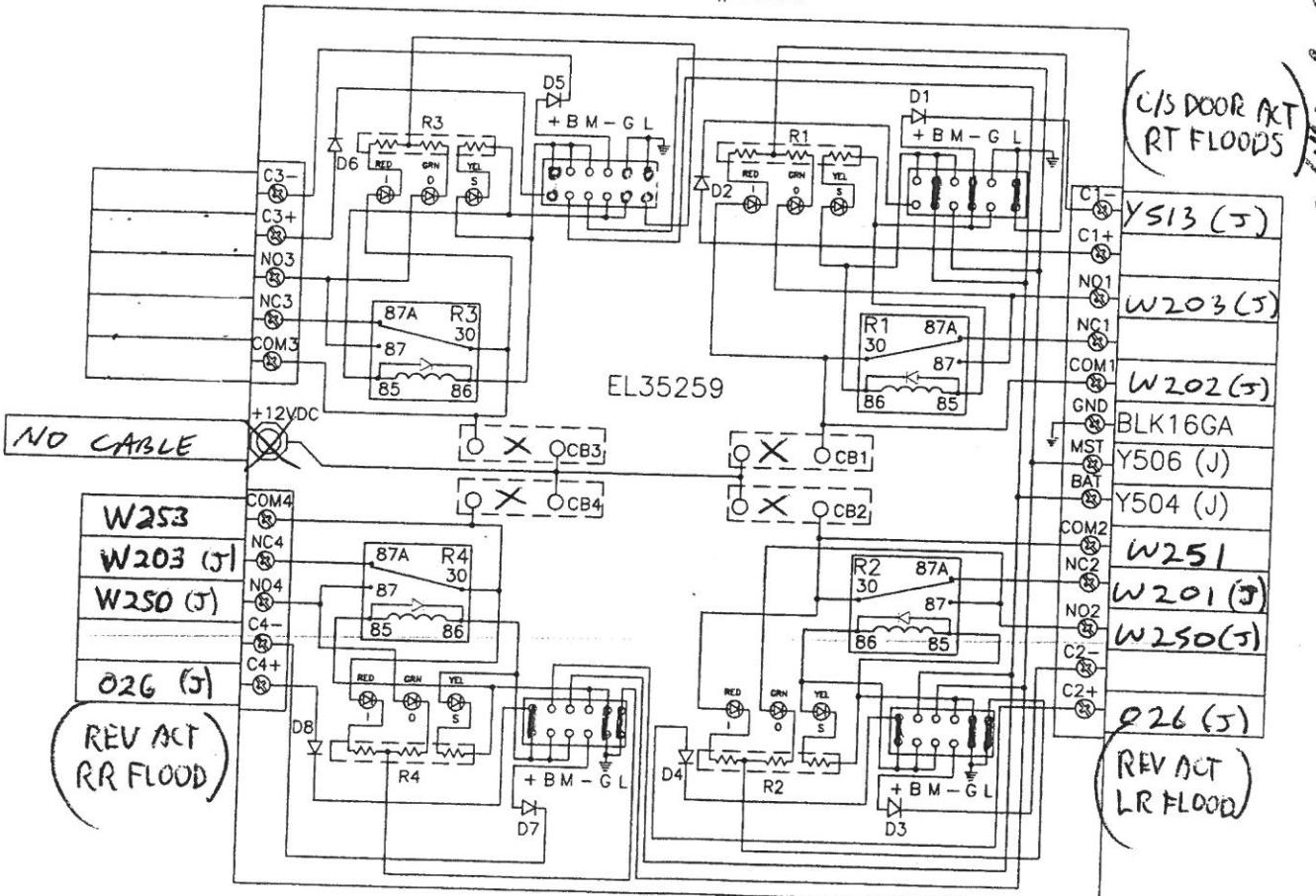


William Oberender  
Electrical Designer

4/20/15

353031


OPTION BOARD # 7

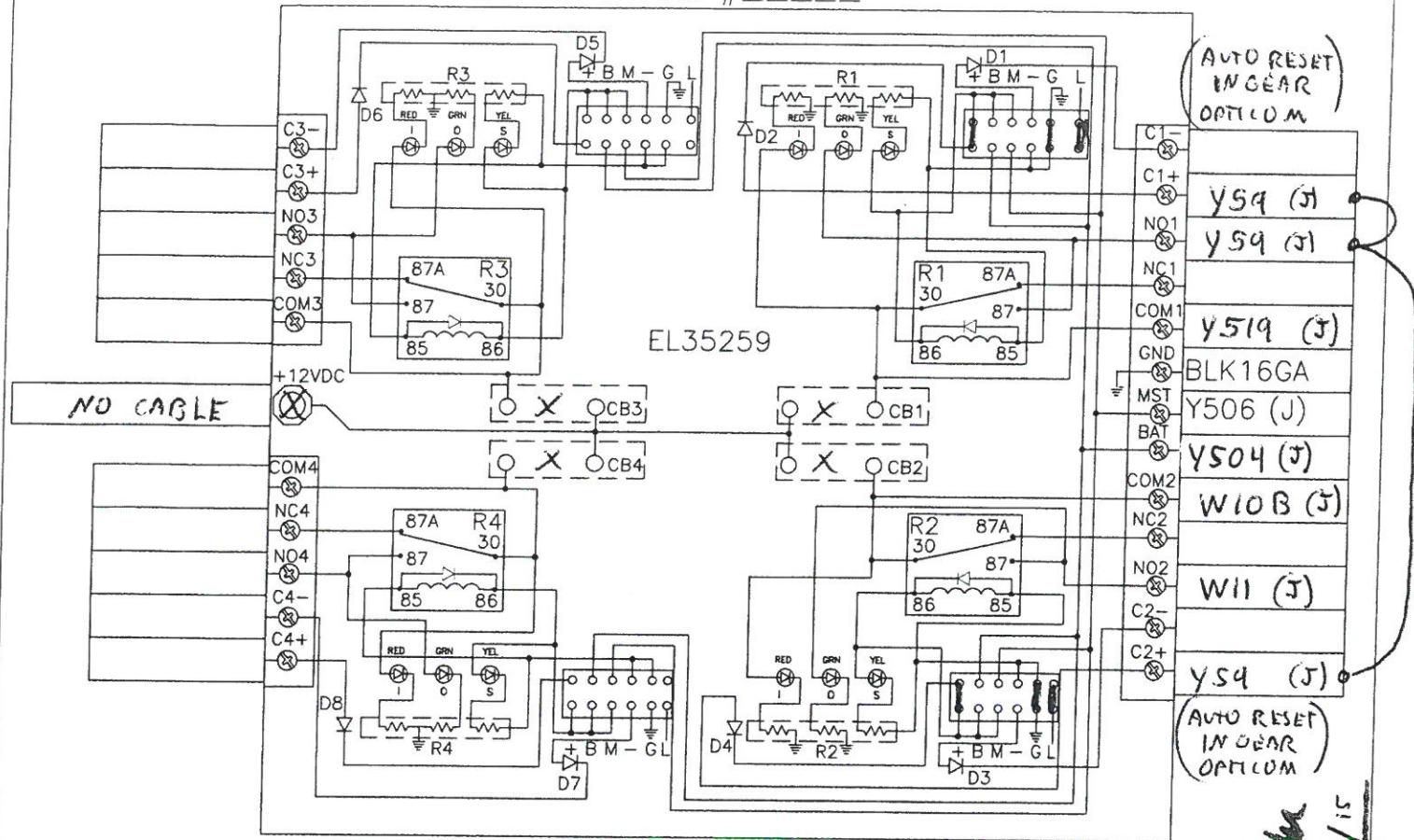



4 / 20 / 15  
Date

William Oberndorfer  
Electrical Designer

353031

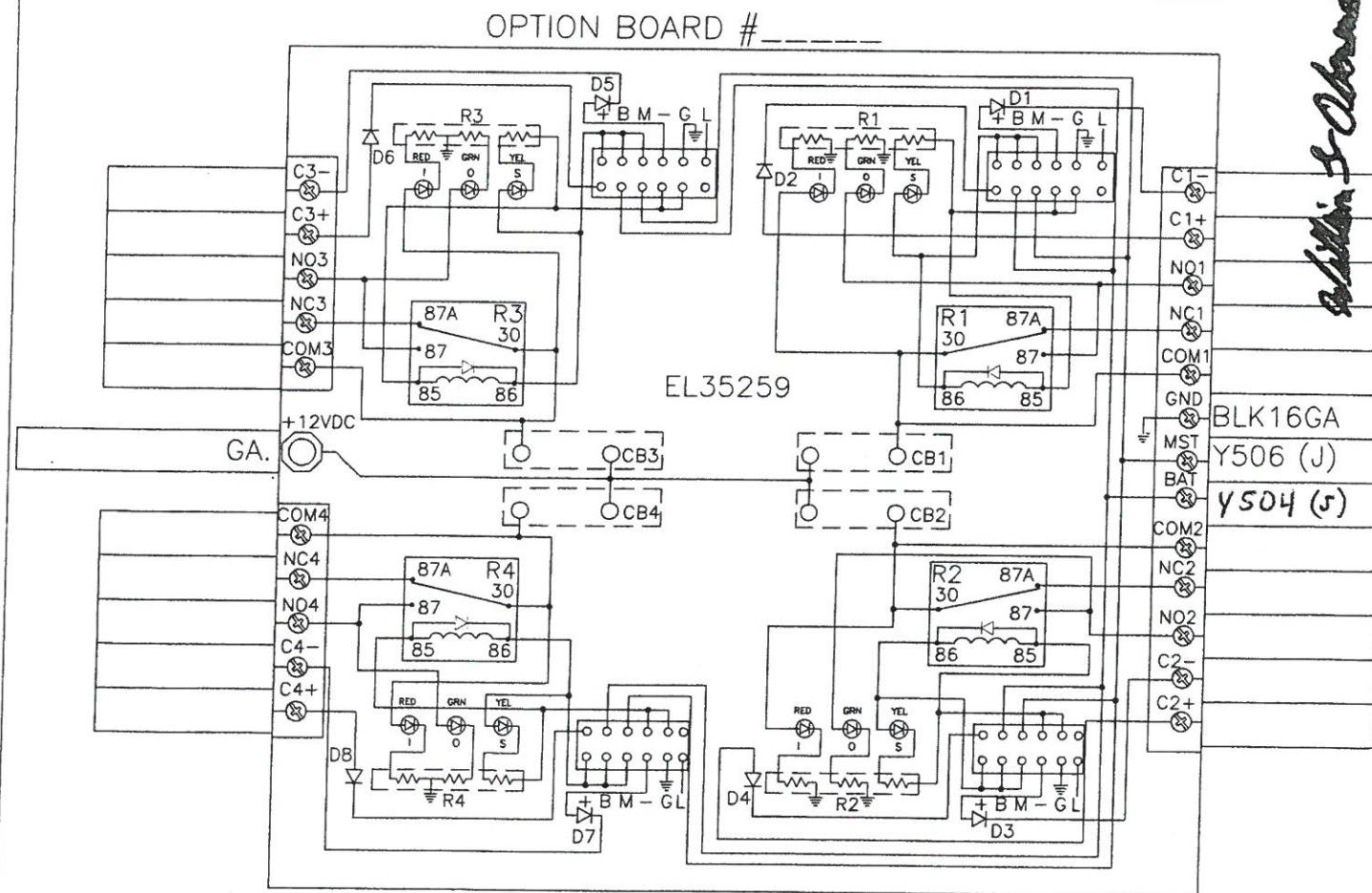

REV. LEVEL	REVISIONS	DESCRIPTION
		5M1

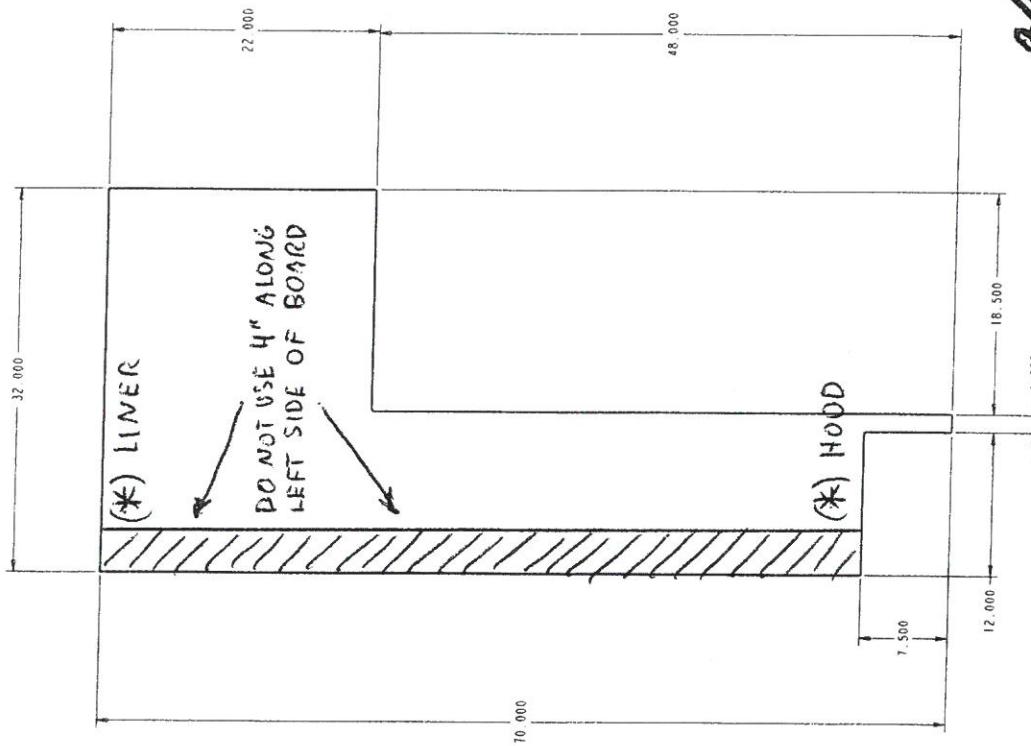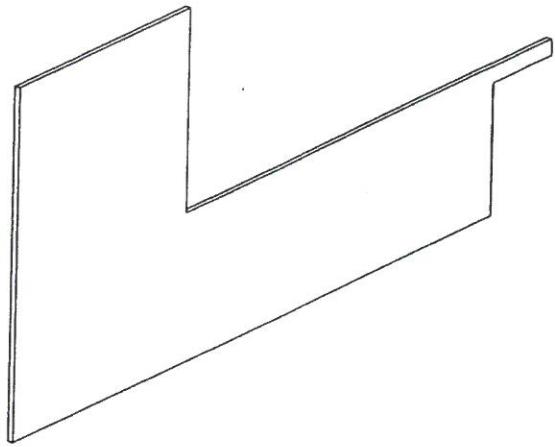

TYPE III / VII / CUTAWAY  
TYPE I FHWT  
WITH 72" HEADROOM

PVC PLASTIC  
P/N: M3MCX


NOTES:  
1. DRAWING CREATED FOR U.S. AMBULANCE PER INC 1786.

ITEM NO.	PL.18226	112 PVC	32.000	32.000	0.864	5.5	13.2
ITEM NO.	PART NO.	DESCRIP/TION	LENGTH	WIDTH	THICK	INCH	INCH
0	1-1111 PVC CIRCUIT BOARD	Wheeled Cart	34.28710				

BY: J. JORGE  
DATE: 17-Jan-14  
TIME: 0 200


353031

PARTS & F. FILENAME: M3MCX

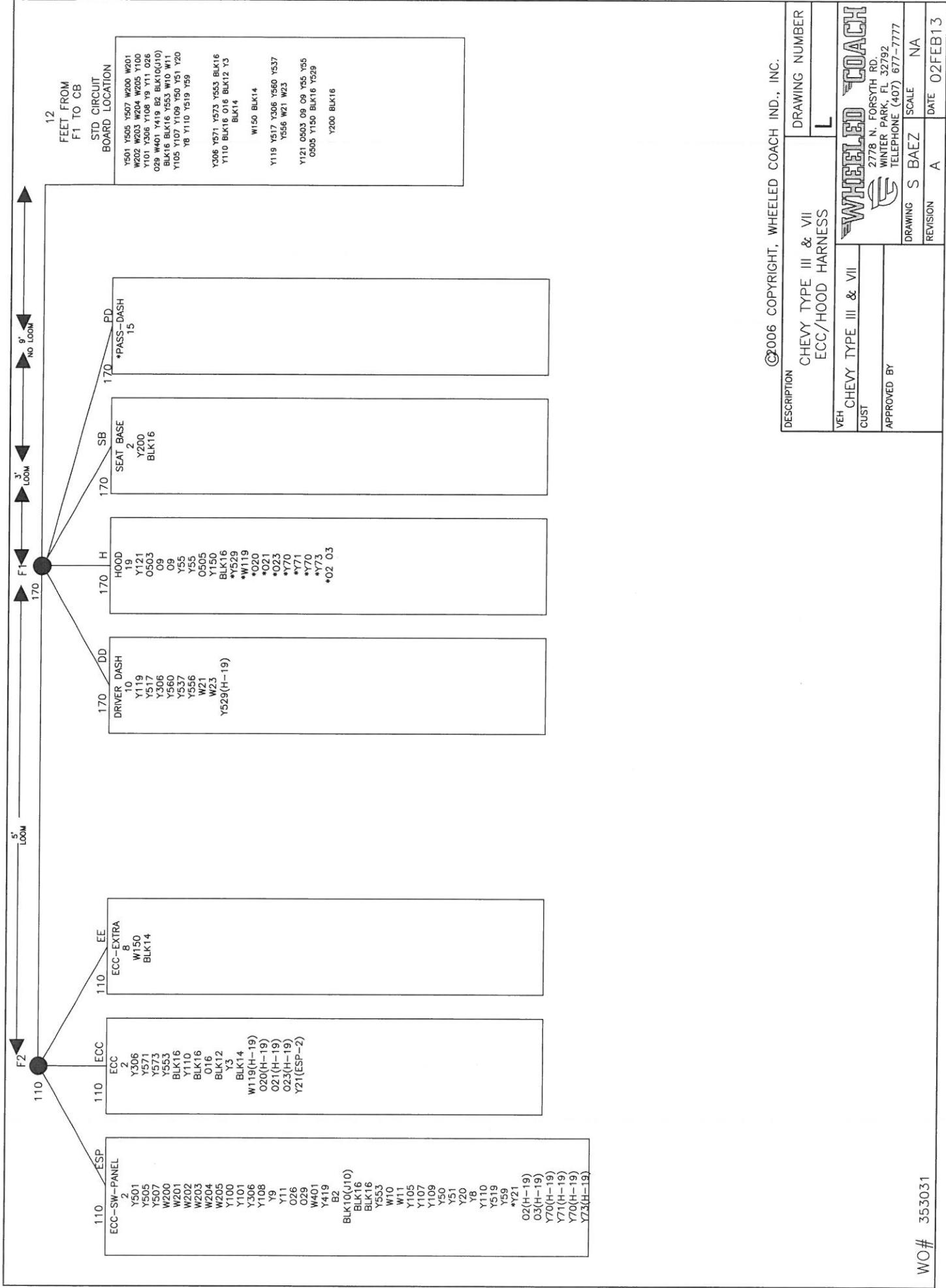

WO# 353031

12

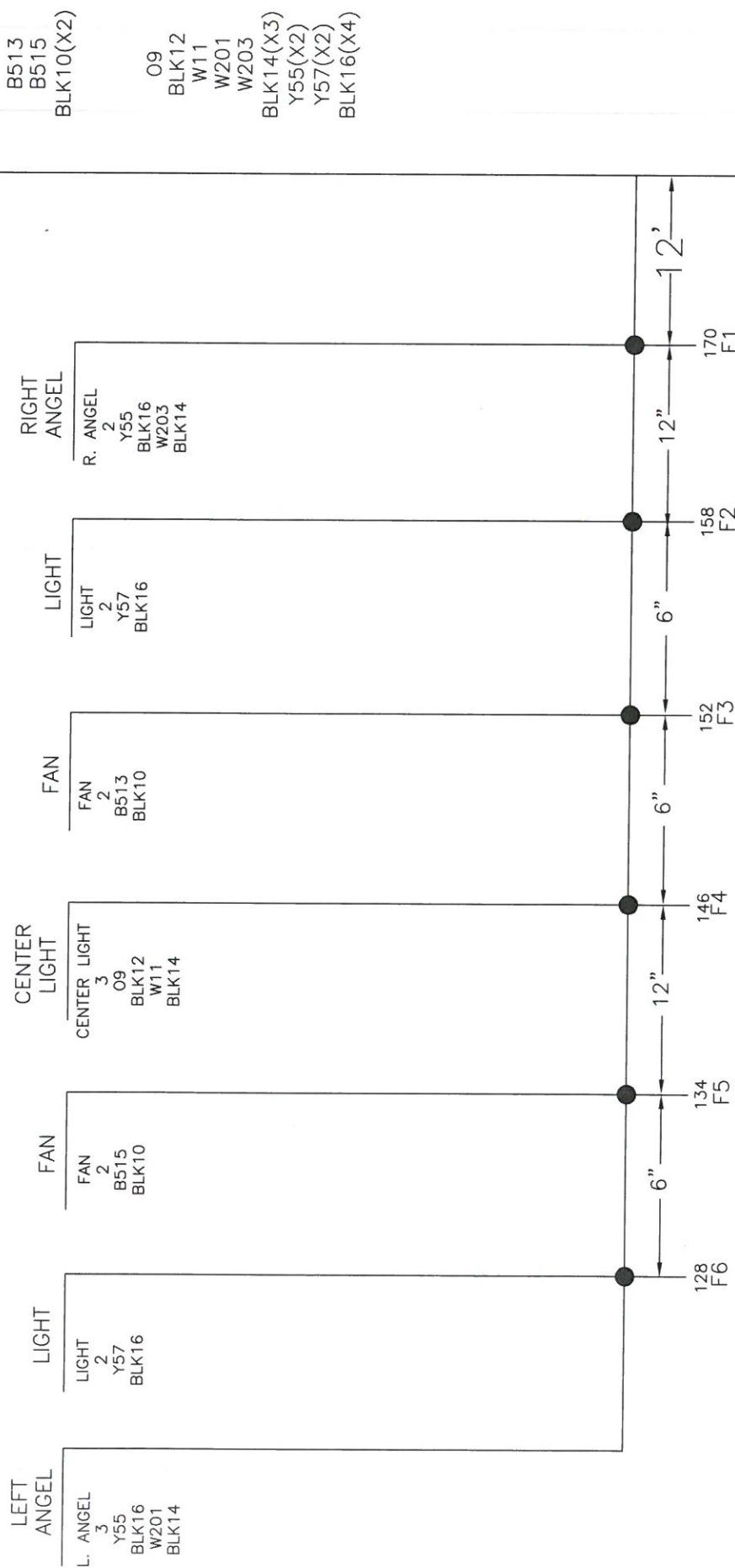FEET FROM  
F1 TO CB

## STANDARD COOLBAR HARNESS

© COPYRIGHT, WHEELED COACH IND., INC.

DRAWING NUMBER  
L 336000

DESCRIPTION COOLBAR HARNESS STD VBH ALL CLUST	WHEELED COACH	DATE 6/27/2013
ADDRESS 21776 FOREST RD. MURFREESBORO, TN 37132 TELEPHONE (401) 572-7777	AMBIANCE TYPE ALL	SCALE N/A
REVISION S	DRAWING BACZ	DATE 6/27/2013
This drawing contains information proprietary to Wheeled Coach Ind. Any reproduction, disclosure or use of the drawing is expressly prohibited without the written permission of Wheeled Coach Ind.		

353031

Chevy Type III and VII Diesel  
1 Batt Under Hood 1 Batt RF

Cable Gauge	Length	Lugs	Color Code	Function
3/0	10"	1/2 - 5/16	Red/Green	J7 to B100
3/0	1'	3/8 - 3/8	Black	Box to Frame Ground
3/0	3' 6"	3/8 - 3/8	Black	Hood Battery (-) Ground
3/0	16"	3/8 - 1/2	Orange	Hood Battery (+) to J7
3/0	5'	3/8 - 3/8	Black	RF Battery Ground to Frame
3/0	5'	3/8 - 3/8	Black	Frame to Frame Ground
3/0	12' 6"	1/2 - 3/8	Blue	J7 to RF Battery
3/0	15'	5/16 - 1/2	Red/Green	B100 to J11
3/0	2'	3/8 - 3/8	Black	Engine to Frame Ground
1/0	1' 6"	3/8 - 3/8	Black	Cab to Frame Ground
1/0	15'	1/4 - 3/8	Black	Liner Ground
2Ga Red	2'	1/2 - 5/16	Red	J7 to OEM Junction
2Ga Red	6'	5/16 - 1/2	Red/Red	Upgrade Alternator to J7
2Ga Red	7' 6"	5/16 - X	Yellow/Red	B99 to Inverter (+)
2Ga Blk	9'	3/8 - X	Black	Frame Ground to Inverter (-)
2Ga Blk	5'	3/8 - 5/16	Black	Frame Ground to J10
2Ga Blk	6'	3/8 - 1/4	Black	Frame Ground to T10
8Ga Blk	11'	3/8 - X	Black	GFI and Case Ground
8Ga Blk	10'	3/8 - X		Both in One Lug
3/0	5"	3/8 - 3/8	Blue	Batt 2 to Batt 3 (+) Jumper
3/0	5"	3/8 - 3/8	Black	Batt 2 to Batt 3 (-) Jumper



William S Oberender

William Oberender  
Electrical Designer

4/29/15  
Date

353031

# LOAD MANAGER / SEQUENCER

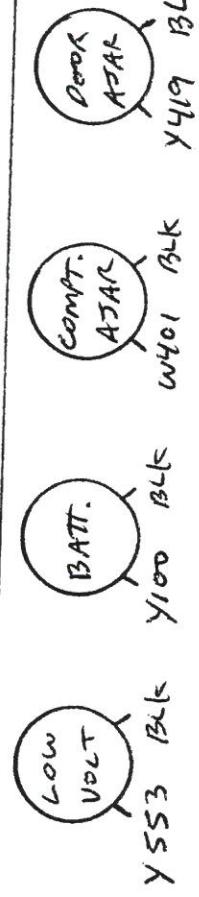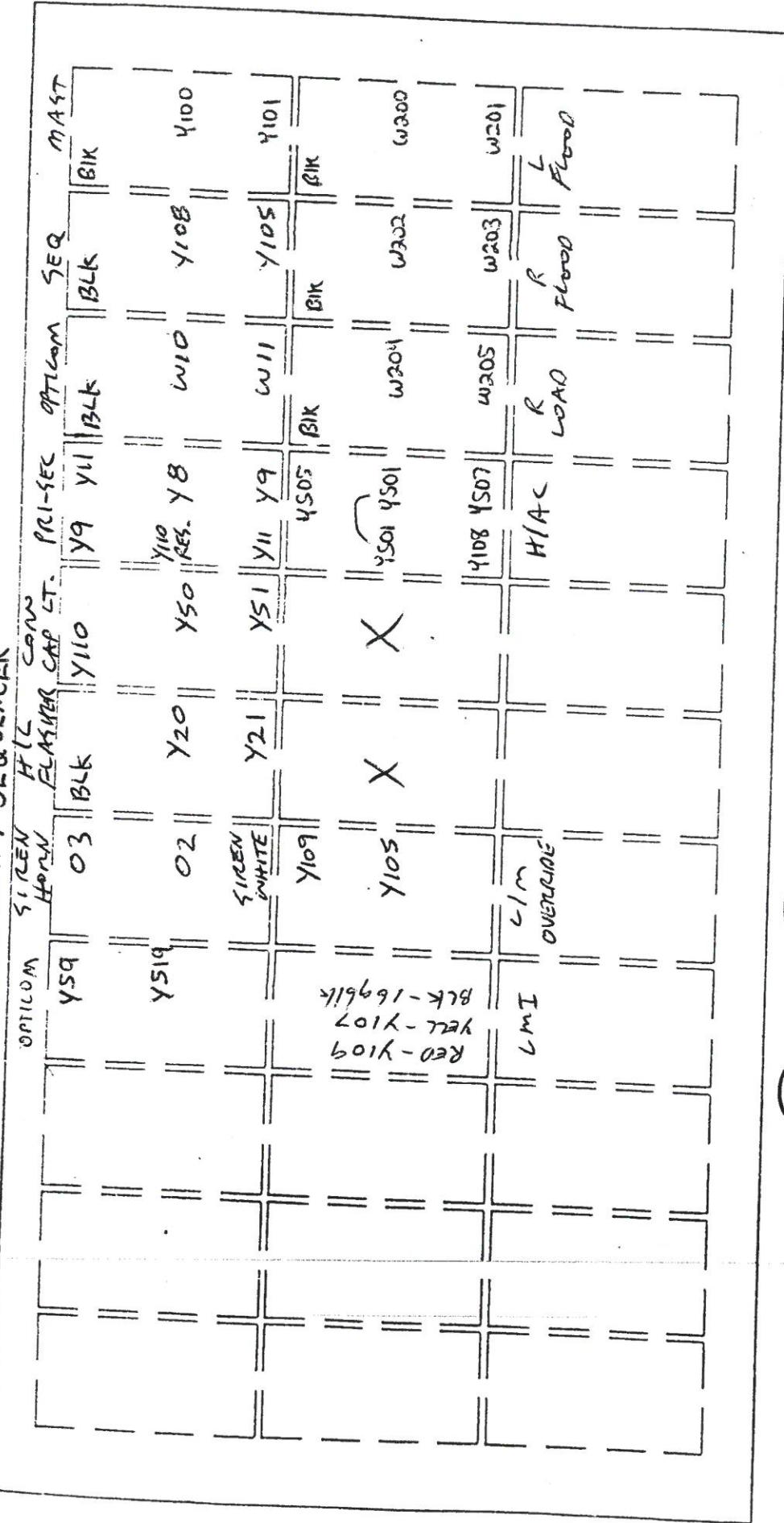

William S Oberender


William S Oberender  
Electrical Designer  
4/27/15  
Date


35303

MATERIAL — .050 BLACK VINYL CLAD ALUMINUM

M320654 — CHEVY DIESEL ECC  
CHEVY


353031

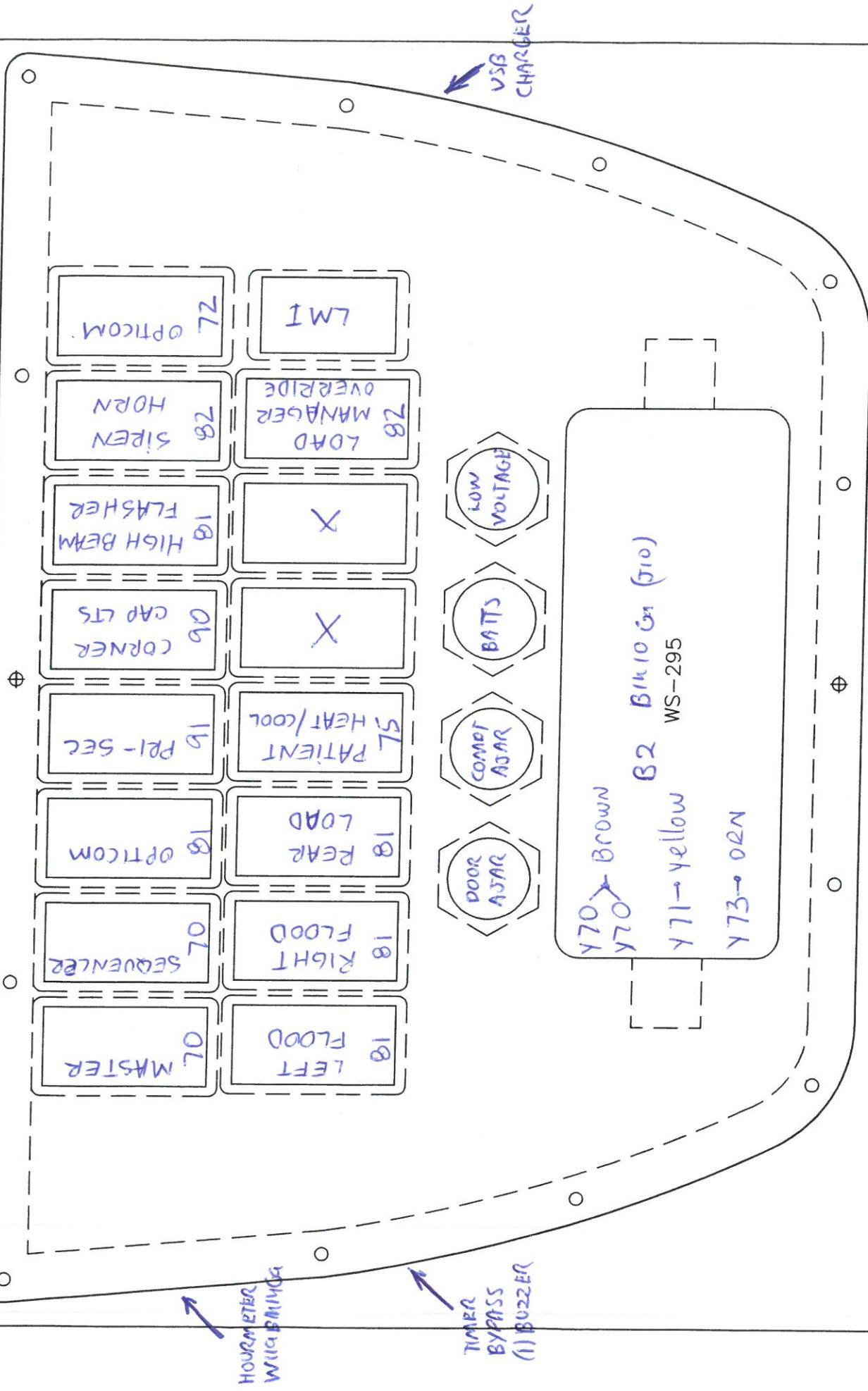

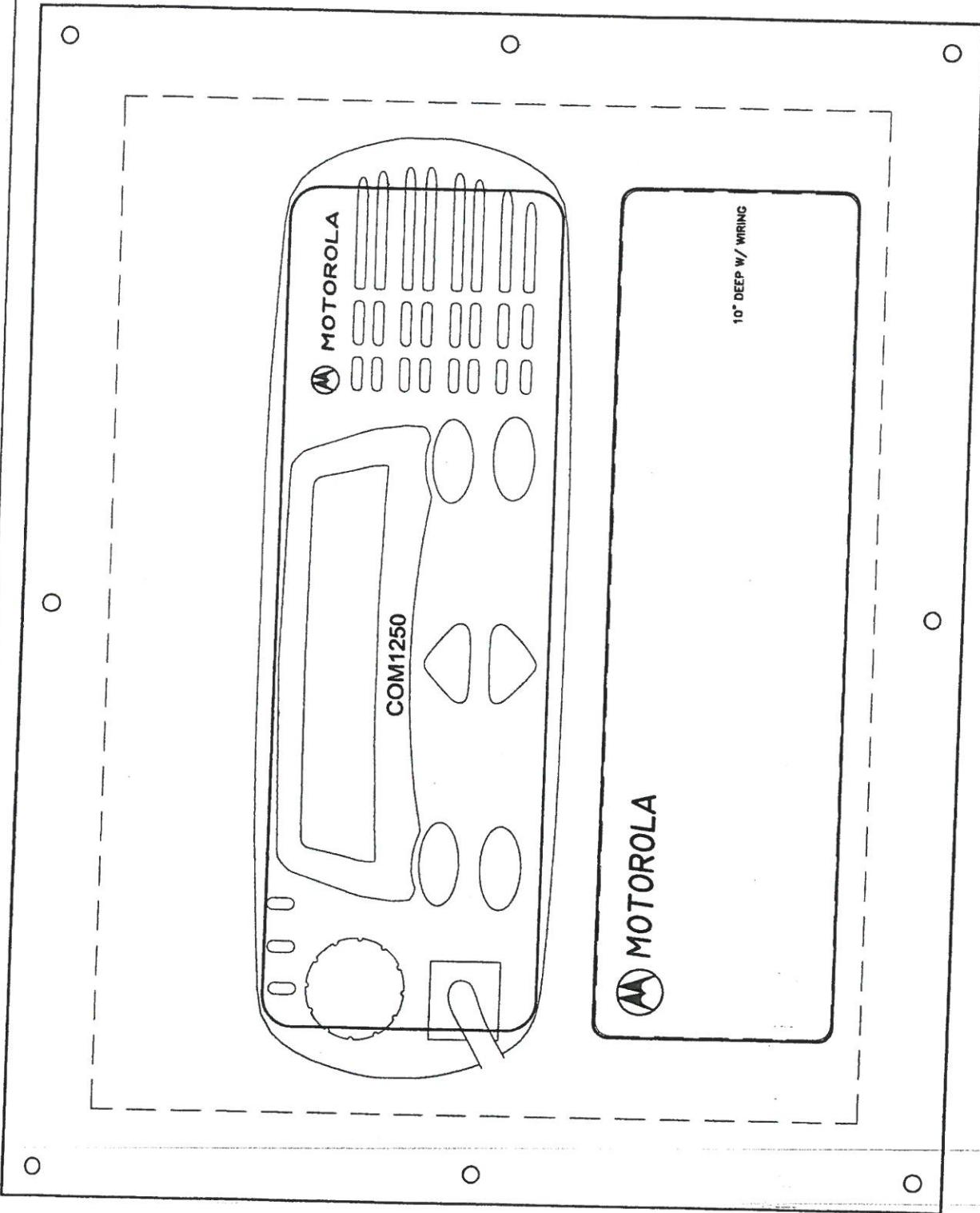

MATERIAL - .050 BLACK VINYL CLAD ALUMINUM  
M290094 - BEV - PNL - BEVERLY HILLS CONSOLE  
10.0 X 8.0  
353031

William Oberender  
Electrical Designer  
Date  
4/29/15

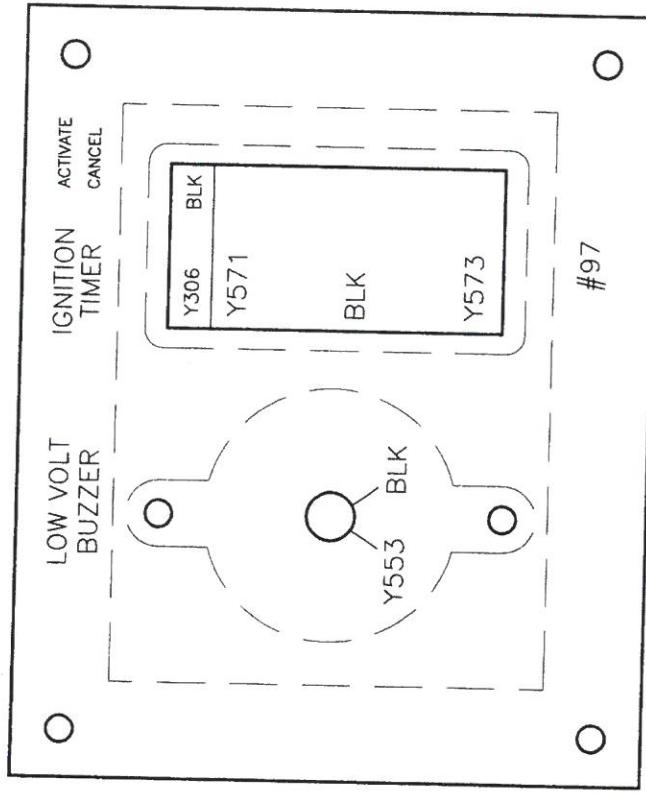

MOUNT ON DRIVERS SIDE OF ECC

MATERIAL - .050 BLACK VINYL CLAD ALUMINUM

353031

M290013 - ACCY - PNL  
4.00 X 3.25

William S. Oberender  
Electrical Designer  
Date

4/29/05

William L. Oberndorfer

4/29/15 Date




*The World's Largest Manufacturer of Ambulances*

REV	DATE	DESCRIPTION	STD A-A SW PNL WIRING	DRAWING NUMBER
REV	DATE	VEH. TYPE	OPTIONAL	
13FEB07	13FEB07	②		-99

150-3

W/O #-


DESCRIPTION

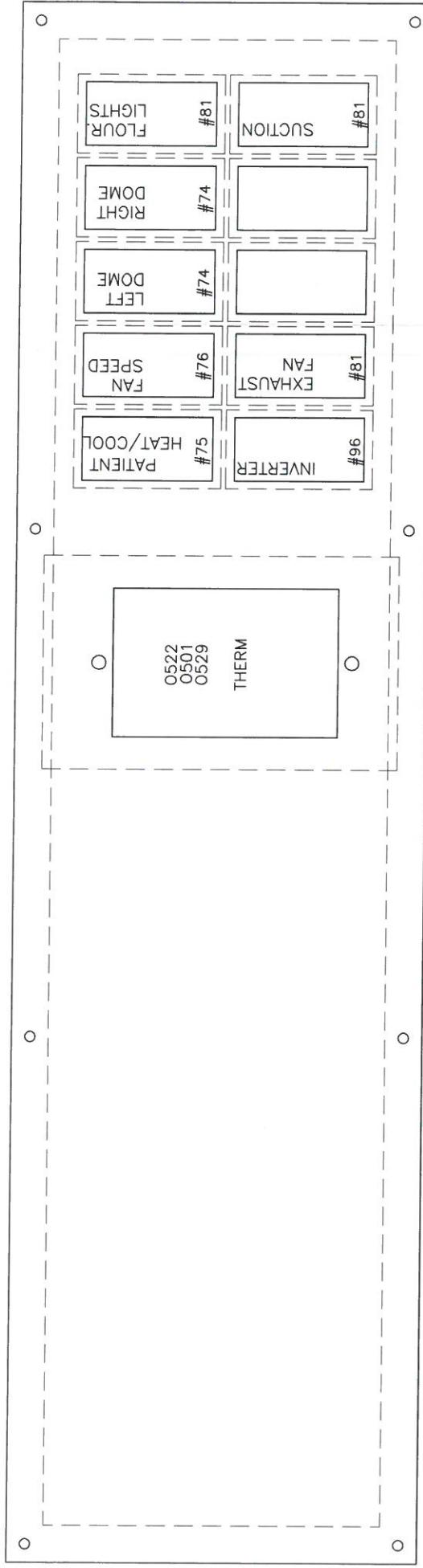

*MATERIAL - .050 BLACK VINYL CLAD ALUMINUM*

L240020- AAPNL MATERIAL .050 BLACK VINYL CLAD ALUMINUM  
SIZE - 21 X 5.5

353031

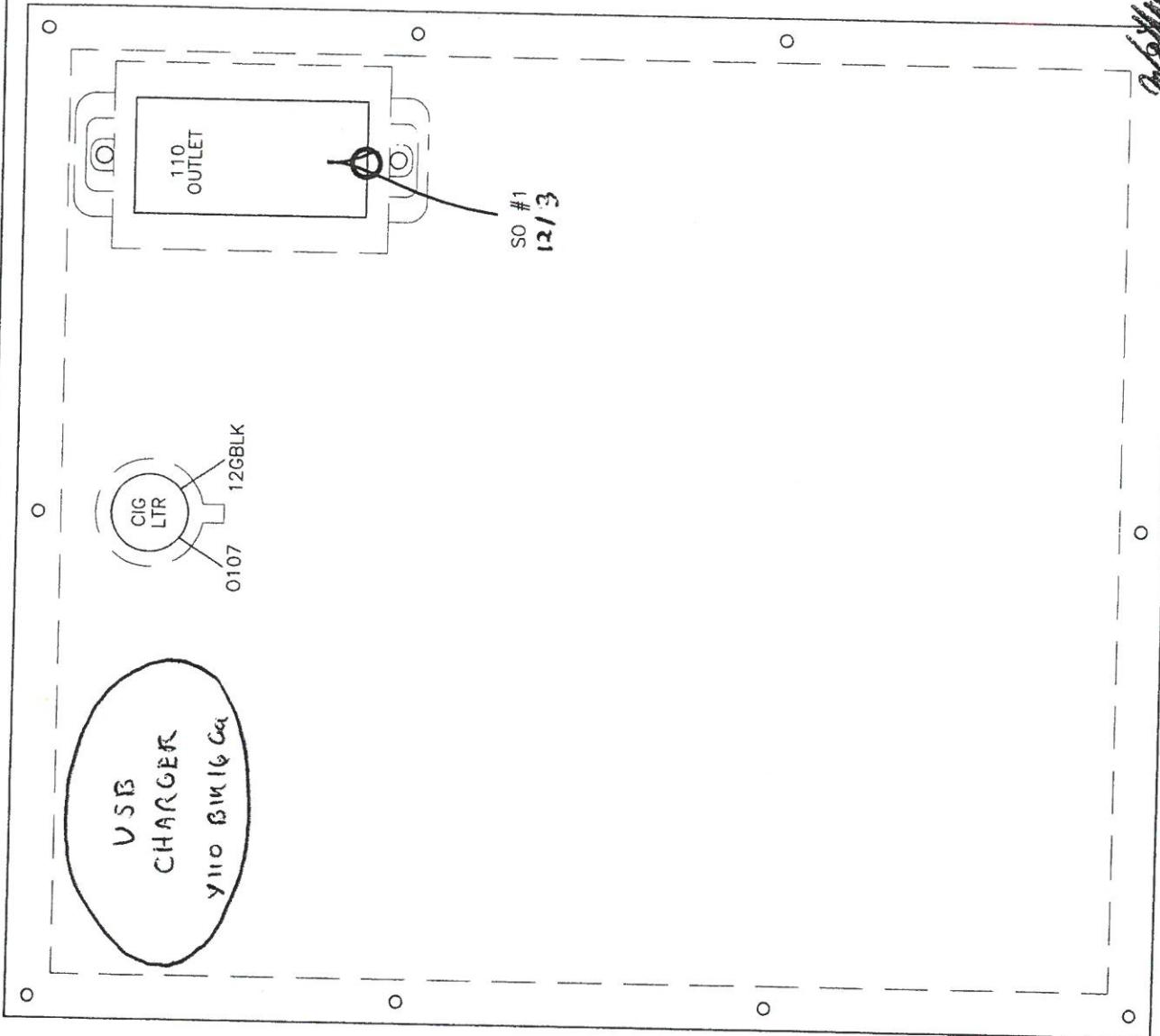


35303

L250104 - LP - (FOR AUTOEJECT)  
SIZE - 11.5 X 13

William Oberender  
Electrical Designer


4/29/15  
Date

William S Oberender

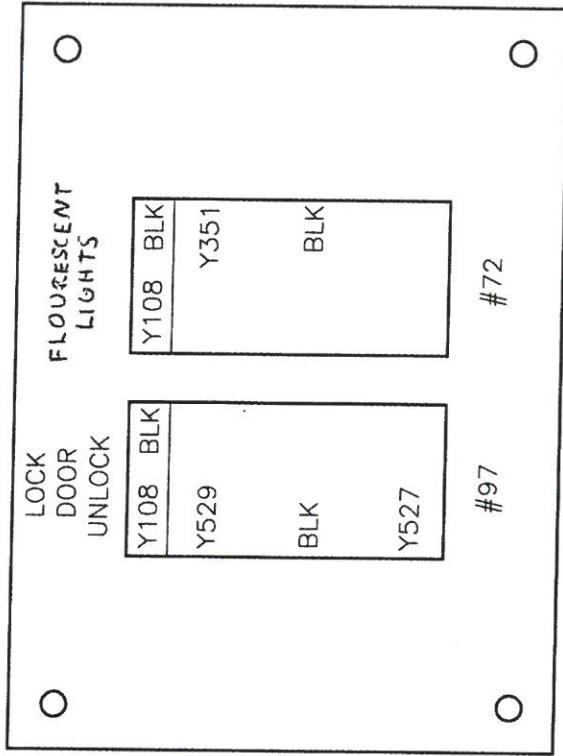

MOUNT ON C/S WALL DOWN LOW NEAR  
SIDE ENTRY DOOR

M81624  
4.125 X 3

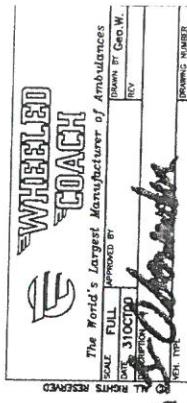

*William Oberender*

353031


William Oberender 4/29/15  
Electrical Designer Date


353031

M80705 - CIG. LIGHTER PANEL  
SIZE - 3 X 3

William S Oberender  
Electrical Designer 4/29/15  
Date

MOUNT IN R/F CABINET PER W/O

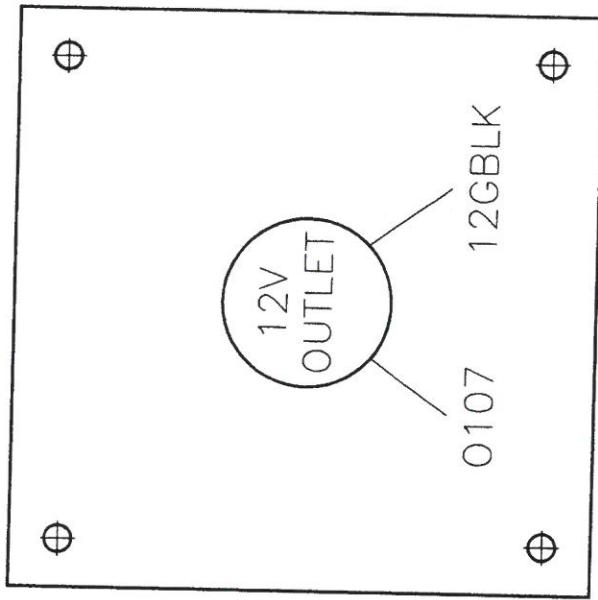

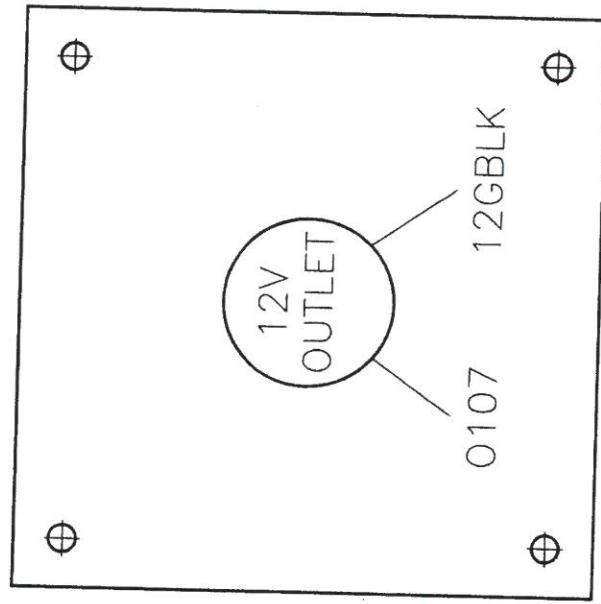


MOUNT IN EXTENDED A/A PER W/O

353031

M80705 - CIG. LIGHTER PANEL  
SIZE - 3 X 3

William Oberender  
Electrical Designer


4/29/15  
Date

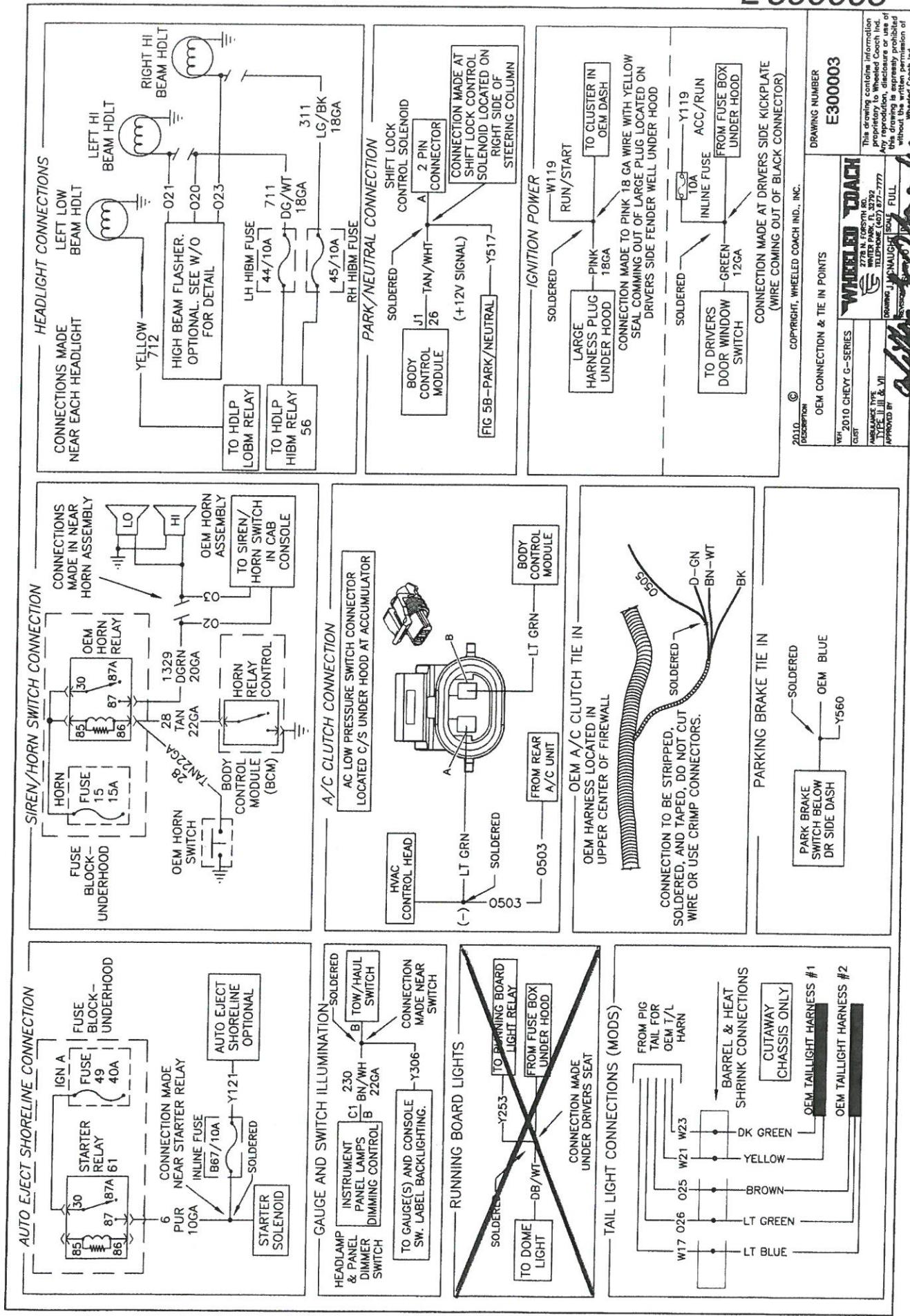

353031

William Oberender  
Electrical Designer

4/29/15  
Date

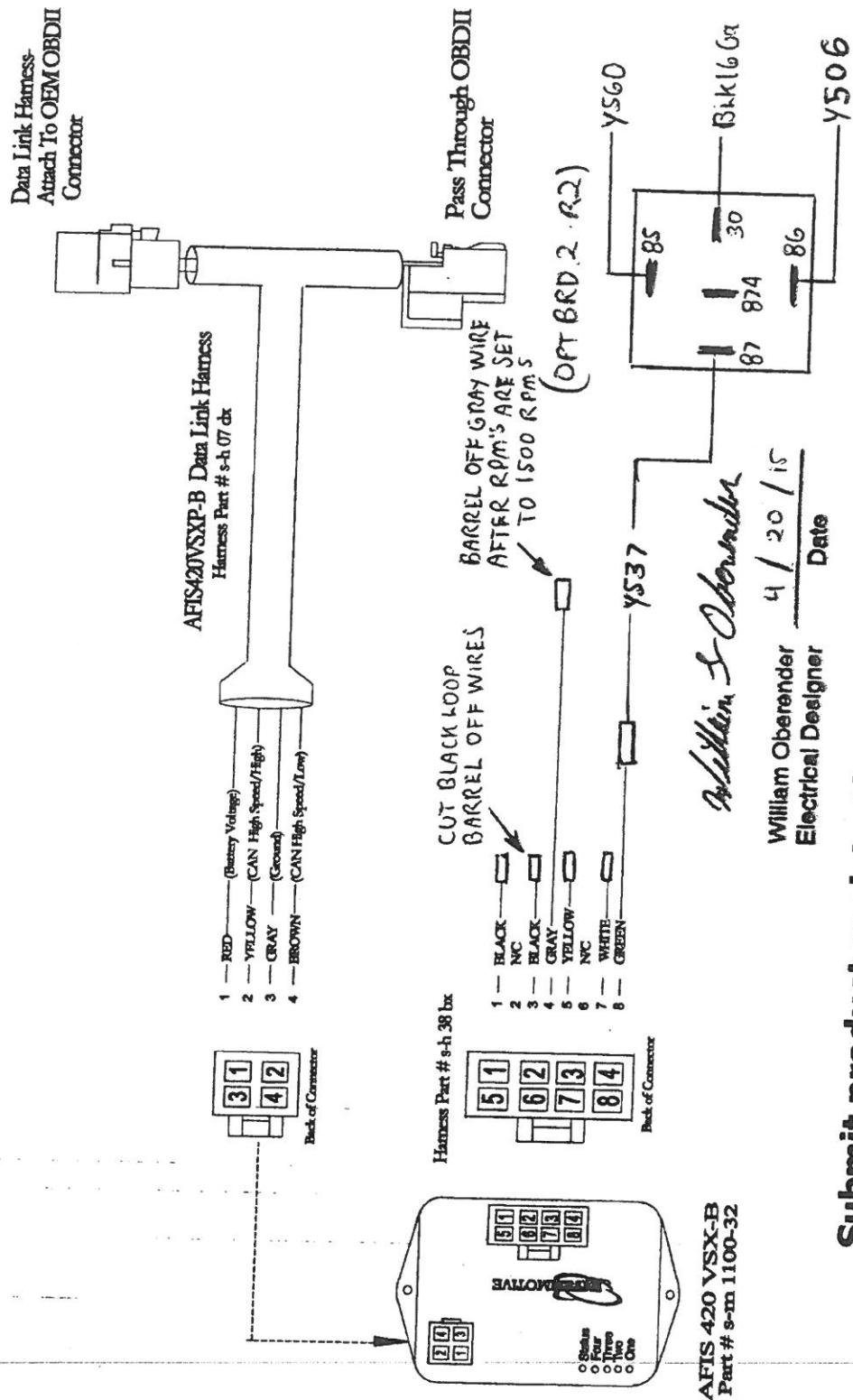


**Submit product registration at [www.intermotive.net](http://www.intermotive.net)**

If the AFIS420VSX-B fails any step in the Post Installation Test, review the installation instructions and check all connections.  
If necessary, call

353031

**InterMotive technical support @ (530) 823-1048.**

AFIS420VSX-B-02-CAD

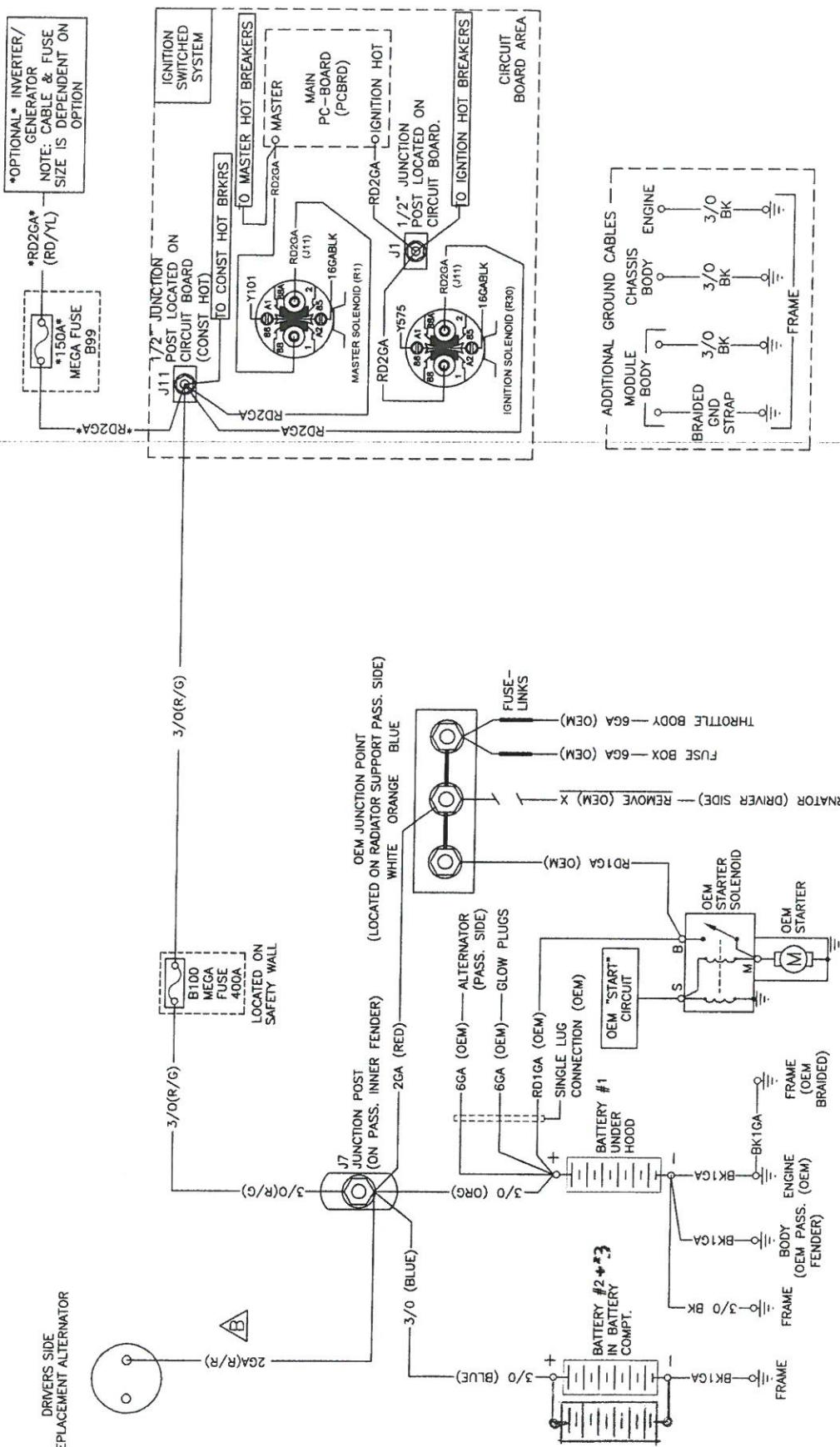

© 2010 COPYRIGHT, WHEELED COACH IND., INC.

DRAWING NUMBER  
E 3000

DRAWING NUMBER  
**E 3000**  
This drawing contains  
proprietary to Wheeler  
Any reproduction,  
of this drawing is  
prohibited without  
permission of Wheeler

**DRAWING NUMBER**

10

E3000

L'Uomo

This document contains

Proportionality to Wholes

Any reproduction, disclosure  
or distribution is illegal.

prohibited without

Copyright of Wheeler

1

10


三


Oberender - 4/29/15 - Dm

William Oberender  
Electrical Designer

333031

E 300099

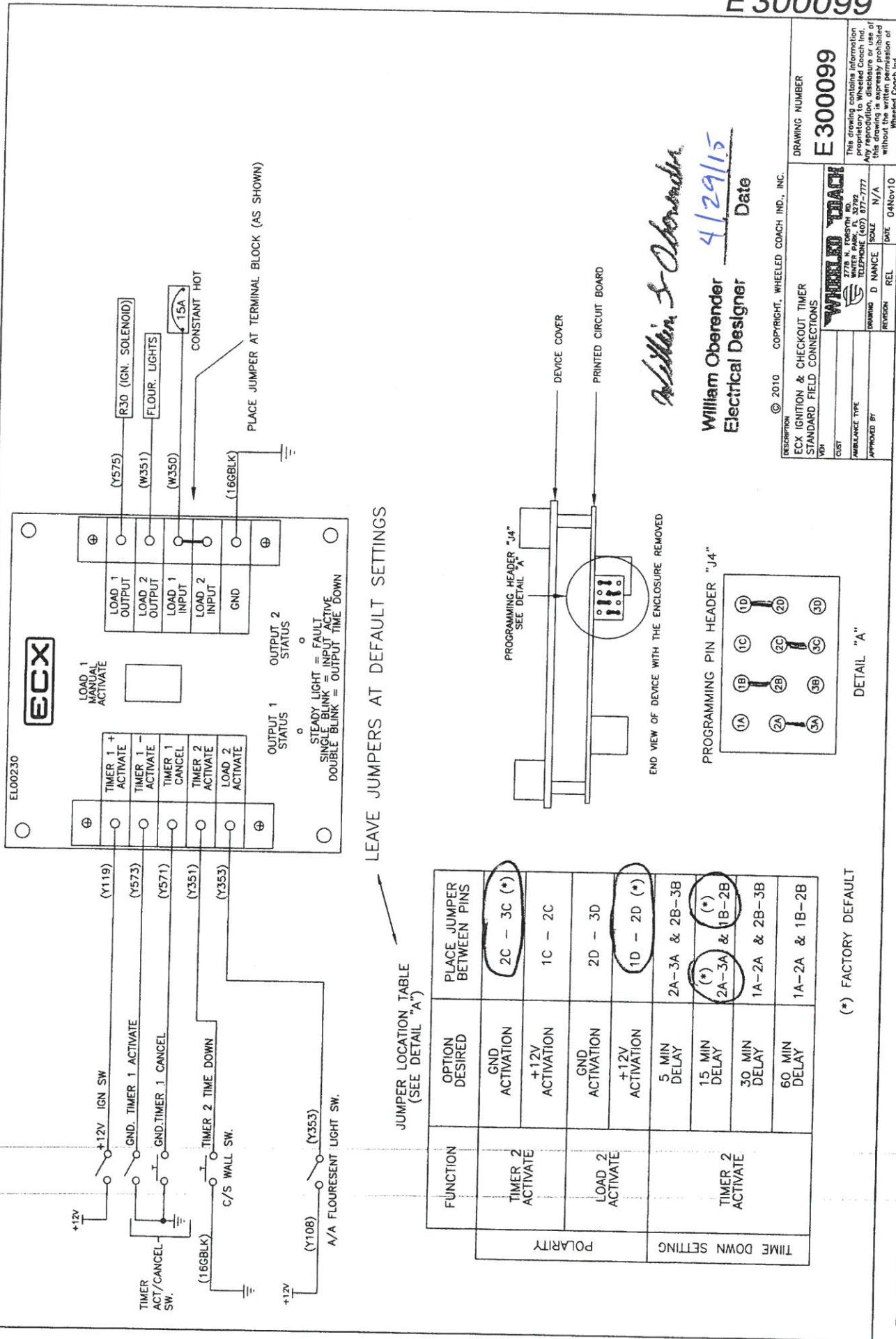

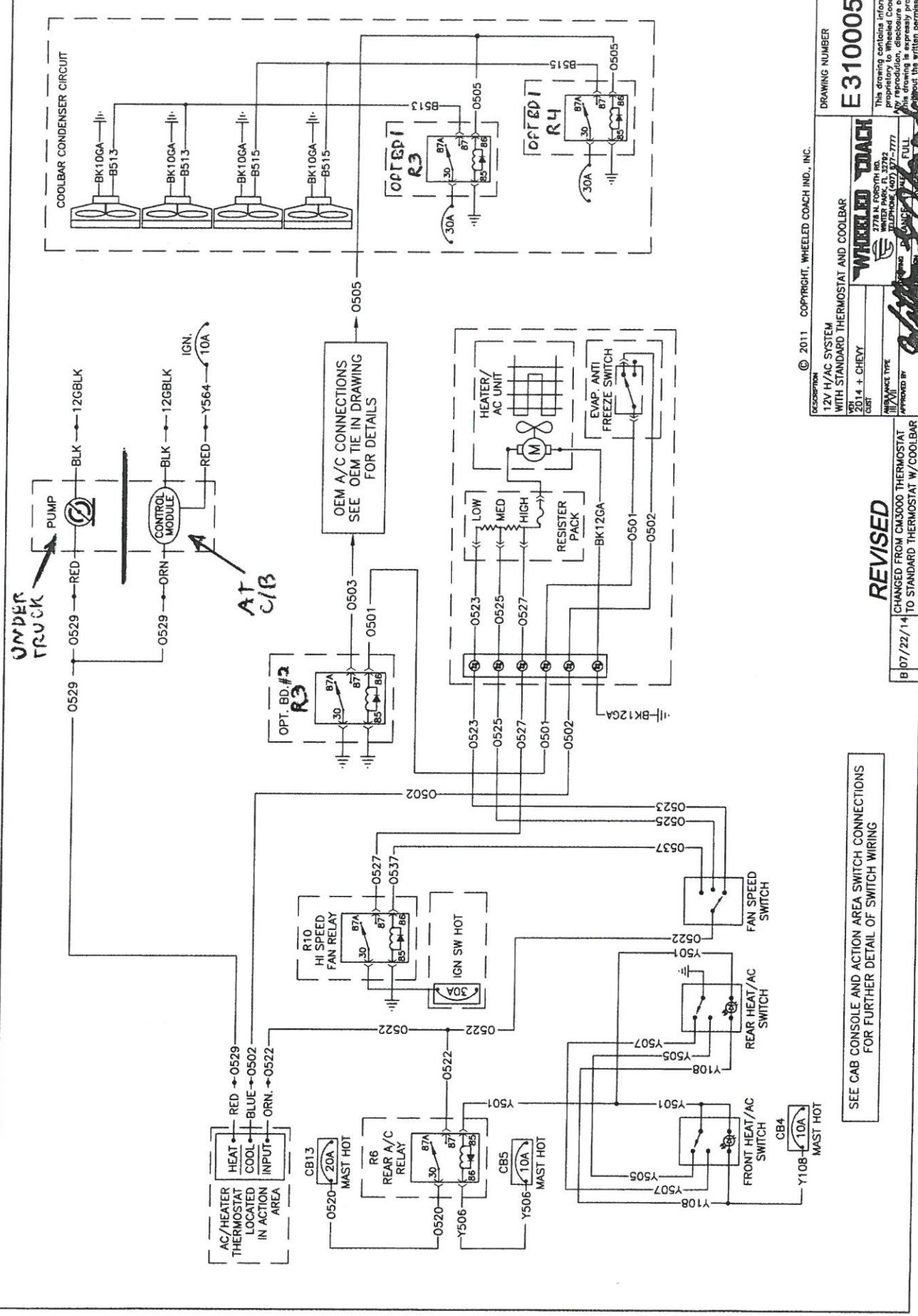


© 2011 COPYRIGHT, WHEELED COACH IND., INC.

E 310005B  
WHEELED COACH  
12V H/AC SYSTEM  
WITH STANDARD THERMOSTAT AND COOLBAR  
2014 + CHEVY  
CLST

Proprietary information of Wheeled Coach Inc.  
No reproduction, disclosure or use of this drawing is expressly prohibited  
without the written permission of Wheeled Coach Inc.

William Oberender  
Electrical Designer  
4/29/15  
Date

353031

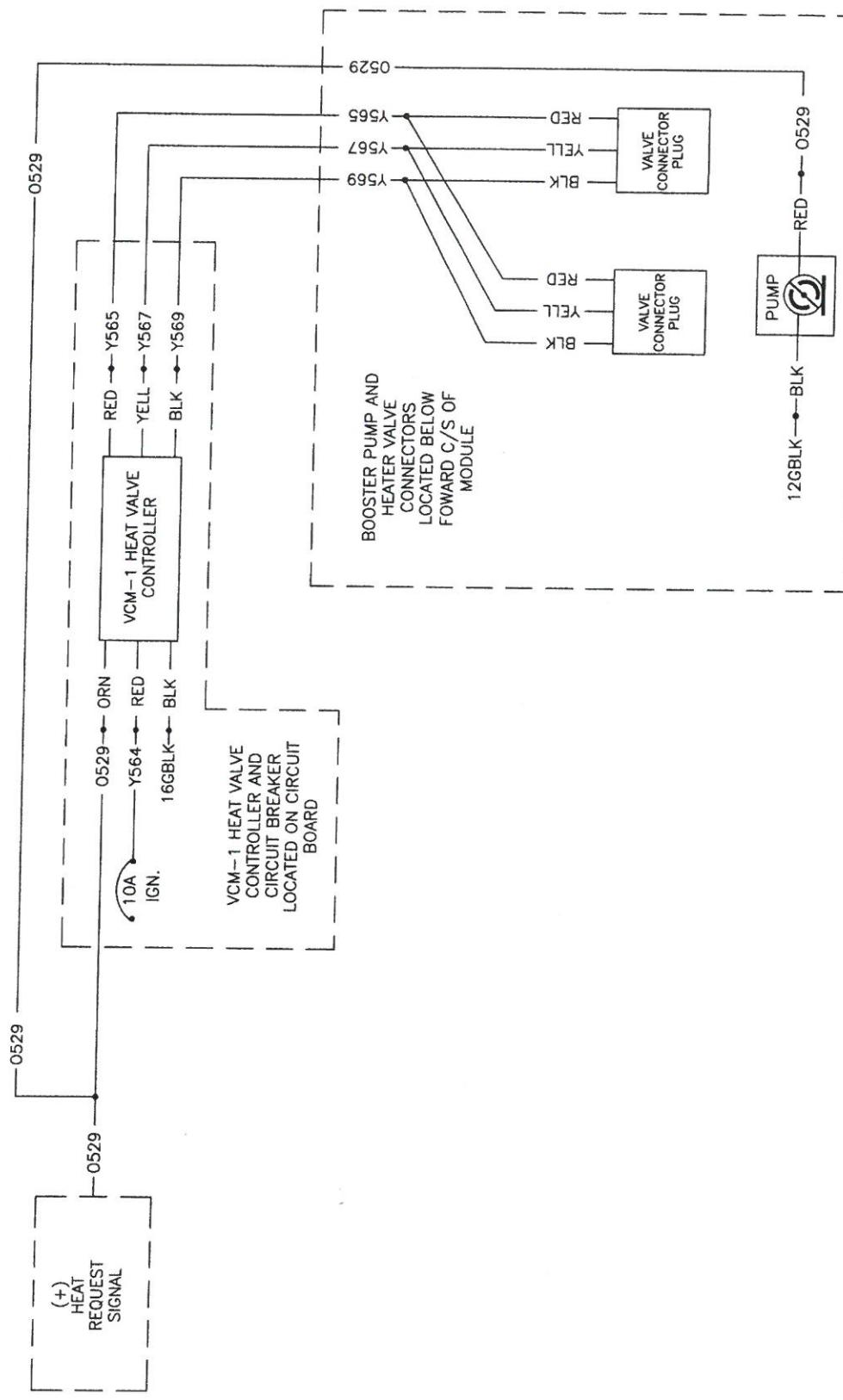

© 2014 COPYRIGHT, WHEELED COACH IND., INC.

DRAWING NUMBER

E 350012

This drawing contains information  
pertaining to Wheeled Coach Inc.  
Proprietary to Wheeled Coach Inc.  
Unauthorized disclosure or use of  
this drawing is expressly prohibited  
by law. All rights reserved.  
The drawing is a complete drawing  
without the written permission of  
Wheeled Coach Inc.



WHEELED COACH

2778 N. FORESTON RD.

WHEELER PARK, FL 32792

TELEPHONE (407) 877-7777

FAX (407) 877-7777

FULL

DRAWING DATE


APPROVED BY

INITIALS

4/29/15

Date


353031

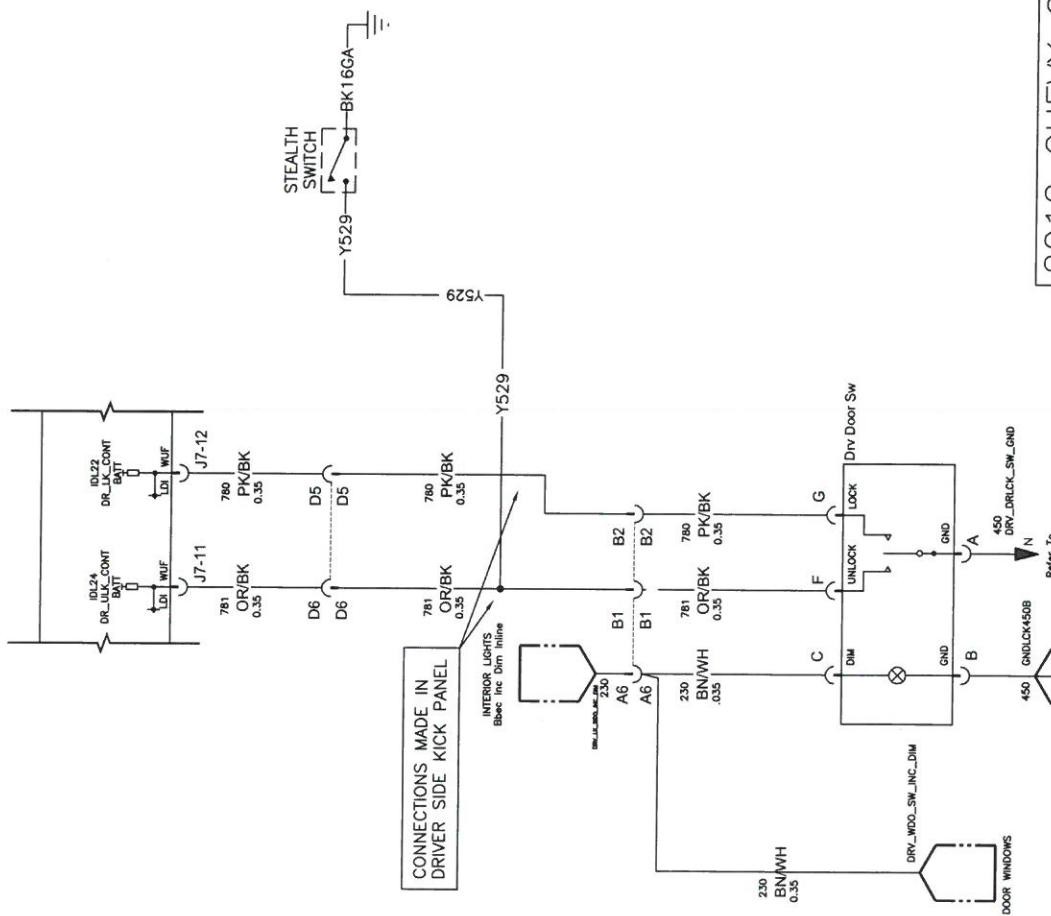

© 2012 COPYRIGHT, WHEELED COACH IND., INC.	DRAWING NUMBER E320017
DESCRIPTION ELECT. DOOR LOCK SYSTEM WITH KEYLESS REMOTE	
WHEELED COACH	2778 N. FOREST RD. WINTER PARK, FL 32792 TELEPHONE: (407) 677-7777
DR. 2012 CHEVY CUTAWAY CAB	AMBIANCE TYPE III
APPROVED BY	SCALE FULL
REVISION	DATE 07 Aug 12

2012 CHEVY CUTAWAY CHASSIS

353031  
GND DIST GND  
OEM DOOR LOCK (DRIVER SIDE)

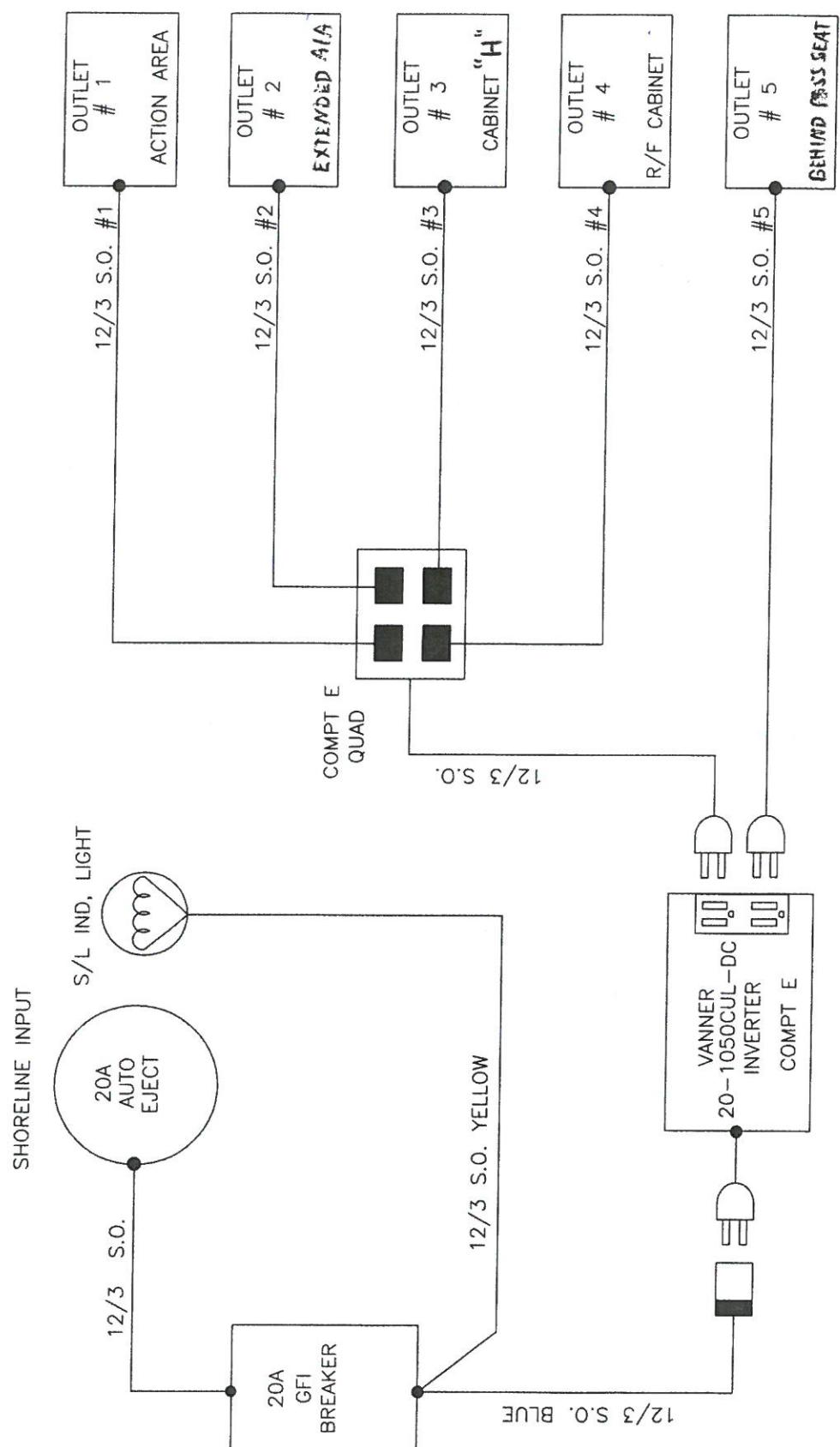


*William L Obernder*

William L Obernder  
Electrical Designer

4/29/15  
Date

353031

E 290378A

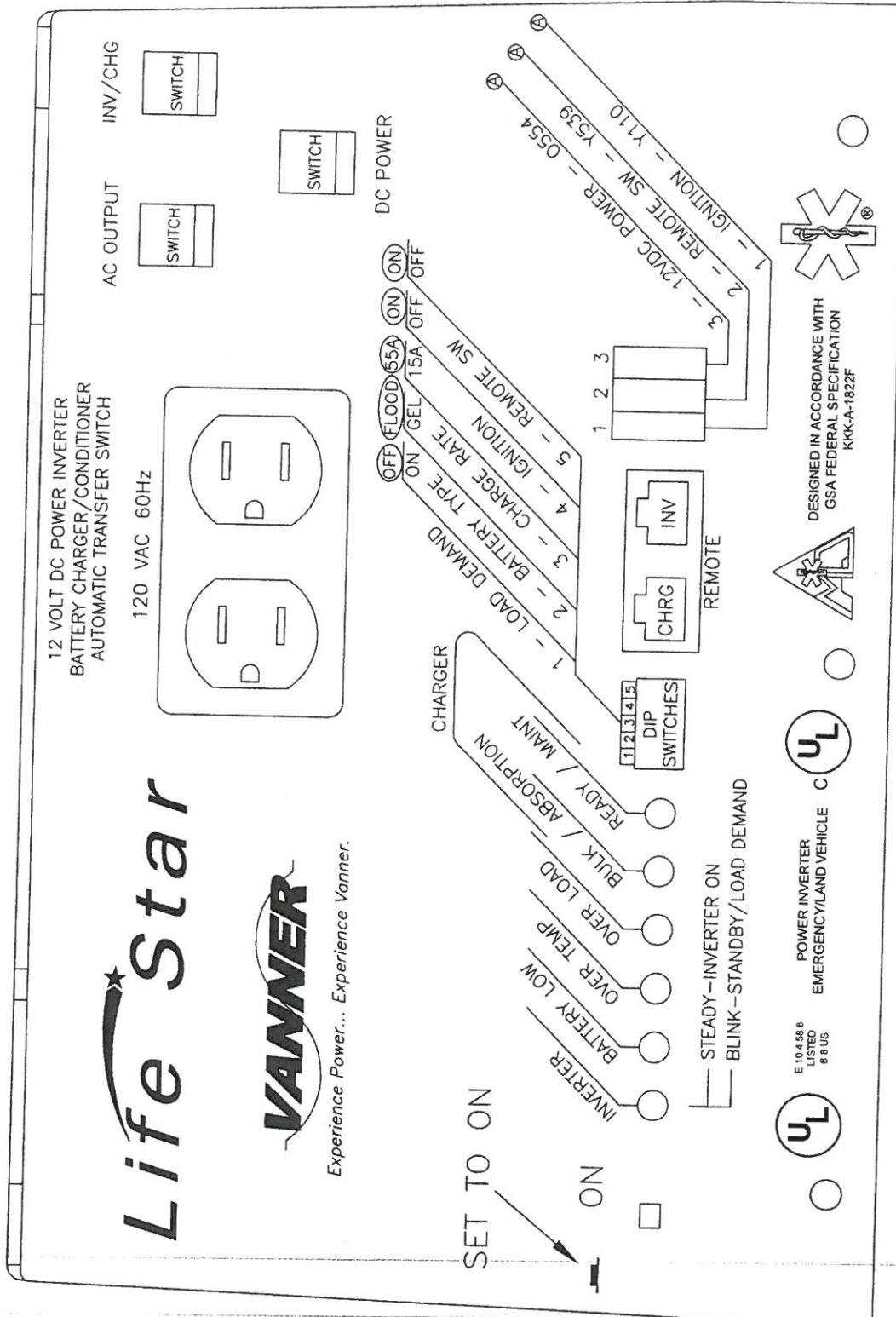

© 2010	COPYRIGHT, WHEELED COACH IND., INC.	DRAWING NUMBER
WHEELED COACH		E 290378A
This drawing contains information proprietary to Wheeled Coach Inc. Any reproduction, disclosure or use of this drawing is expressly prohibited without the written permission of Wheeled Coach Inc.		
DESCRIPTION	LIFE STAR VANNER 20-1050 INVERTER/BATTERY CHARGER VEN	VEN
CLASS	2778 N. MICHIGAN RD. TELEPHONE: 404-677-7777	TELEPHONE: 404-677-7777
MANUFACTURE TYPE	DRAWING C. ELLIS	SCALE: N/A
DISPOSITION	FOR INFORMATION	REL DATE 01 March 10


**REVISED**

④ 03/30/10	CHANGED THE LABELING OF THE REMOTE CONNECTOR (C) <i>John</i>
------------	---

353031


William Oberender  
Electrical Designer

9/29/15  
Date

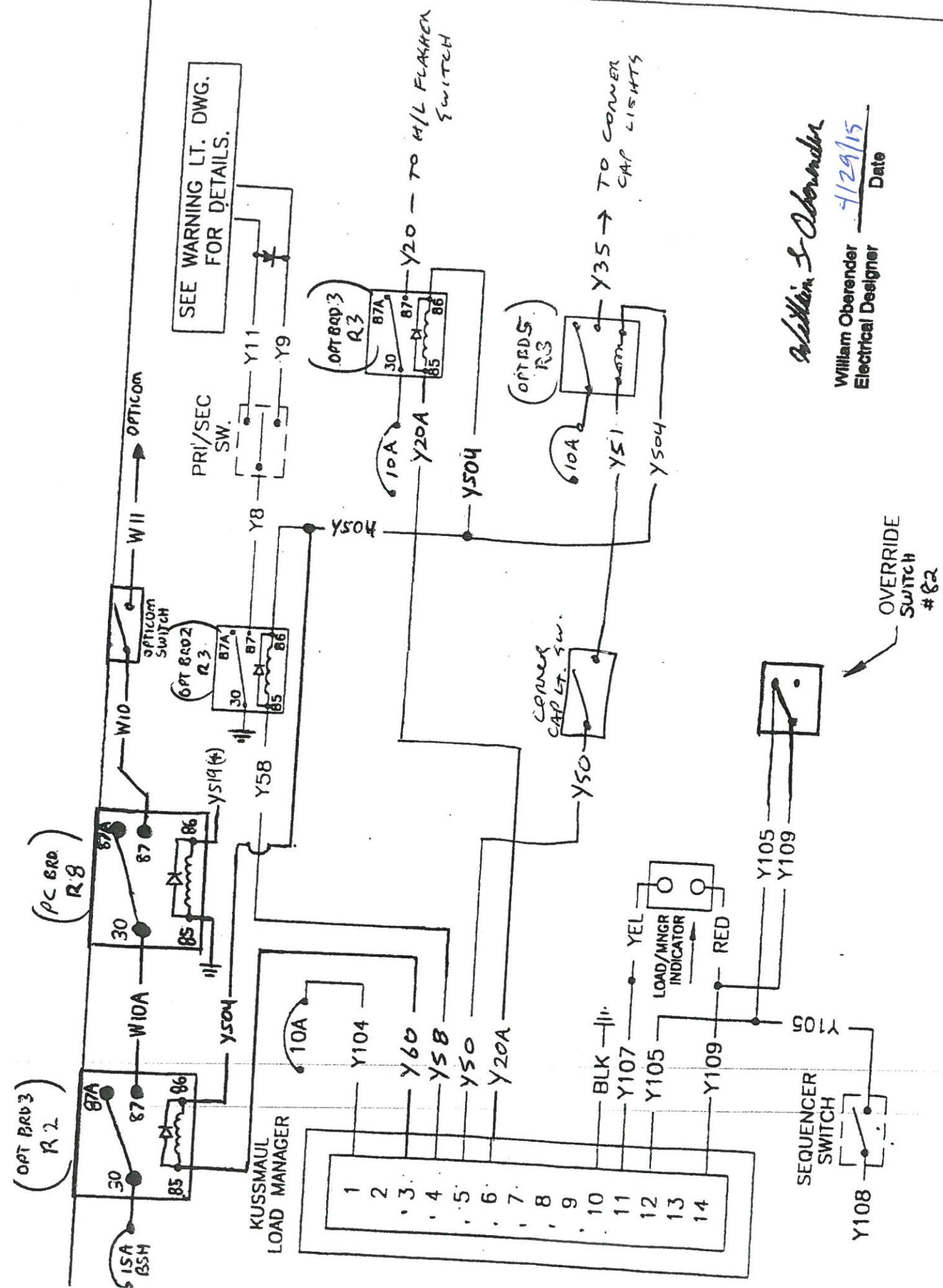

53031

GND FIRED PRI / SEC SWITCHING - (WITH OPTICOM)

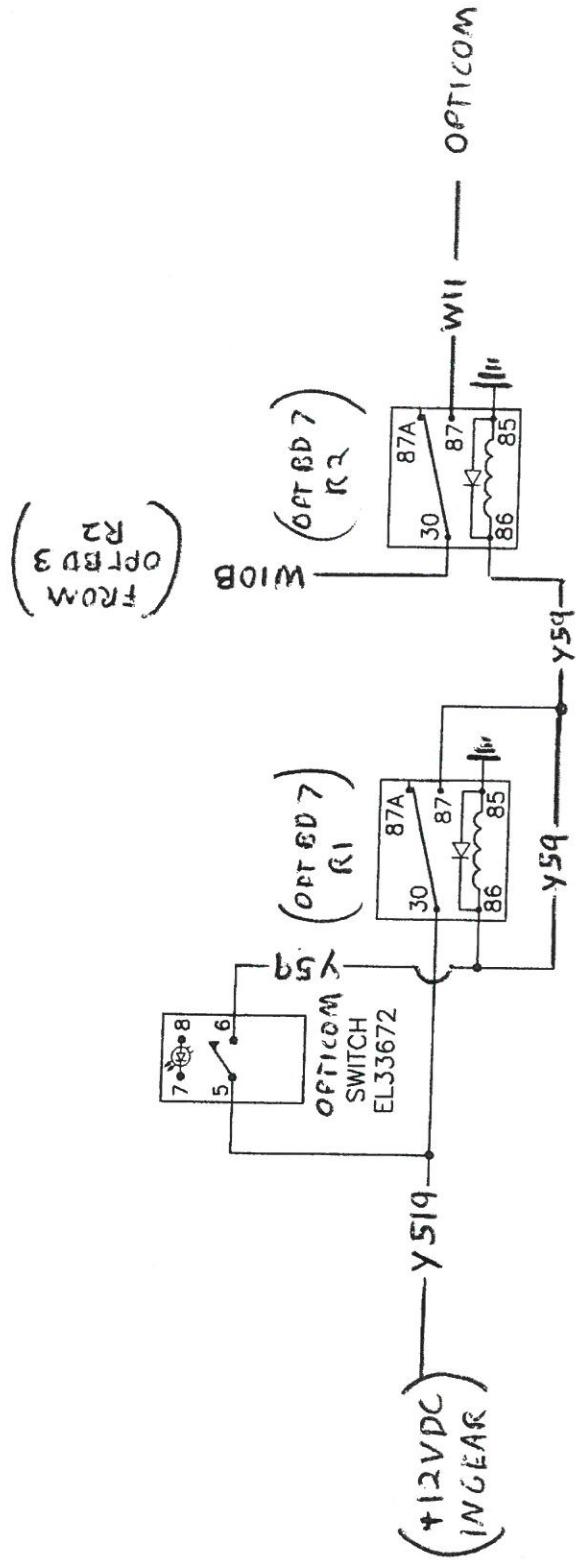


IGNITION HOT (1) CYCLE OPTICOM CIRCUIT  
 OPTICOM WILL BE SHUT OFF IN P/N.  
 SWITCH WILL TURN OPTICOM ON IN GEAR FOR (1) USE  
 WILL REMAIN ON UNTIL SHIFTING INTO PARK

(NOTE: REFERENCE ONLY. WIRE NUMBERS  
 REFLECT MODIFICATIONS TO A CERTAIN VEH.)

*William L. Oberender*

353031


William Oberender 4/29/15  
 Electrical Designer  
 Date

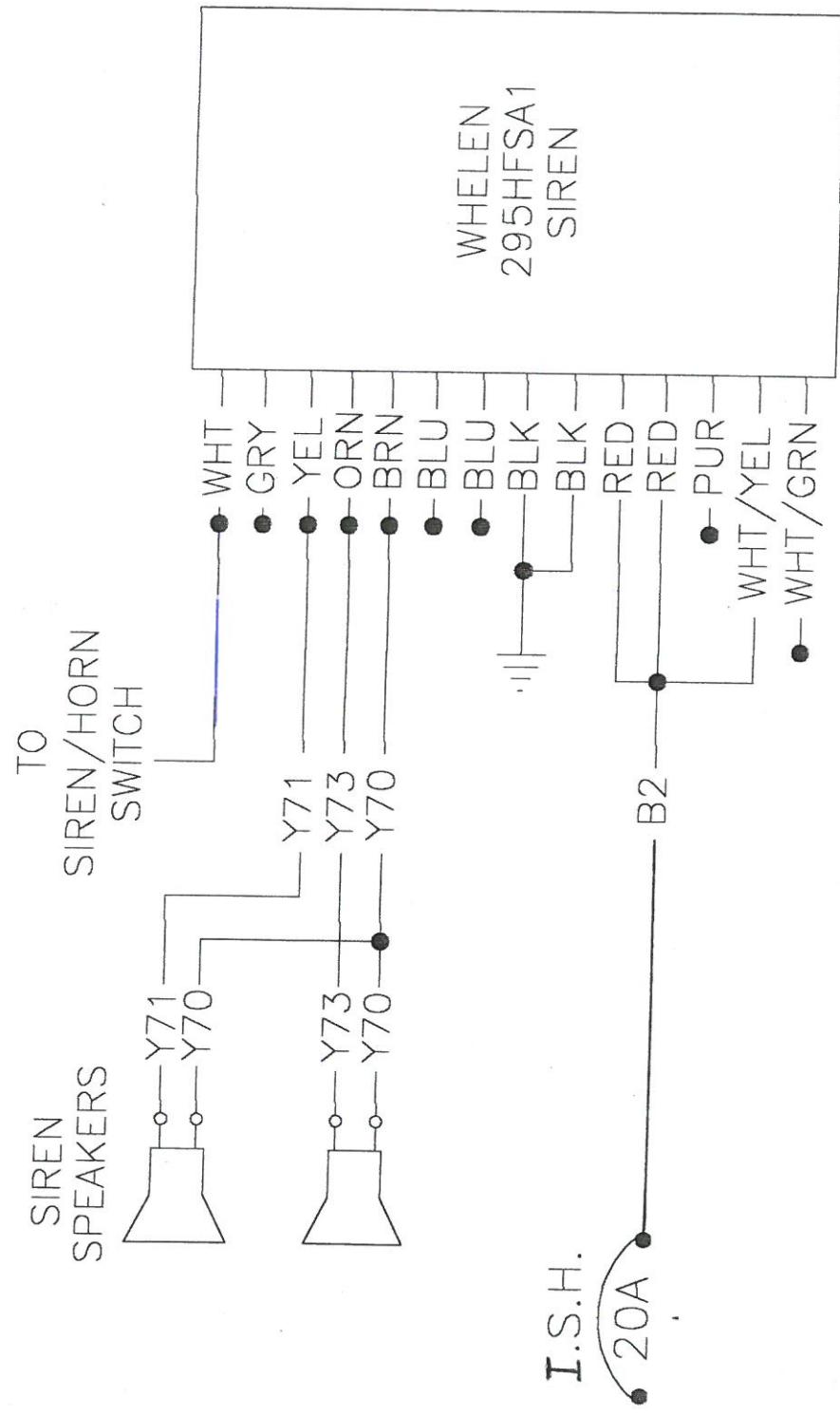

393031

STANDARD WS295-HFSA1 CONNECTIONS

CONNECTIONS MADE  
NEAR EACH HEADLIGHT

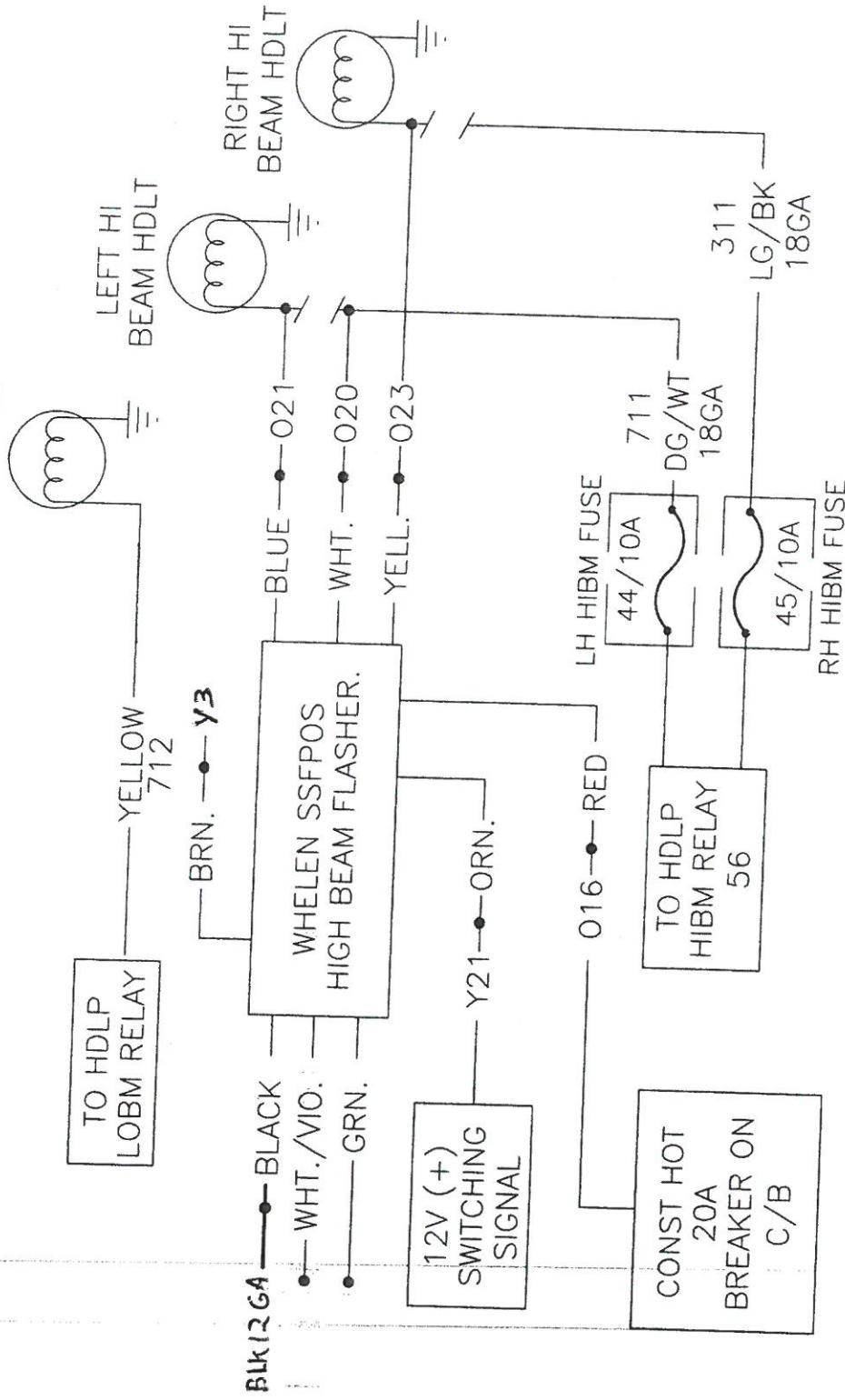

© 2012 COPYRIGHT, WHEELED COACH IND., INC.

E 320016

WHEELED COACH


DESCRIPTION	ITEM NO.	QUANTITY	UNIT	DRAWING NUMBER
WHELEN SSFPOS HEAD LIGHT FLASHER	2012 CHEVY	1	PCB	E 320016
COST				
MANUFACTURE TYPE	II AND III			
APPROVED BY	C. ELLIS	N/A		
REVISION	REL	DATE		
		28Jun12		

This drawing contains information  
proprietary to Wheeled Coach Ind.  
Any reproduction, disclosure or use  
of this drawing is expressly prohibited  
without written permission of  
Wheeled Coach Ind.

353031

William Oberndorfer  
Electrical Designer 412915  
Date


E290266

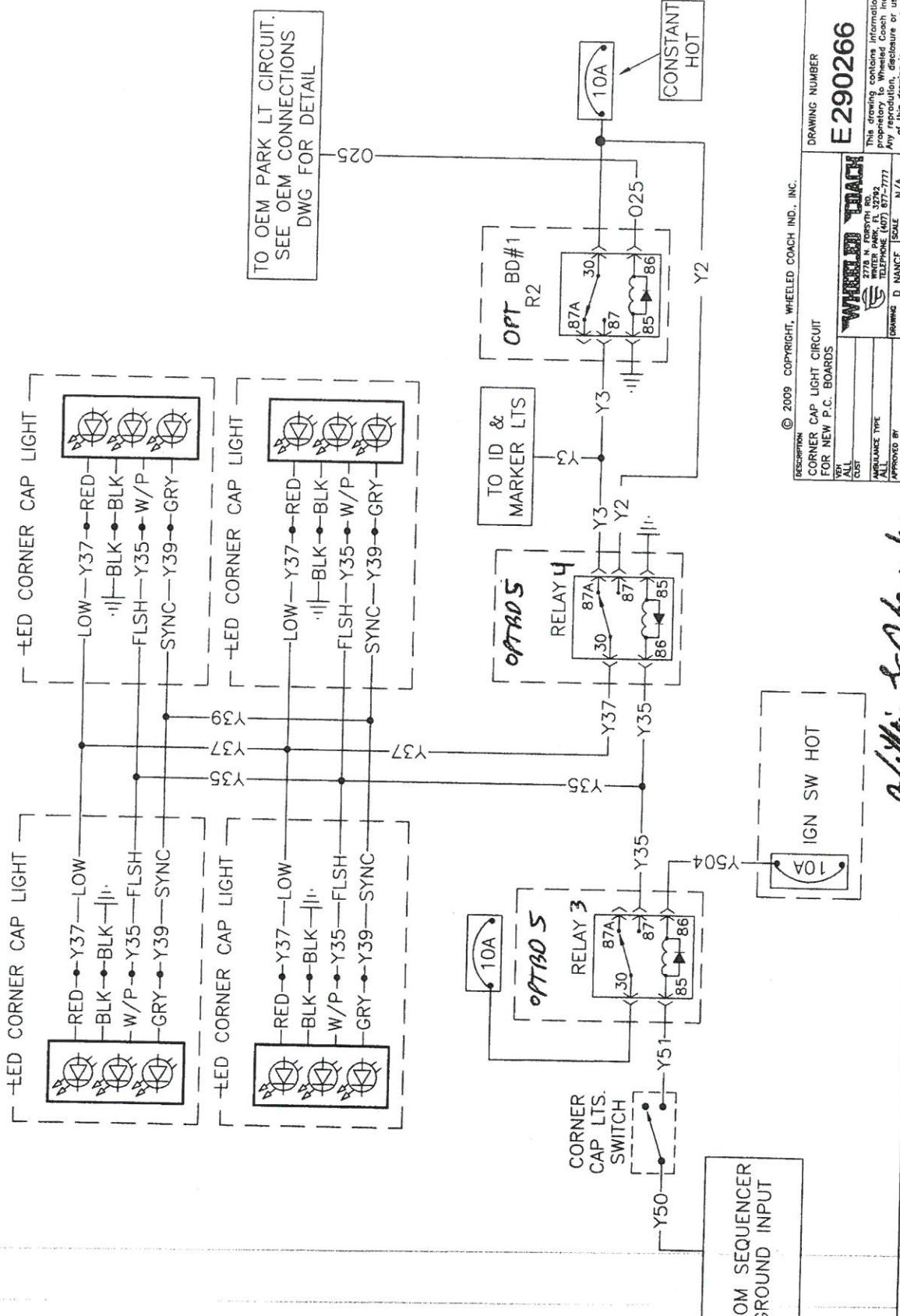

© 2009 COPYRIGHT, WHEELED COACH IND., INC.	DRAWING NUMBER E290266
This drawing contains information proprietary to Wheeled Coach Ind. Any reproduction, disclosure or use of this drawing is expressly prohibited without the written permission of Wheeled Coach Ind.	
DESCRIPTION CORNER CAP LIGHT CIRCUIT FOR NEW P.C. BOARDS TER ALL CLASS	WHEELED COACH 2778 N STORTH RD. WINTER PARK, FL 32792 TELEPHONE (407) 877-7777
APPLIANCE TYPE ALL	DRAWING D NANCE SCALE N/A
APPROVED BY	REVISION REL
	DATE 05/05/09


353031

William S Oberender  
Electrical Designer  
Date  
4/29/15

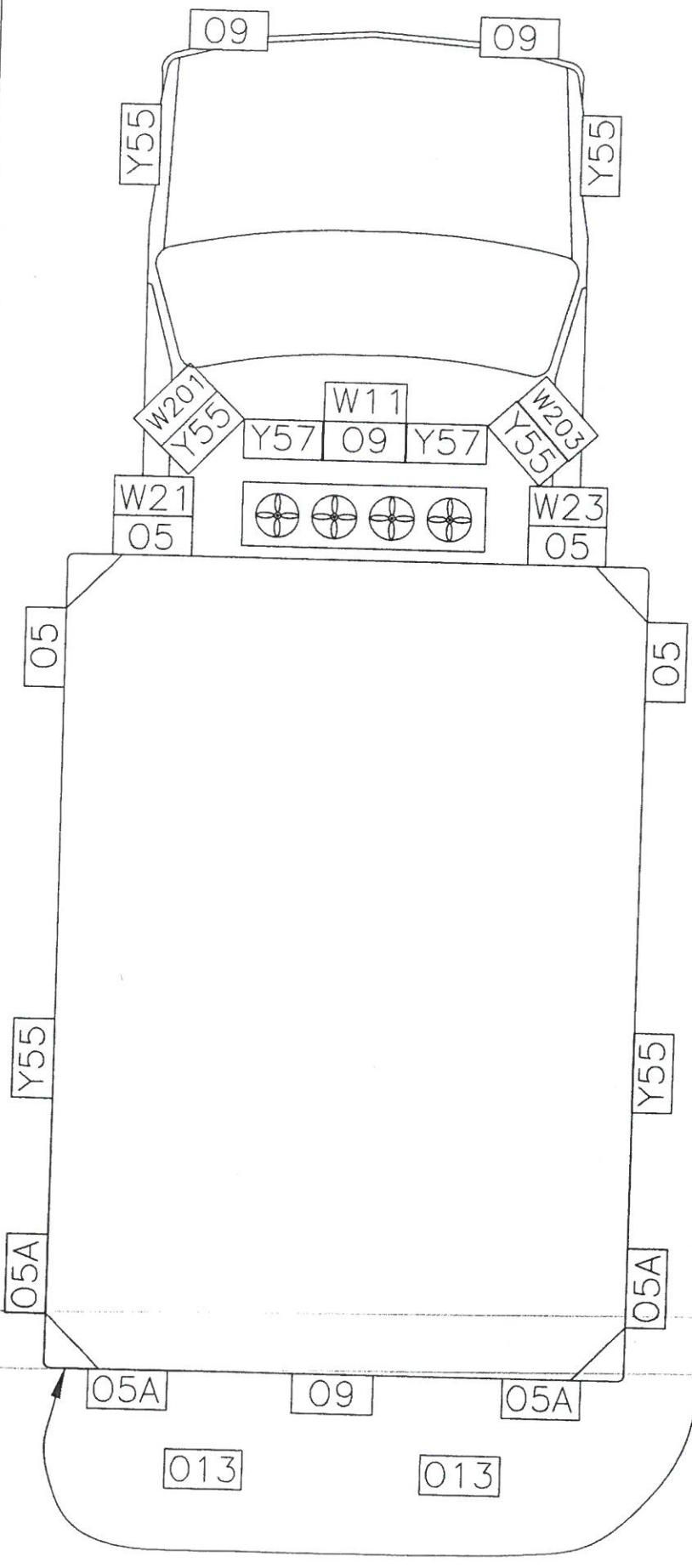

### OEM TAILLIGHT CONNECTIONS FOR CHEVY CUTAWAY

FROM PIG  
TAIL FOR  
OEM T/L  
HARN

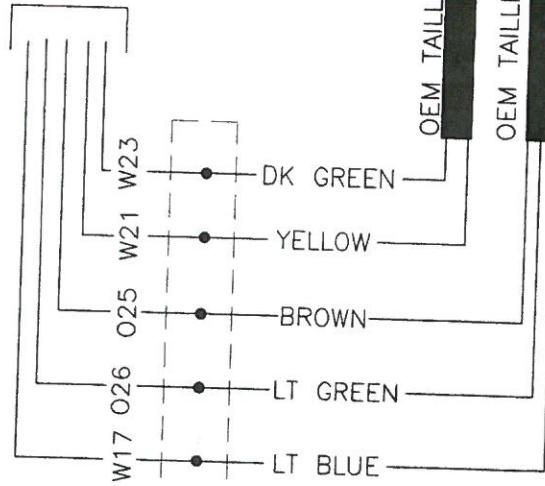

BARREL CONNECTION

35303 |

*William S. Oberender*  
 William Oberender  
 Electrical Designer  
 4/29/15  
 Date